

# EUROPEAN SOUTHERN OBSERVATORY

Organisation Européenne pour des Recherches Astronomiques dans l'Hémisphère Austral

Europäische Organisation für astronomische Forschung in der südlichen Hemisphäre

# VERY LARGE TELESCOPE

# FLAMES-UVES Pipeline User Manual

VLT-MAN-ESO-19500-3016

Issue 5.0

Date 2008-11-28

Prepared:	Andrea Modigliani and Jonas Møller Larsen	2008-11-27
	Name	Date
		Signature

Approved: P.Ballester

---

Name Date Signature

Released:	M. Peron		
	Name	Date	Signature

This page was intentionally left blank

<b>ESO</b>	<b>FLAMES-UVES Pipeline User Manual</b>	Doc:	VLT-MAN-ESO-19500-3016
		Issue:	Issue 5.0
		Date:	Date 2008-11-28
		Page:	3 of <a href="#">110</a>

#### Change record

Issue/Rev.	Date	Section/Parag. affected	Reason/Initiation/Documents/Remarks
5.0	28/11/2008		Documentation describing CPL based release 4.3.0

This page was intentionally left blank

<b>ESO</b>	<b>FLAMES-UVES Pipeline User Manual</b>	Doc:	VLT-MAN-ESO-19500-3016
		Issue:	Issue 5.0
		Date:	Date 2008-11-28
		Page:	5 of 110

## Contents

<b>1</b>	<b>Introduction</b>	<b>11</b>
1.1	Purpose . . . . .	11
1.2	Acknowledgements . . . . .	11
1.3	Scope . . . . .	11
1.4	Reference documents . . . . .	12
<b>2</b>	<b>Overview</b>	<b>13</b>
<b>3</b>	<b>What's new in pipeline release 4.3.0</b>	<b>14</b>
<b>4</b>	<b>FLAMES Instrument Description</b>	<b>15</b>
4.1	Instrument overview . . . . .	15
<b>5</b>	<b>Quick start</b>	<b>18</b>
5.1	FLAMES-UVES pipeline recipes . . . . .	18
5.2	An introduction to Gasgano and EsoRex . . . . .	18
5.2.1	Using Gasgano . . . . .	19
5.2.2	Using EsoRex . . . . .	24
5.3	Example of data reduction using EsoRex . . . . .	26
<b>6</b>	<b>Known problems</b>	<b>39</b>
<b>7</b>	<b>Instrument Data Description</b>	<b>40</b>
7.1	Bias frames . . . . .	40
7.2	Dark frames . . . . .	40
7.3	Slit flat-field frames . . . . .	40
7.4	Formatcheck frames . . . . .	41
7.5	Single fibre order definition frames . . . . .	41
7.6	Fibre Flat-Field (odd-even-all) frames . . . . .	41
7.7	Wavelength calibration frames . . . . .	42
7.8	Supported raw frames (keyword identifiers) . . . . .	43

<b>ESO</b>	<b>FLAMES-UVES Pipeline User Manual</b>	Doc:	VLT-MAN-ESO-19500-3016
		Issue:	Issue 5.0
		Date:	Date 2008-11-28
		Page:	6 of 110

<b>8</b>	<b>Static Calibration Data</b>	<b>50</b>
8.1	Line reference table . . . . .	50
<b>9</b>	<b>Data Reduction</b>	<b>51</b>
9.1	Data reduction overview . . . . .	51
9.2	Required input data . . . . .	51
9.3	Reduction cascade . . . . .	52
9.3.1	Data reduction peculiarities . . . . .	53
<b>10</b>	<b>Pipeline Recipes Interfaces</b>	<b>55</b>
10.1	uves_cal_mbias . . . . .	56
10.1.1	Input . . . . .	56
10.1.2	Output . . . . .	56
10.1.3	Quality control . . . . .	56
10.1.4	Parameters . . . . .	57
10.2	uves_cal_mdark . . . . .	57
10.2.1	Input . . . . .	57
10.2.2	Output . . . . .	57
10.2.3	Quality control . . . . .	57
10.2.4	Parameters . . . . .	58
10.3	flames_cal_mkmaster . . . . .	58
10.3.1	Input . . . . .	58
10.3.2	Output . . . . .	59
10.3.3	Quality control . . . . .	59
10.3.4	Parameters . . . . .	59
10.4	flames_cal_predict . . . . .	59
10.4.1	Input . . . . .	60
10.4.2	Output . . . . .	60
10.4.3	Quality control . . . . .	61
10.4.4	Parameters . . . . .	61
10.5	flames_cal_orderpos . . . . .	62

<b>ESO</b>	<b>FLAMES-UVES Pipeline User Manual</b>	Doc:	VLT-MAN-ESO-19500-3016
		Issue:	Issue 5.0
		Date:	Date 2008-11-28
		Page:	7 of 110

10.5.1	Input	62
10.5.2	Output	62
10.5.3	Quality control	63
10.5.4	Parameters	63
10.6	flames_cal_prep_sff_ofpos	64
10.6.1	Input	64
10.6.2	Output	65
10.6.3	Quality control	65
10.6.4	Parameters	65
10.7	flames_cal_wavecald	66
10.7.1	Input	66
10.7.2	Output	66
10.7.3	Quality control	68
10.7.4	Parameters	69
10.8	flames_obs_scired	71
10.8.1	Input	71
10.8.2	Output	72
10.8.3	Quality control	72
10.8.4	Quality control	72
10.8.5	Parameters	73
10.9	flames_obs_redchain	73
10.9.1	Input	73
10.9.2	Output	74
10.9.3	Quality control	74
10.9.4	Parameters	74
<b>11</b>	<b>Algorithms and recipe details</b>	<b>81</b>
11.1	Algorithms	81
11.1.1	Error model	81
11.1.2	Physical model offset setting	81
11.1.3	Fibre frames preparation	81

<b>ESO</b>	<b>FLAMES-UVES Pipeline User Manual</b>	Doc:	VLT-MAN-ESO-19500-3016
		Issue:	Issue 5.0
		Date:	Date 2008-11-28
		Page:	8 of 110

11.1.4	Slit frames preparation . . . . .	82
11.1.5	Lit fibres flag setting . . . . .	82
11.1.6	Inter-order background determination . . . . .	82
11.1.7	Fibre shift determination . . . . .	83
11.1.8	Generation of UVES related fibre information . . . . .	84
11.1.9	Frames preparation . . . . .	84
11.1.10	Bad pixel map preparation . . . . .	85
11.1.11	Fibre order tracing . . . . .	85
11.1.12	Iter order background determination . . . . .	86
11.1.13	Measure of flat part size of Slit frames . . . . .	86
11.1.14	Slit frames preparation . . . . .	86
11.1.15	Fiber frames preparation . . . . .	86
11.1.16	Fiber frames normalization . . . . .	87
11.1.17	Fiber frames extraction . . . . .	88
11.1.18	Fiber frames merging . . . . .	88
11.1.19	Interorder background subtraction . . . . .	89
11.1.20	Hough transform . . . . .	89
11.1.21	Order tracing . . . . .	90
11.1.22	Line Search . . . . .	90
11.1.23	Wavelength calibration first solution determination . . . . .	91
11.1.24	Line identification . . . . .	91
11.2	Recipes . . . . .	92
11.2.1	Master bias generation . . . . .	92
11.2.2	Master dark generation . . . . .	92
11.2.3	Spectral format guess . . . . .	92
11.2.4	Single fibre order table generation . . . . .	92
11.2.5	Master slit flat field generation . . . . .	92
11.2.6	Slit flats preparation, odd-even-all fibre flats preparation, normalization and fibre-order table determination . . . . .	93
11.2.7	Wavelength Calibration . . . . .	93
11.2.8	Science data reduction . . . . .	94



<b>ESO</b>	<b>FLAMES-UVES Pipeline User Manual</b>	Doc:	VLT-MAN-ESO-19500-3016
		Issue:	Issue 5.0
		Date:	Date 2008-11-28
		Page:	9 of 110

11.2.9 Full UVES-FIBRE data reduction . . . . .	95
<b>A Installation . . . . .</b>	<b>96</b>
A.1 Supported platforms . . . . .	96
A.2 Requirements . . . . .	96
A.3 Building the UVES pipeline . . . . .	97
<b>B FLAMES-UVES pipeline data reduction parameters . . . . .</b>	<b>99</b>
<b>C Full list of changes . . . . .</b>	<b>106</b>
<b>D Abbreviations and acronyms . . . . .</b>	<b>108</b>

<b>ESO</b>	FLAMES-UVES Pipeline User Manual	Doc:	VLT-MAN-ESO-19500-3016
		Issue:	Issue 5.0
		Date:	Date 2008-11-28
		Page:	10 of <a href="#">110</a>

<b>ESO</b>	<b>FLAMES-UVES Pipeline User Manual</b>	Doc:	VLT-MAN-ESO-19500-3016
		Issue:	Issue 5.0
		Date:	Date 2008-11-28
		Page:	11 of 110

# 1 Introduction

## 1.1 Purpose

The FLAMES-UVES pipeline is a subsystem of the *VLT Data Flow System* (DFS). Its target user is ESO *Data Flow Operations* (DFO) in the generation of master calibration data, in the reduction of scientific exposures, and in the data quality control. It should also serve as a quick look tool for *Paranal Science Operations* (PSO). Additionally, the FLAMES-UVES pipeline recipes are made public to the user community, to allow a more personalised processing of the data from the instrument.

This manual is a complete description of the data reduction recipes implemented by the CPL based FLAMES-UVES pipeline, reflecting the status of the FLAMES-UVES pipeline as of Nov 28, 2008 (version 4.3.0). Release 4.3.0 supports the reduction of UVES frames obtained when the UVES slit is fed by the fiber link to FLAMES.

## 1.2 Acknowledgements

The FLAMES-UVES pipeline has been initially (2001-2003) developed as a MIDAS package by Andrea Modigliani, as a wrapper of the FLAMES-UVES Data Reduction Software developed by Giacomo Mulas, Ignazio Porceddu and Francesco Damiani of the Ital-FLAMES consortium. Starting in winter 2006, to standardize implementation across all ESO pipelines, the FLAMES-UVES pipeline has been ported to CPL by Andrea Modigliani and Jonas Møller Larsen. We would like to thank for John Pritchard, from the Data Products Department <sup>2</sup>, and Domenic Naef from Paranal Science Operations for the feedback provided to improve this release.

## 1.3 Scope

This document describes the CPL based FLAMES-UVES pipeline used at ESO-Garching and ESO-Paranal for the purpose of data assessment and data quality control.

Updated versions of the present document may be found on [1]. For general information about the current instrument pipeline status we remind the user of [2]. Quality Control information are at [3].

Additional information on QFITS, the Common Pipeline Library (CPL) and EsoRex can be found respectively at [4], [5], [6]. The Gasgano front end is described in [7]. Description of the instrument are in [8,9,10]. The FLAMES-UVES instrument user manual is in [10]. The FLAMES-UVES calibration plan is in [11], while results of Science Verifications (SV) are at [12]. Additional information on the DFS and VLT data interfaces are in [13,14] and [15]. A clear and detailed description of the FLAMES-UVES data reduction software design is in [16]. First result obtained with the FLAMES-UVES data reduction software are described in [17]. The UVES user manual can be found at [18]. The UVES pipeline user manual can be found at [19]. A clear and compact description of the FLAMES-UVES pipeline is in [20].

---

<sup>1</sup>While the pipeline will generally cope with data acquired in non-standard modes ESO offers NO quality control over such data, and therefore it is responsibility of the user to verify the scientific efficacy of any calibration and/or science products produced by the pipeline for non-standard modes

<sup>2</sup>formerly Data Flow Operation Department

<b>ESO</b>	<b>FLAMES-UVES Pipeline User Manual</b>	Doc:	VLT-MAN-ESO-19500-3016
		Issue:	Issue 5.0
		Date:	Date 2008-11-28
		Page:	12 of 110

## 1.4 Reference documents

- [1] FLAMES-UVES Pipeline Users' Manual VLT-MAN-ESO-19500-3016  
[www.eso.org/pipelines](http://www.eso.org/pipelines)
- [2] Current pipeline status  
[www.eso.org/observing/dfo/quality/pipeline-status.html](http://www.eso.org/observing/dfo/quality/pipeline-status.html)
- [3] ESO-Data Flow Operation home page [www.eso.org/observing/dfo/quality/](http://www.eso.org/observing/dfo/quality/)
- [4] QFITS home page [www.eso.org/projects/aot/qfits/](http://www.eso.org/projects/aot/qfits/)
- [5] CPL home page [www.eso.org/cpl](http://www.eso.org/cpl)
- [6] EsoRex home page [www.eso.org/cpl/esorex.html](http://www.eso.org/cpl/esorex.html)
- [7] Gasgano User's Manual VLT-PRO-ESO-19000-1932
- [8] FLAMES home page [www.eso.org/instruments/flames](http://www.eso.org/instruments/flames)
- [9] *Installation and commissioning of FLAMES, the VLT Multifibre Facility* The Messenger, **110**, 1, 2002.
- [10] VLT Paranal Science Operations  
FLAMES User Manual VLT-MAN-ESO-13700-2994  
[www.eso.org/instruments/flames/doc](http://www.eso.org/instruments/flames/doc)
- [11] VLT FLAMES Calibration Plan VLT-PLA-ESO-13200-1567  
[www.eso.org/instruments/flames/doc](http://www.eso.org/instruments/flames/doc)
- [12] FLAMES SV home page  
[www.eso.org/sci/activities/vltsv/flamessv](http://www.eso.org/sci/activities/vltsv/flamessv)
- [13] VLT Data Flow System Specifications for Pipeline  
and Quality Control VLT-SPE-ESO-19600-1233
- [14] DFS Pipeline & Quality Control – User Manual VLT-MAN-ESO-19500-1619
- [15] ESO DICB – Data Interface Control Document GEN-SPE-ESO-00000-0794
- [16] FLAMES-UVES Data Reduction Software VLT-TRE-ITA-13750-0002  
Design and Architecture Report  
[www.eso.org/org/dmd/DFS-related-papers.html](http://www.eso.org/org/dmd/DFS-related-papers.html)
- [17] Automatic data reduction in support of the FLAMES-UVES VLT facility  
Proceedings of the SPIE, Volume 4844, pp. 310-320 (2002).  
[www.eso.org/org/dmd/DFS-related-papers.html](http://www.eso.org/org/dmd/DFS-related-papers.html)
- [18] UVES User Manual VLT-MAN-ESO-13200-1825  
[www.eso.org/instruments/uves/doc](http://www.eso.org/instruments/uves/doc)
- [19] UVES Pipeline User Manual VLT-MAN-ESO-19500-2964  
[www.eso.org/pipelines](http://www.eso.org/pipelines)
- [20] The FLAMES-UVES Pipeline The ESO Messenger, 118, 8.

<b>ESO</b>	<b>FLAMES-UVES Pipeline User Manual</b>	Doc:	VLT-MAN-ESO-19500-3016
		Issue:	Issue 5.0
		Date:	Date 2008-11-28
		Page:	13 of 110

## 2 Overview

In collaboration with instrument consortia, the Pipeline Systems Department (PSD) of the Software Development Division is implementing data reduction pipelines for the most commonly used VLT/VLTI instrument modes. These data reduction pipelines have the following three main purposes:

**Data quality control:** pipelines are used to produce the quantitative information necessary to monitor instrument performance.

**Master calibration product creation:** pipelines are used to produce master calibration products (*e.g.*, combined bias frames, super-flats, wavelength dispersion solutions).

**Science product creation:** using pipeline-generated master calibration products, science products are produced for the supported instrument modes (*e.g.*, combined ISAAC jitter stacks; bias-corrected, flat-fielded FORS images, wavelength-calibrated UVES spectra). The accuracy of the science products is limited by the quality of the available master calibration products and by the algorithmic implementation of the pipelines themselves. **In particular, adopted automatic reduction strategies may not be suitable or optimal for all scientific goals.**

Instrument pipelines consist of a set of data processing modules that can be called from the command line, from the automatic data management tools available on Paranal or from Gasgano.

ESO offers two front-end applications for launching pipeline recipes, *Gasgano* [14] and *EsoRex*, both included in the pipeline distribution (see Appendix A, page 96). These applications can also be downloaded separately from [www.eso.org/gasgano](http://www.eso.org/gasgano) and [www.eso.org/cpl/esorex.html](http://www.eso.org/cpl/esorex.html). An illustrated introduction to Gasgano is provided in Section 5.

The FLAMES facility and the different types of FLAMES-UVES raw frames and auxilliary data are described in Sections 4, 7, and 8.

A brief introduction to the usage of the available reduction recipes using Gasgano or EsoRex is presented in Section 5. In section 6 we advise the user about known data reduction problems.

An overview of the data reduction, the input data, and the recipes involved in the calibration cascade is provided in section 9.

More details on inputs, products, quality control measured quantities, and controlling parameters of each recipe is given in section 10.

More detailed descriptions of the data reduction algorithms used by the individual pipeline recipes can be found in Section 11.

In Appendix A the installation of the FLAMES-UVES pipeline recipes is described, in Appendix B the full list of relevant data reduction parameters is provided, in Appendix C the full list of changes is given, and in Appendix D a list of used abbreviations and acronyms is given.

<b>ESO</b>	<b>FLAMES-UVES Pipeline User Manual</b>	Doc:	VLT-MAN-ESO-19500-3016
		Issue:	Issue 5.0
		Date:	Date 2008-11-28
		Page:	14 of <a href="#">110</a>

### 3 What's new in pipeline release 4.3.0

The following major changes has been realised on the FLAMES-UVES pipeline:

- The pipeline has been ported to CPL.

For a more detailed list of changes see Appendix [C](#).

<b>ESO</b>	<b>FLAMES-UVES Pipeline User Manual</b>	Doc:	VLT-MAN-ESO-19500-3016
		Issue:	Issue 5.0
		Date:	Date 2008-11-28
		Page:	15 of 110

## 4 FLAMES Instrument Description

FLAMES is the multi-object, intermediate and high resolution spectrograph of the VLT. Mounted at UT2, FLAMES can access targets over a field of view 25 arcmin in diameter. FLAMES feeds two different spectrographs covering the whole visual spectral range: GIRAFFE and UVES. GIRAFFE allows the observation of up to 130 targets at the time or to do integral field spectroscopy, with intermediate resolution (either  $R \approx 25000$  or  $R \approx 10000$ ). UVES provides the maximum possible resolution ( $R=47000$ ) but can access only up to 8 objects at the time. The instrument has been made available to the community and started operations in Paranal on April 1<sup>st</sup>, 2000.

In this chapter a brief description of the FLAMES-UVES fiber link and the UVES instrument is given. A more complete documentation can be found in the FLAMES and UVES User Manuals [10,18].



Figure 4.0.1: A photo of FLAMES mounted at the Nasmyth A focus of Kuyen (VLT-UT2).

### 4.1 Instrument overview

FLAMES consists of three main components:

- A Fibre Positioner (OzPoz) hosting two plates: while one plate is observing the other positions the fibres for the subsequent observations, therefore limiting the dead time between one observation and the next to less than 15 minutes, including the telescope preset and the acquisition of the next field.

<b>ESO</b>	<b>FLAMES-UVES Pipeline User Manual</b>	Doc:	VLT-MAN-ESO-19500-3016
		Issue:	Issue 5.0
		Date:	Date 2008-11-28
		Page:	16 of 110

- A medium-high resolution optical spectrograph, GIRAFFE, with three types of fibre-fed systems : MEDUSA (132 individually deployable 1.2 arcsec entrance aperture fibres), IFU (15 deployable 2x3 arcmin aperture IFUs), ARGUS (an 11.5 or 6.6x4.2 arcsec Integral Field Unit).
- A link to the UVES spectrograph (Red Arm) via 8 single fibres of 1 arcsec entrance aperture.

Special observing software (FLAMES OS) coordinates the operation of the different subsystems, also allowing simultaneous acquisition of UVES and GIRAFFE observations. For combined observations, the exposure times for UVES and GIRAFFE do not need to be the same. Note that it is not possible to observe simultaneously in two GIRAFFE modes, or to observe the same target simultaneously with the two spectrographs.

UVES is the high resolution spectrograph at UT2 of the VLT (see Section 6.4). It was designed to work in long slit mode but it has been possible to add a fibre mode (6 or 8 fibres, depending on setup and/or mode) fed by the FLAMES positioner to its Red Arm only. Only three standard UVES Red setups are offered, with central wavelength of 520, 580 and 860 nm respectively (see the manual for details). The standard readout mode of FLAMES-UVES is 225 kHz (unbinned) which ensures low readout noise. As of P76 a high-speed readout mode (625 kHz, unbinned, low gain) with increased readout noise but less overheads is offered in visitor mode only. No pipeline support is available in this mode. With an aperture on the sky of 1 arcsec, the fibres project onto a resolution spot of diameter 5 UVES pixels giving a resolving power of 47000. For faint objects and depending on the spectral region, one or more fibres can be devoted to recording the sky contribution. In addition, for the 580 nm setup only, a separate calibration fibre is available to acquire simultaneous ThAr calibration spectra. This allows very accurate radial velocity determinations. In this configuration, 7 fibres remain available for targets on sky.

Spectrograph	Mode	N. of Objects	Aperture (arcsecs)	Resolving power	Spectral Band [nm]	Central wavelength	Observing Mode
UVES	Red arm	8	1.0	47000	200	520/580/860	OzPoz
UVES	Red arm	7+1(calibration)	1.0	47000	200	580	SimCal

Table 4.1: This table summarize main features of UVES fibres observations



<b>ESO</b>	<b>FLAMES-UVES Pipeline User Manual</b>	Doc:	VLT-MAN-ESO-19500-3016
		Issue:	Issue 5.0
		Date:	Date 2008-11-28
		Page:	17 of 110

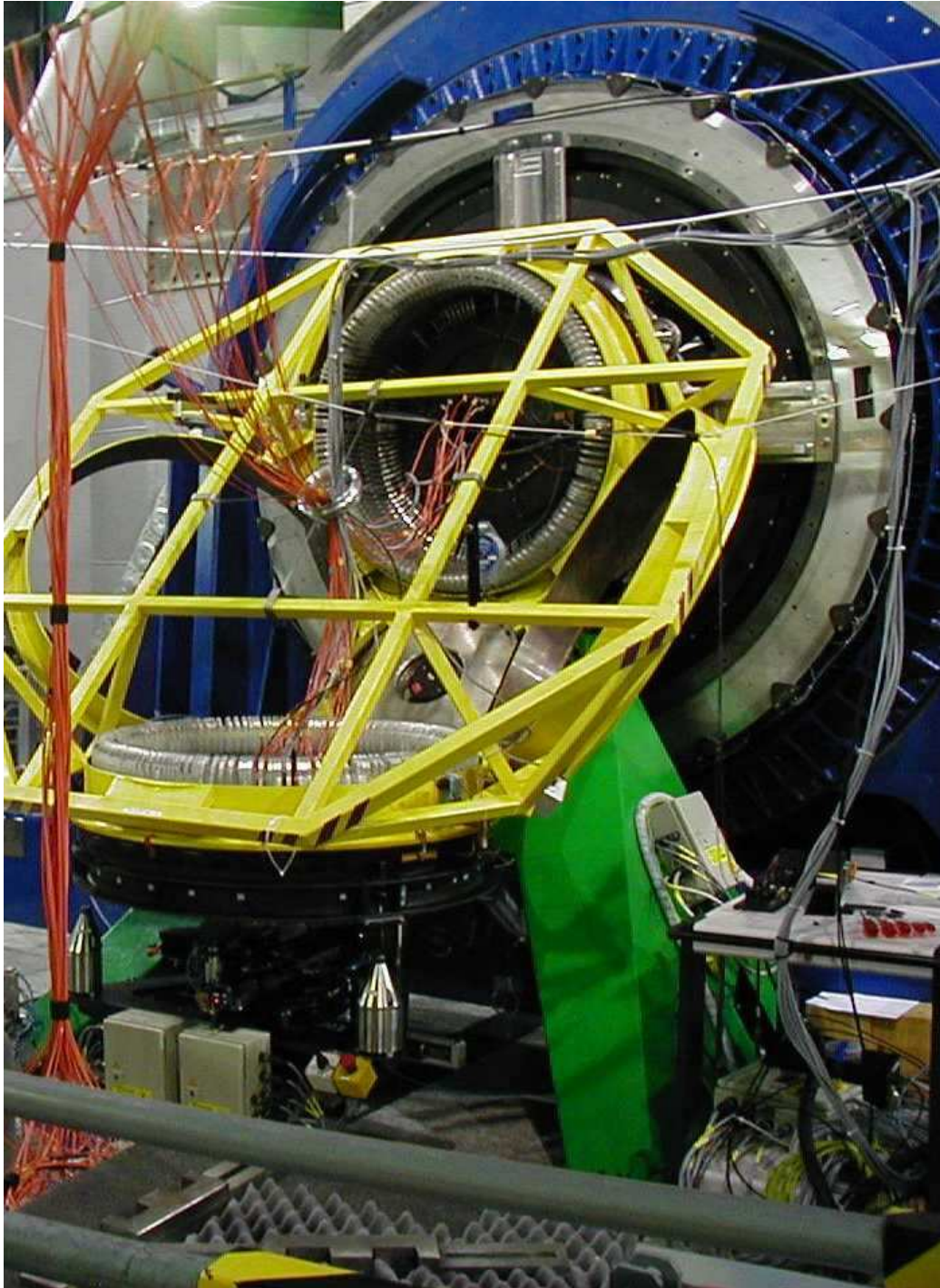


Figure 4.1.1: A photo of the OzPoz fibre positioner. While one plate is observing, the other one is positioning the fibres for the subsequent observations.

<b>ESO</b>	<b>FLAMES-UVES Pipeline User Manual</b>	Doc:	VLT-MAN-ESO-19500-3016
		Issue:	Issue 5.0
		Date:	Date 2008-11-28
		Page:	18 of 110

## 5 Quick start

This section describes the most immediate usage of the FLAMES-UVES pipeline recipes.

### 5.1 FLAMES-UVES pipeline recipes

The current FLAMES-UVES pipeline is based on a set of 7 stand-alone recipes plus one utility recip (and uses 2 recipes from the UVES pipeline to combine biases and darks) involved in the data reduction cascade:

**uves\_cal\_mbias** creates a master bias frame.

**uves\_cal\_mdark** creates a master dark frame.

**flames\_cal\_predict** implements the UVES physical model (and applies a shift to reflect the offset between the calibration fibre and the order centre).

**flames\_cal\_orderpos** defines the calibration fibre order positions.

**flames\_cal\_mkmaster** creates master slit flat-field frames.

**flames\_cal\_prep\_sff\_ofpos** process master slit flats and odd-even-all fibre flats to determine the a table tracing all fibres and creates additional frames, needed to extract the science data.

**flames\_cal\_wavecal** performs the wavelength calibration.

**flames\_obs\_scired** reduces a science frame.

**flames\_utl\_unpack** unpack products relative to each fibre from the products usually delivered by the Quality Control Group of ESO, Garching.

Alternatively, the full data reduction chain can be run using the recipe **flames\_obs\_redchain**.

### 5.2 An introduction to Gasgano and EsoRex

Before being able to call pipeline recipes on a set of data, the data must be properly classified, and associated with the appropriate calibrations. *Data Classification* consists of tasks such as: "What kind of data am I?", *e.g.*, BIAS, "to which group do I belong?", *e.g.*, to a particular Observation Block or template. *Data Association* is the process of selecting appropriate calibration data for the reduction of a set of raw science frames. Typically, a set of frames can be associated if they share a number of properties, such as instrument and detector configuration. As all the required information is stored in the FITS headers, data association is based on a set of keywords (called "association keywords") which, generally, are specific to each type of calibration.

The process of data classification and association is known as data organisation.

Service Mode users also receive, as part of their Service Mode Data Package, files called Association Blocks (or ABs) which, amongst other things, specify the "optimum" (normally in the sense of closest in time certified by QC calibration) associations as determined by the Data Product and Quality Control Group of ESO, Garching.

<b>ESO</b>	<b>FLAMES-UVES Pipeline User Manual</b>	Doc:	VLT-MAN-ESO-19500-3016
		Issue:	Issue 5.0
		Date:	Date 2008-11-28
		Page:	19 of 110

An instrument pipeline consists of a set of data processing modules that can be called from different host applications, either from the command line with *EsoRex*, from the automatic data management tools available at Paranal, or from the graphical *Gasgano* tool.

*Gasgano* is a data management tool that simplifies the data organisation process, offering automatic data classification and making the data association easier (*even if automatic association of frames is not yet provided*). *Gasgano* determines the classification of a file by applying an instrument specific rule, while users must provide this information to the recipes when they are executed manually using *EsoRex* from the command line. In addition, *Gasgano* allows the user to execute directly the pipeline recipes on a set of selected files.

### 5.2.1 Using Gasgano

To get familiar with the FLAMES-UVES pipeline recipes and their usage, it is advisable to begin with *Gasgano*, because it provides a complete graphic interface for data browsing, classification and association, and offers several other utilities such as easy access to recipe documentation and preferred data display tools.

*Gasgano* can be started from the system prompt in the following way:

```
gasgano $HOME/gasgano/config/UVES.prefs &
```

where we have passed as first optional argument explicitly the UVES preference file which defines proper defaults for UVES data reduction (for both ECHELLE and FIBER modes). The user may like to alias the *gasgano* command to the previous command. The *Gasgano* main window will appear. In Figure 5.2.1 (next page), a view on a set of FLAMES-UVES data is shown as an example. *Gasgano* can be pointed to the directories where the data to be handled are located using the navigation panels accessible via the *Add/Remove Files* entry of the *File* menu (shown on the upper left of the figure).

The data are hierarchically organised as preferred by the user. After each file name are shown the classification and the values of the following FITS keywords (we omit the prefix HIERARCH.ESO):

Keyword name	Purpose
CLASSIFICATION	Data classification
OBS.TARG.NAME	Observation Block target name
EXPTIME	Exposure time
DATE	Observing date
INS.GRAT2.WLEN	Instrument setting central wavelength (for red arm)
DET.READ.SPEED	Readout speed
INS.OBSPLATE	Observation plate
OCS.SIMCAL	Simultaneous calibration

Additionally *Gasgano* uses other keywords needed to properly classify UVES data taken using the slit mode.

The **CLASSIFICATION** field provides either the value of the **PRO.CATG**, for pipeline products; or a user defined file classification, if provided, defined in the classification rule file, which can be accessed by *Gasgano* from the Tools → Classification rules... tab; or the default value “UNDEFINED”. File classification rules are selection rules which assign to a FITS file a classification based on the value of a few FITS keywords, usually the

<b>ESO</b>	<b>FLAMES-UVES Pipeline User Manual</b>	Doc:	VLT-MAN-ESO-19500-3016
		Issue:	Issue 5.0
		Date:	Date 2008-11-28
		Page:	20 of 110

DPR.TYPE, DPR.TECH, DPR.CATG values, which respectively define the file data type, acquisition technique and category.

Keyword name	Purpose
DET.CHIPS	# of chips in detector array
INS.MODE	Instrument mode used
INS.GRAT1.NAME	Instrument grating name (for blue arm)
INS.GRAT1.WLEN	Instrument setting central wavelength (for blue arm)
INS.SLIT2.WID	Instrument slit width (for blue arm)
INS.GRAT2.NAME	Instrument grating name (for red arm)
INS.SLIT3.WID	Instrument slit width (for red arm)
DET.WIN1.BINX	Binning factor along X
DET.WIN1.BINY	Binning factor along Y

Additional relevant keywords for FLAMES-UVES data are INS.GRATi.WLEN, INS.OBSPLATE, OCS.SIMCAL. In order to properly classify data taken in slit mode it is useful also to display the keyword values of DET.CHIPS, INS.SLITi.NAME, DET.WIN1.BINX/Y, INS.MODE, INS.GRATi.NAME. Those relevant keywords are indicated by Gasgano either in the file section, or by selecting each file, in the section which shows the FITS file header content.

More information about a single frame can be obtained by clicking on its name: the corresponding FITS file header will be displayed on the bottom panel, where specific keywords can be filtered and searched. Images and tables may be easily displayed using the viewers specified in the appropriate *Preferences* fields. Such a field also allows to set the file filter, which should point to the \$HOME/gasgano/config/UVES.rul. This rule file provides simple filtering rules to select FLAMES-UVES data corresponding to a given standard data reduction setting.

Frames can be selected from the main window for by the appropriate recipe: in Figure 5.2.2, the Wavelength calibration frame (CLASSIFICATION=FIB\_ARC\_LAMP\_RED), previously produced master bias and master flat frames, together with line and order tables, a table with the reference standard star spectra and one with the atmospheric dispersion are all selected and sent to the *flames\_cal\_wavecal* recipe. This will open a *Gasgano* recipe execution window (see Figure 5.2.3), having all the specified files listed in its *Input Frames* panel.

Help for the recipe is available from the *Help* menu. Before launching the recipe, its configuration may be modified on the *Parameters* panel (on top). The window contents might be saved for later use by selecting the *Save Current Settings* entry from the *File* menu, as shown in figure.

At this point the recipe can be launched by pressing the *Execute* button. Messages from the running recipe will appear on the *Log Messages* panel at the bottom, and in case of successful completion the products will be listed on the *Output Frames* panel, where they can be easily viewed and located back on the Gasgano main window. To produce useful plots the user needs to set to 'gnuplot -persist' the **plotter** recipe parameter value (and have a valid installation of gnuplot package, and the gnuplot command available in the PATH). This basic plotting functionality is supported by EsoRex but not by gasgano.

Please refer to the *Gasgano User's Manual* [7] for a more complete description of the *Gasgano* interface.



ESO	FLAMES-UVES Pipeline User Manual	Doc:	VLT-MAN-ESO-19500-3016
		Issue:	Issue 5.0
		Date:	Date 2008-11-28
		Page:	21 of 110

GASGANO Version: 2.2.7 esomidas / Linux

File Selected files Tools Help

Default grouping collapse Find entry: find 1

Displaying 89 files Unfiltered.

File	CLASSIFICATION	INS.SETUP.ID	INS.OPT11.NAME	TPL.EXPNO	TPL.NEXP	OBS.TARG.NAME	DPR.TYPE
073.D-0211(A) UVES UNKNOWN							
155910 N0288HR11b							
UVES.2004-07-26T05:24:13.263.fits	FIB_SCL_RED	OUT		2	2	NGC288	OBJECT_O2Poz
60.A-9022(A) UVES UNKNOWN							
200006102 red_225kHz_1x1_low_cal							
mbias_l.fits	MASTER_BIAS_REDL			1	5	Calibration-ob	BIAS
mbias_u.fits	MASTER_BIAS_REDU			1	5	Calibration-ob	BIAS
mbsre1x1.fits	MASTER_BIAS_REDL			1	5	Calibration-ob	BIAS
mbsrm1x1.fits	MASTER_BIAS_REDU			1	5	Calibration-ob	BIAS
60.A-9022(B) UVES							
200117229 Calibration							
UVES.2003-02-02T11:08:23.305.fits	SFLAT_RED	OUT		1	9		LAMP,SFLAT
UVES.2003-02-02T11:09:47.871.fits	SFLAT_RED	OUT		2	9		LAMP,SFLAT
UVES.2003-02-02T11:11:12.282.fits	SFLAT_RED	OUT		3	9		LAMP,SFLAT
UVES.2003-02-02T11:14:04.807.fits	SFLAT_RED	OUT		5	9		LAMP,SFLAT
UVES.2003-02-02T11:15:28.876.fits	SFLAT_RED	OUT		6	9		LAMP,SFLAT
UVES.2003-02-02T11:17:03.161.fits	SFLAT_RED	OUT		7	9		LAMP,SFLAT
UVES.2003-02-02T11:18:27.022.fits	SFLAT_RED	OUT		8	9		LAMP,SFLAT
UVES.2003-02-02T11:19:52.398.fits	SFLAT_RED	OUT		9	9		LAMP,SFLAT
UVES.2004-07-25T10:58:31.580.fits	FIB_FF_ODD_RED	OUT		1	3		LAMP,FLAT,O...
UVES.2004-07-25T11:06:43.602.fits	FIB_FF_EVEN_RED	OUT		2	3		LAMP,FLAT,EV...
UVES.2004-07-25T11:15:18.016.fits	FIB_FF_ALL_RED	OUT		3	3		LAMP,FLAT,AL...
UVES.2004-07-25T11:24:23.142.fits	FIB_ARC_LAMP_RED	OUT		1	1		LAMP,WAVE,O...
UVES.2004-07-25T13:08:16.774.fits	FIB_ORDEF_RED	OUT		1	1		LAMP,ORDER...
UVES.2004-07-25T13:09:39.070.fits	FIB_ARC_LAMP_FORM_RED	OUT		1	1		LAMP,FMTCH...
b_set_1_mf580_1x1_s10_l.fits	MASTER_SFLAT_REDL	OUT		1	9		LAMP,SFLAT
b_set_1_mf580_1x1_s10_u.fits	MASTER_SFLAT_REDU	OUT		1	9		LAMP,SFLAT
b_set_2_mf580_1x1_s10_l.fits	MASTER_SFLAT_REDL	OUT		4	9		LAMP,SFLAT
b_set_2_mf580_1x1_s10_u.fits	MASTER_SFLAT_REDU	OUT		4	9		LAMP,SFLAT
b_set_3_mf580_1x1_s10_l.fits	MASTER_SFLAT_REDL	OUT		7	9		LAMP,SFLAT
b_set_3_mf580_1x1_s10_u.fits	MASTER_SFLAT_REDU	OUT		7	9		LAMP,SFLAT

/media/disk/flames/flames\_uvcs\_demo/raw/mos/UVES.2004-07-25T11:24:23.142.fits FLAMES\_UVES\_WAVE207\_0001.fits FIB\_ARC\_LAMP\_RED

Extension: HEADER Find in header: find Load Filter Filter Auto Display

Keyword	Value
SIMPLE	T
BITPIX	16
NAXIS	0
EXTEND	T
ORIGIN	ESO
DATE	2004-07-25T11:25:15.713
TELESCOP	ESO-VLT-U2
INSTRUME	UVES
OBJECT	LAMP,WAVE,O2Poz
EXPTIME	51.6637
MJD-OBS	53211.47526785
DATE-OBS	2004-07-25T11:24:23.142
FL_COI	FLAMES Operation Terms

Figure 5.2.1: The Gasgano main window.

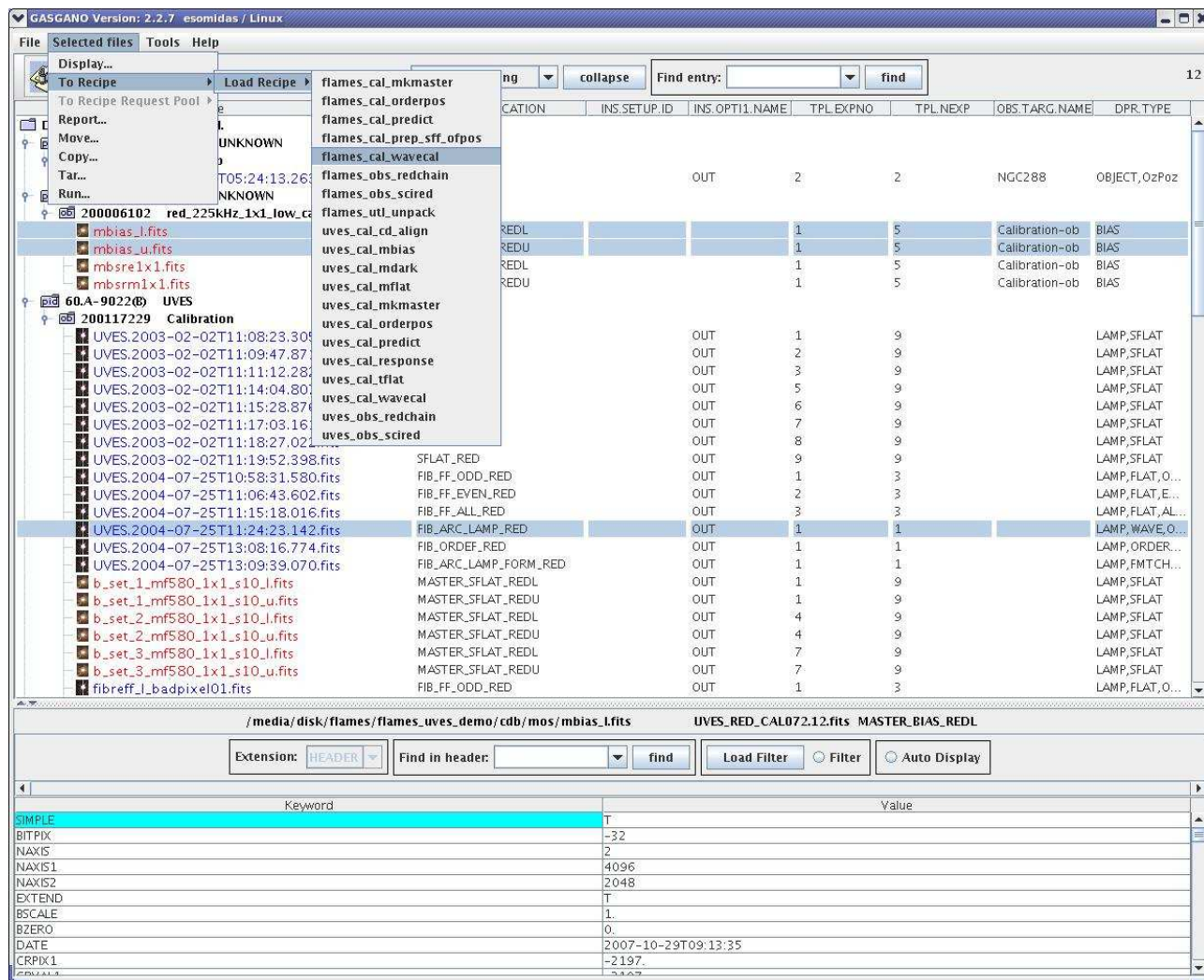


Figure 5.2.2: *Selecting files to be processed by a FLAMES-UVES pipeline recipe.*

ESO	FLAMES-UVES Pipeline User Manual	Doc:	VLT-MAN-ESO-19500-3016
		Issue:	Issue 5.0
		Date:	Date 2008-11-28
		Page:	23 of 110

Flames\_cal\_wavecal v40203

File Help

Current Queued Executing

**Parameters**

Name	Value	Default	Range
uves.debug	<input type="checkbox"/>	<input type="checkbox"/>	
uves.plotter	no	no	
uves.process.chip	both	both	
flames_cal_wavecal.nwindows	1	1	1..2147483647
flames_cal_wavecal.length	7.0	7.0	-2.0..1.797693134...
flames_cal_wavecal.offset	0.0	0.0	0.0..1.7976931348...

**Input Frames**

Include	Filename	Classification	Locate	Display
<input checked="" type="checkbox"/>	mbias_l.fits	MASTER_BIAS_REDL	Locate	Display
<input checked="" type="checkbox"/>	mbias_u.fits	MASTER_BIAS_REDU	Locate	Display
<input checked="" type="checkbox"/>	UVES.2004-07-25T11:24:23.142.fits	FIB_ARC_LAMP_RED	Locate	Display
<input checked="" type="checkbox"/>	lineguesstable_redl.fits	FIB_LIN_GUE_REDL	Locate	Display
<input checked="" type="checkbox"/>	lineguesstable_redu.fits	FIB_LIN_GUE_REDU	Locate	Display
<input checked="" type="checkbox"/>	order_def_redl.fits	FIB_ORDEF_REDL	Locate	Display

**Product Naming**

Product Root Directory: /diska/home/esomidas/tmp  Naming Scheme: Numeric

**Output Frames**

Filename	Classification	Locate	Display
linetable_redl_0000.fits	FIB_LINE_TABLE_REDL	Locate	Display
linetable_redu_0000.fits	FIB_LINE_TABLE_REDU	Locate	Display

**Log Messages**

[ INFO ] Wavecal: Line table for trace 7, window #1 saved to extensions 19-21 of 'linetable\_redu.fits'

[ INFO ] Wavecal: Line table for trace 8, window #1 saved to extensions 22-24 of 'linetable\_redu.fits'

[WARNING] Wavecal: Recipe produced 22 warnings (excluding this one)

Product frames:

/diska/home/esomidas/tmp/linetable\_redl\_0000.fits

/diska/home/esomidas/tmp/linetable\_redu\_0000.fits

Completion status: SUCCESS

Figure 5.2.3: The Gasgano recipe execution window.

ESO	FLAMES-UVES Pipeline User Manual	Doc:	VLT-MAN-ESO-19500-3016
		Issue:	Issue 5.0
		Date:	Date 2008-11-28
		Page:	24 of 110

### 5.2.2 Using EsoRex

*EsoRex* is a command line utility for running pipeline recipes. It can be used in data reduction scripts for the automation of processing tasks. Unlike when using *Gasgano*, the user must classify and associate the data using the information contained in the FITS header keywords (see Section 7, page 40)<sup>3</sup>. The user should also take care of defining the input set-of-frames and the appropriate configuration parameters for each recipe run:

**The set-of-frames:** Each pipeline recipe is run on a set of input FITS data files. When using *EsoRex* the filenames must be listed together with their DO category<sup>4</sup> in an ASCII file, the *set-of-frames* (SOF), that is required when launching a recipe.<sup>5</sup> DO categories for the supported FLAMES-UVES input raw frames are indicated in section 7.8.

```

/path_raw/UVES.2004-08-14T10:20:56.497.fits FIB_ARC_LAMP_RED
/path_pro/ordertable_redl.fits FIB_ORDEF_TABLE_REDL
/path_pro/lineguess_redl.fits FIB_LIN_GUE_REDL
/path_pro/ordertable_redu.fits FIB_ORDEF_TABLE_REDU
/path_pro/lineguess_redu.fits FIB_LIN_GUE_REDU
/path_pro/thargood_3.fits LINE_REFER_TABLE
/path_pro/masterbias_redl.fits MASTER_BIAS_REDL (optional)
/path_pro/masterbias_redu.fits MASTER_BIAS_REDU (optional)

```

It contains for each input frame the full path file name and its DO category. The pipeline recipe will access the listed files when required by the reduction algorithm.

Note that the FLAMES-UVES pipeline recipes do not verify in any way the correctness of the classification tags specified by the user in the SOF. In the above example, the recipe *flames\_cal\_wavecal* will treat the frame `/path_raw/UVES.2004-08-14T10:20:56.497.fits` as an `FIB_ARC_LAMP_RED`, the frame `/path_pro/ordertable_redl.fits` as a `FIB_ORDEF_TABLE_REDL`, etc., even when they do not contain this type of data. The recipe will also assume that all frames are associated correctly, *i.e.*, that they all come from the same arm, dichroic and bin setting.

The reason for this lack of control is that the FLAMES-UVES recipes are just the DRS component of the complete pipeline running on Paranal, where the task of data classification and association is carried out by separate applications. Moreover, using *Gasgano* as an interface to the pipeline recipes will always ensure a correct classification of all the data frames, assigning the appropriate DO category to each one of them (see Section 5.2.1, page 19). Also this lack of control allows the user to reduce e.g. an arc lamp frame pretending it is a science frame.

A recipe handling an incorrect SOF may stop or display unclear error messages at best. In the worst cases, the recipe would apparently run without any problem, producing results that may look reasonable, but are actually flawed.

<sup>3</sup>Service Mode users could also use the association information contained in the Association Blocks (or ABs) included in their Service Mode Data Packages

<sup>4</sup>The indicated *DO category* is a label assigned to any data type after it has been classified, which is then used to identify the frames listed in the *set-of-frames*

<sup>5</sup>The set-of-frames corresponds to the *Input Frames* panel of the *Gasgano* recipe execution window (see Figure 5.2.3, page 23).



<b>ESO</b>	<b>FLAMES-UVES Pipeline User Manual</b>	Doc:	VLT-MAN-ESO-19500-3016
		Issue:	Issue 5.0
		Date:	Date 2008-11-28
		Page:	25 of 110

**EsoRex syntax:** The basic syntax to use EsoRex is the following:

**esorex [esorex\_options] recipe\_name [recipe\_options] set\_of\_frames**

To get more information on how to customise EsoRex (see also [6]) run the command:

**esorex - -help**

To generate a configuration file `esorex.rc` in the directory `$HOME/.esorex` run the command:

**esorex - -create-config**

A list of all available recipes, each with a one-line description, can be obtained using the command:

**esorex - -recipes**

EsoRex searches for recipes in the directory specified by the option

**esorex - -recipe-dir=installation\_directory**

All recipe parameters (aliases) and their default values can be displayed by the command

**esorex - -params recipe\_name**

To get a brief description of each parameter meaning execute the command:

**esorex - -help recipe\_name**

To get more details about the given recipe give the command at the shell prompt:

**esorex - -man-page recipe\_name**

**Recipe configuration:** Each pipeline recipe may be assigned an *EsoRex* configuration file, containing the default values of the parameters related to that recipe.<sup>6</sup> The configuration files are normally generated in the directory `$HOME/.esorex`, and have the same name as the recipe to which they are related, with the filename extension `.rc`. For instance, the recipe *flames\_cal\_wavecal* has its *EsoRex* generated configuration file named `flames_cal_wavecal.rc`, and is generated with the command:

**esorex - -create-config flames\_cal\_wavecal**

The definition of one parameter of a recipe may look like this:

```
# --tolerance
# Tolerance of fit. If positive, 'tolerance' is in pixel units. If negative,
# abs('tolerance') is in wavelength units. Lines with residuals worse than
# the tolerance are excluded from the final fit. Unlike in previous versions,
# this parameter is not corrected for CCD binning.
flames_cal_wavecal.calibrate.tolerance=0.6
```

In this example, the parameter `flames_cal_wavecal.calibrate.tolerance` is set to the value `0.6`. In the configuration file generated by *EsoRex*, one or more comment lines are added containing information about the possible values of the parameter, and an alias that could be used as a command line option.

The command

**esorex - -create-config recipe\_name**

---

<sup>6</sup>The *EsoRex* recipe configuration file corresponds to the *Parameters* panel of the *Gasgano* recipe execution window (see Figure 5.2.3, page 23).

<b>ESO</b>	<b>FLAMES-UVES Pipeline User Manual</b>	Doc:	VLT-MAN-ESO-19500-3016
		Issue:	Issue 5.0
		Date:	Date 2008-11-28
		Page:	26 of 110

generates a default configuration file **recipe\_name.rc** in the directory **\$HOME/.esorex**<sup>7</sup>.

A recipe configuration file different from the default one can be specified on the command line:

**esorex - -recipe-config=my\_alternative\_recipe\_config**

Recipe parameters are provided in section 10 and their role is described in Section 11.

More than one configuration file may be maintained for the same recipe but, in order to be used, a configuration file not located under `$HOME/.esorex`, or having a name different from the recipe name, should be explicitly specified when launching a recipe.

**Recipe execution:** A recipe can be run by giving its name to *EsoRex*, together with the name of a set-of-frames. For instance, the following command line would be used to run the recipe *flames\_cal\_wavecal* for processing the files specified in the set-of-frames *flames\_cal\_wavecal.sof*:

**esorex flames\_cal\_wavecal flames\_cal\_wavecal.sof**

The recipe parameters can be modified either by editing directly the used configuration file, or by specifying new parameter values on the command line using the command line options defined for this purpose. Such command line options should be inserted after the recipe name and before the SOF name, and they will supersede the system defaults and/or the configuration file settings. For instance, to set the *flames\_cal\_wavecal* recipe *tolerance* parameter to 0.07, the following should be typed:

**esorex flames\_cal\_wavecal - -tolerance=0.07 flames\_cal\_wavecal.sof**

Every recipe provides a parameter *debug* which – when enabled – causes intermediate results to be saved to the local directory. This allows more detailed inspection of the recipe processing

**esorex flames\_cal\_wavecal - -debug flames\_cal\_wavecal.sof**

Basic plotting functionality can be enabled on systems which have the gnuplot tool:

**esorex flames\_cal\_wavecal - -plotter='gnuplot -persist' flames\_cal\_wavecal.sof**

For more advanced visualisation a dedicated FITS viewer should be used with the FITS files produced by the recipe.

For more information on *EsoRex*, see [www.eso.org/cpl/esorex.html](http://www.eso.org/cpl/esorex.html).

### 5.3 Example of data reduction using EsoRex

Here we provide an example of data reduction for data obtained with the FLAMES-UVES.

In Table 5.3 we list the category Data Products (DPR) keywords for the currently supported data.

The plate No. is identified through the keyword ESO.INS.OBSPLATE. Plate No. 1 corresponds to a keyword value equal to 1. Plate No. 2 corresponds to a keyword value equal to 2.

The keyword ESO.OCS.SIMCAL will be used to associate the proper master calibration frame to the input raw frame separating the cases in which the simultaneous calibration fibre is involved (value equal to 1) or not (value equal to 0).

<sup>7</sup>If a number of recipe parameters are specified on the command line, the given values will be used in the created configuration file.

<b>ESO</b>	<b>FLAMES-UVES Pipeline User Manual</b>	Doc:	VLT-MAN-ESO-19500-3016
		Issue:	Issue 5.0
		Date:	Date 2008-11-28
		Page:	27 of 110

Frame Type	TYPE	CATG	TECH	Note
Science frames				
Science in OzPoz mode	OBJECT,OzPoz	SCIENCE	MOS	8 OzPoz fibres
Science in SimCal mode	OBJECT,SimCal	SCIENCE	MOS	7 OzPoz, 1 SimCal
Calibration fibre fed frames				
Bias	BIAS	CALIB	IMAGE	Same as UVES
Dark	DARK	CALIB	IMAGE	Same as UVES
Nasmyth Screen FF	LAMP,FLAT,NASMYTH	CALIB	MOS	Nasmyth Screen
Wavelength calibration	LAMP,WAVE,OzPoz	CALIB	MOS	8 OzPoz fed by ThAr
Wavelength calibration	LAMP,WAVE,SimCal	CALIB	MOS	7 OzPoz, 1 SimCal fed by ThAr
Odd FF	LAMP,FLAT,ODD,OzPoz	CALIB	MOS	Odd Fiber FF, 4 OzPoz
Even FF	LAMP,FLAT,EVEN,OzPoz	CALIB	MOS	Even Fiber FF, 4 OzPoz
All FF	LAMP,FLAT,ALL,OzPoz	CALIB	MOS	All Fiber FF, 8 OzPoz
Odd FF	LAMP,FLAT,ODD,SimCal	CALIB	MOS	Odd fibre FF, 4 OzPoz
Even FF	LAMP,FLAT,EVEN,SimCal	CALIB	MOS	Even Fiber FF, 4 OzPoz
All FF	LAMP,FLAT,ALL,SimCal	CALIB	MOS	All Fiber FF, 7 OzPoz
Formatcheck	LAMP,FMTCHK,SimCal	CALIB	MOS	1 SimCal fed by ThAr
Order Trace	LAMP,ORDERDEF,SimCal	CALIB	MOS	1 SimCal fed by Flat
Special frames in echelle mode				
Slit Flat Fields	LAMP,SFLAT	CALIB	ECHELLE	3 sets of 3 slit FF each each set taken at diff Y pos

Table 5.3: Required DPR keywords

When an accuracy better than 2% is required for the absolute fluxes, one should take flat fields using the Nasmyth screen and “attached templates”.

In Table 5.4 we list the DO\_CLASSIFICATION values for the currently supported data.

<b>ESO</b>	<b>FLAMES-UVES Pipeline User Manual</b>	Doc:	VLT-MAN-ESO-19500-3016
		Issue:	Issue 5.0
		Date:	Date 2008-11-28
		Page:	28 of 110

Frame Type	DO_CLASSIFICATION
Raw Frames	
ThAr line reference table	LINE_REFER_TABLE
Formatcheck frame	FMTCHK_CHIP
Reference Formatcheck frame	ARC_LAMP_FORM_CHIP
Single Order trace frame	FIB_ORDEF_CHIP
Bias	BIAS_CHIP
Master Bias	MASTER_CHIP
Dark	DARK_CHIP
Slit Flat Field	SFLAT_CHIP
Nasmyth Screen Flat Field	FIB_FF_NASMYTH_CHIP
Wavelength Calibration frame	FIB_ARC_LAMP_CHIP
Odd Flat Field	FIB_FF_ODD_CHIP
Even Flat Field	FIB_FF_EVEN_CHIP
All Flat Field	FIB_FF_ALL_CHIP
Science OzPoz	FIB_SCI_CHIP
Science SimCal	FIB_SCI_SIM_CHIP
Science COMBINED	FIB_SCI_COM_CHIP
Data Products	
MASTER Slit Flat Field	MASTER_SFLAT_CHIP
DRS Setup Table	FIB_DRS_CHIP
Guess Order Table	FIB_ORD_GUE_CHIP
Guess Line Table	FIB_LIN_GUE_CHIP
Background Table	FIB_BKG_CHIP
Order-Fibre Table	FIB_ORDEF_TABLE_CHIP
Line-Fibre Table	FIB_LINE_TABLE_CHIP
Fibre data* frames	FIB_FF_D0n_CHIP
Fibre common frames	FIB_FF_COM_CHIP
Fibre norm frames	FIB_FF_NOR_CHIP
Fibre nsigma frames	FIB_FF_NSIGMA_CHIP
Fibre sigma frames	FIB_FF_S0n_CHIP
Fibre badpixel frames	FIB_FF_BPn_CHIP
Slit FF data frames	SLIT_FF_D0n_CHIP
Slit FF common frames	SLIT_FF_COM_CHIP
Slit FF norm frames	SLIT_FF_NOR_CHIP
Slit FF sigma frames	SLIT_FF_S0n_CHIP
Slit FF badpixel frames	SLIT_FF_BPn_CHIP
Slit FF bound frames	SLIT_FF_BNDn_CHIP

Table 5.4: List of frame types and corresponding tag values. \*Note: CHIP= REDL or REDU; n=1,2,...nflats where in case of slit flat fields NFLATS is a FITS keyword indicating the (automatically determined) minimum number of equalized flat fields which are necessary to cover the total Y span, while in case of the fibre flat fields it indicates the number of basic (uncontaminated) fibre flat fields, typically (but not necessarily) 2 (odd/even).

<b>ESO</b>	<b>FLAMES-UVES Pipeline User Manual</b>	Doc:	VLT-MAN-ESO-19500-3016
		Issue:	Issue 5.0
		Date:	Date 2008-11-28
		Page:	29 of 110

The simplest and least interactive way to reduce the data is to create a SOF (set-of-frames) file from all raw calibrations, master calibrations and raw science frames:

```

/path_raw/uves_bias_red1.fits      BIAS_RED
/path_raw/uves_bias_red2.fits      BIAS_RED
/path_raw/uves_bias_red3.fits      BIAS_RED
/path_raw/uves_bias_red4.fits      BIAS_RED
/path_raw/uves_bias_red5.fits      BIAS_RED

/path_raw/uves_dark_red1.fits      DARK_RED
/path_raw/uves_dark_red2.fits      DARK_RED
/path_raw/uves_dark_red3.fits      DARK_RED

/path_raw/uves_flat_set1_red1.fits  SFLAT_RED
/path_raw/uves_flat_set1_red2.fits  SFLAT_RED
/path_raw/uves_flat_set1_red3.fits  SFLAT_RED
/path_raw/uves_flat_set2_red1.fits  SFLAT_RED
/path_raw/uves_flat_set2_red2.fits  SFLAT_RED
/path_raw/uves_flat_set2_red3.fits  SFLAT_RED
/path_raw/uves_flat_set3_red1.fits  SFLAT_RED
/path_raw/uves_flat_set3_red2.fits  SFLAT_RED
/path_raw/uves_flat_set3_red3.fits  SFLAT_RED

/path_raw/flames_uves_odd_red.fits  FIB_FF_ODD_RED
/path_raw/flames_uves_even_red.fits FIB_FF_EVEN_RED
/path_raw/flames_uves_all_red.fits  FIB_FF_ALL_RED

/path_raw/flames_uves_arc_form_red.fits FIB_ARC_FORM_RED

/path_raw/flames_uves_order_flat_red.fits  FIB_ORDER_FLAT_RED

/path_raw/flames_uves_arc_lamp_red.fits FIB_ARC_LAMP_RED

/path_raw/flames_uves_science_red.fits  FIB_SCI_RED
/path_raw/flames_uves_science_red.fits  FIB_SCI_SIM_RED
/path_raw/flames_uves_science_red.fits  FIB_SCI_COM_RED

/path_ref/thargood_3.fits           LINE_REFER_TABLE
/path_ref/flxstd.fits               FLUX_STD_TABLE
/path_ref/atmoexan.fits             EXTCOEFF_TABLE

```

Then run

**esorex flames\_obs\_redchain flames\_obs\_redchain.sof**

which will execute the necessary recipes and create the pipeline products listed in following section.

<b>ESO</b>	<b>FLAMES-UVES Pipeline User Manual</b>	Doc:	VLT-MAN-ESO-19500-3016
		Issue:	Issue 5.0
		Date:	Date 2008-11-28
		Page:	30 of 110

In this example five raw bias frames are provided; therefore the `uves_cal_mbias` recipe will be executed in order to produce the master bias frame. For corresponding reasons the `uves_cal_mdark` and `flames_cal_mflat` will be executed. If a master frame (master bias, master dark or master slit flat) is provided instead of the corresponding raw frames, the corresponding master creation recipe (`uves_cal_mbias`, `uves_cal_mdark`, `flames_cal_mflat`) will not be executed.

In the following a typical step-by-step data reduction procedure is described. <sup>8</sup> Figure 5.3.1 gives an overview.

`/path_ref` indicates the full path to the source tree directory containing reference ancillary data, `/path_pro` indicates the full path to the source tree directory containing product data.

Bias frames: these frames are characterized by `DPR.TYPE='BIAS'`

```
/path_raw/uves_bias_red1.fits  BIAS_RED
/path_raw/uves_bias_red2.fits  BIAS_RED
/path_raw/uves_bias_red3.fits  BIAS_RED
/path_raw/uves_bias_red4.fits  BIAS_RED
/path_raw/uves_bias_red5.fits  BIAS_RED
```

Dark frames: these frames are characterized by `DPR.TYPE='DARK'`

```
/path_raw/uves_dark_red1.fits  DARK_RED
/path_raw/uves_dark_red2.fits  DARK_RED
/path_raw/uves_dark_red3.fits  DARK_RED
```

Slit Flat-field lamp frames: these frames are characterized by `DPR.TYPE='LAMP,SFLAT'`

```
/path_raw/uves_flat_set1_red1.fits SFLAT_RED
/path_raw/uves_flat_set1_red2.fits SFLAT_RED
/path_raw/uves_flat_set1_red3.fits SFLAT_RED
/path_raw/uves_flat_set2_red1.fits SFLAT_RED
/path_raw/uves_flat_set2_red2.fits SFLAT_RED
/path_raw/uves_flat_set2_red3.fits SFLAT_RED
/path_raw/uves_flat_set3_red1.fits SFLAT_RED
/path_raw/uves_flat_set3_red2.fits SFLAT_RED
/path_raw/uves_flat_set3_red3.fits SFLAT_RED
```

Formatcheck: these frames are characterized by `DPR.TYPE='LAMP,FMTCHK'`,

```
/path_raw/flames_uves_arc_lamp_form_blue.fits FIB_ARC_LAMP_FORM_RED
```

Single fibre order tracing flat-field frames: these frames are characterized by `DPR.TYPE='LAMP,ORDERDEF,SimCal'`

---

<sup>8</sup>The procedure using *Gasgano* is conceptually identical.

<b>ESO</b>	<b>FLAMES-UVES Pipeline User Manual</b>	Doc:	VLT-MAN-ESO-19500-3016
		Issue:	Issue 5.0
		Date:	Date 2008-11-28
		Page:	31 of 110

```
/path_raw/flames_uves_order_flat_blue.fits FIB_ORDER_FLAT_RED
```

Arc lamp frames to compute the wavelength calibration: these frames have DPR.TYPE respectively equal to 'LAMP,WAVE,OzPoz'.

```
/path_raw/flames_uves_arc_lamp_red.fits FIB_ARC_LAMP_RED
```

science frames: these frames are characterized by DPR.TYPE='OBJECT'.

```
/path_raw/flames_uves_science_red.fits FIB_SCI_RED
```

In the following examples we assume that pipeline product filenames are the original as set by the pipeline. This corresponds to have the parameter **esorex.caller.suppress-prefix** in the EsoRex configuration file (\$HOME/.esorex/esorex.rc) set to TRUE. Otherwise EsoRex will rename the pipeline products using a common prefix (set by the parameter **esorex.caller.output-prefix**), a four digits increasing number, and terminating the FITS file with the extension “.fits”. We suggest to verify to have the flag *readonly* set to FALSE, if the user would like to run the same recipe several times with EsoRex having standard values for product files. This setting allows the pipeline to overwrite previously generated products <sup>9</sup>.

1. Generated master flat frames. One selects the raw biases and lists them in an ASCII file uves\_cal\_mbias.sof:

```
/path_raw/uves_bias_red1.fits BIAS_RED
/path_raw/uves_bias_red2.fits BIAS_RED
/path_raw/uves_bias_red3.fits BIAS_RED
/path_raw/uves_bias_red4.fits BIAS_RED
/path_raw/uves_bias_red5.fits BIAS_RED
```

The command:

**esorex uves\_cal\_mbias uves\_cal\_mbias.sof**

will generate the following products:

default recipe filename	format	PRO.CATG	short description
masterbias_redl	2d image (pix-pix)	MASTER_BIAS_REDL	master bias
masterbias_redu	2d image (pix-pix)	MASTER_BIAS_REDU	master bias

and generates the ASCII file masterbias\_redl-0.paf, masterbias\_redu-0.paf which logs QC parameters.

**mv \*.fits \*.paf /path\_pro**

---

<sup>9</sup>By default installation in the EsoRex configuration file (\$HOME/.esorex/esorex.rc) the flag *suppress-prefix* is set to FALSE and the flag *readonly* is set to FALSE, a possible combination, in which case pipeline product filenames will have a prefix out\_, an increasing four digit number, and extension .fits.



<b>ESO</b>	<b>FLAMES-UVES Pipeline User Manual</b>	Doc:	VLT-MAN-ESO-19500-3016
		Issue:	Issue 5.0
		Date:	Date 2008-11-28
		Page:	32 of 110

- Then one selects the raw darks and lists them in an ASCII file `uves_cal_mdark.sof`.

```

/path_raw/uves_dark_red1.fits  DARK_RED
/path_raw/uves_dark_red2.fits  DARK_RED
/path_raw/uves_dark_red3.fits  DARK_RED
/path_pro/masterbias_red1.fits MASTER_BIAS_REDL (optional)
/path_pro/masterbias_redu.fits MASTER_BIAS_REDU (optional)

```

The command:

**esorex uves\_cal\_mdark uves\_cal\_mdark.sof**

will generate the following products:

default recipe filename	format	PRO.CATG	short description
masterdark_redl	2d image (pix-pix)	MASTER_DARK_REDL	master dark
masterdark_redu	2d image (pix-pix)	MASTER_DARK_REDU	master dark

and generates the ASCII file `masterdark_redl-0.paf`, `masterdark_redu-0.paf` which logs QC parameters.

**mv \*.fits \*.paf /path\_pro**

- Generate order and line guesses. Formatcheck frames are listed together with the required calibration frames in an ASCII file, `flames_cal_predict.sof`. This file will look like as follows:

```

/path_raw/flames_uves_arc_lamp_form_red.fits FIB_ARC_LAMP_FORM_RED
/path_ref/thargood_3.fits                     LINE_REFER_TABLE
/path_pro/uves_masterbias_red1.fits           MASTER_BIAS_REDL (optional)
/path_pro/uves_masterbias_redu.fits           MASTER_BIAS_REDU (optional)

```

Then the user can generate the order and line guesses with the command

**esorex flames\_cal\_predict flames\_cal\_predict.sof**

This command will generate four files (in the following table FITS files have extension `.fits`):

default recipe filename	format	PRO.CATG	short description
lineguesstable_redl	table	FIB_LIN_GUE_REDL	guess line table
orderguesstable_redl	table	FIB_ORD_GUE_REDL	guess order table
lineguesstable_redu	table	FIB_LIN_GUE_REDU	guess line table
orderguesstable_redu	table	FIB_ORD_GUE_REDU	guess order table

and generates the ASCII files `lineguesstable_redl-0.paf`, and `lineguesstable_redu-0.paf` which logs QC parameters.

**mv \*.fits \*.paf /path\_pro**

- Generate order table and order definition frames. Order definition frames are necessary later in the data reduction to have a full set of flat frames that cover also the simultaneous calibration fibre as the odd-even



<b>ESO</b>	<b>FLAMES-UVES Pipeline User Manual</b>	Doc:	VLT-MAN-ESO-19500-3016
		Issue:	Issue 5.0
		Date:	Date 2008-11-28
		Page:	33 of 110

fibre frames cover only the other fibres.

A set of narrow slit raw flat-field frames may be put in the ASCII file flames\_cal\_orderpos.sof.

```
/path_raw/flames_uves_order_flat_red.fits FIB_ORDEF_RED
/path_pro/orderguesstable_redl.fits FIB_ORD_GUE_REDL
/path_pro/orderguesstable_redu.fits FIB_ORD_GUE_REDU
/path_pro/masterbias_redl.fits      MASTER_BIAS_REDL  (optional)
/path_pro/masterbias_redu.fits      MASTER_BIAS_REDU  (optional)
```

The user can generate order tables and reoriented single fibre order flats with the command:

**esorex flames\_cal\_orderpos flames\_cal\_orderpos.sof**

This command will generate the following products:

default recipe filename	format	PRO.CATG	short description
ordertable_redl	table	FIB_ORD_TAB_REDL	simultaneous calibration fibre order table
order_def_redl	image	FIB_ORDEF_REDL	simultaneous calibration fibre order frame
ordertable_redu	table	FIB_ORD_TAB_REDU	simultaneous calibration fibre order table
order_def_redu	image	FIB_ORDEF_REDU	simultaneous calibration fibre order frame

and generates the ASCII files ordertable\_redl-0.paf, ordertable\_redu-0.paf which logs QC parameters.

**mv \*.fits \*.paf /path\_pro**

- Then one selects the raw slit flat-field frames and lists them in an ASCII file flames\_cal\_mflat.sof together with some master calibrations and previously obtained products:

```
/path_raw/flames_uves_flat_set1_red1.fits SFLAT_RED
/path_raw/flames_uves_flat_set1_red2.fits SFLAT_RED
/path_raw/flames_uves_flat_set1_red3.fits SFLAT_RED
/path_raw/flames_uves_flat_set2_red1.fits SFLAT_RED
/path_raw/flames_uves_flat_set2_red2.fits SFLAT_RED
/path_raw/flames_uves_flat_set2_red3.fits SFLAT_RED
/path_raw/flames_uves_flat_set3_red1.fits SFLAT_RED
/path_raw/flames_uves_flat_set3_red2.fits SFLAT_RED
/path_raw/flames_uves_flat_set3_red3.fits SFLAT_RED
/path_pro/masterbias_redl.fits      MASTER_BIAS_REDL  (optional)
/path_pro/masterbias_redu.fits      MASTER_BIAS_REDU  (optional)
```

The command:

**esorex flames\_cal\_mflat flames\_cal\_mflat.sof**

will generate the following products:

default recipe filename	format	PRO.CATG	Note
-------------------------	--------	----------	------

<b>ESO</b>	<b>FLAMES-UVES Pipeline User Manual</b>	Doc:	VLT-MAN-ESO-19500-3016
		Issue:	Issue 5.0
		Date:	Date 2008-11-28
		Page:	34 of 110

masterflat_set1_redl	2d image (pix-pix)	MASTER_SFLAT_REDL	master flat
masterflat_set2_redl	2d image (pix-pix)	MASTER_SFLAT_REDL	master flat
masterflat_set3_redl	2d image (pix-pix)	MASTER_SFLAT_REDL	master flat
masterflat_set1_redu	2d image (pix-pix)	MASTER_SFLAT_REDU	master flat
masterflat_set2_redu	2d image (pix-pix)	MASTER_SFLAT_REDU	master flat
masterflat_set3_redu	2d image (pix-pix)	MASTER_SFLAT_REDU	master flat

and generates the ASCII file masterflat\_redl-0.paf, masterflat\_redu-0.paf which log QC parameters.

**mv \*.fits \*.paf /path\_pro**

- Then the order table for all fibre traces is determined. A set of three raw fibre frames illuminating respectively the odd, even and all fibre frames are put in an the ASCII file flames\_cal\_prep\_sff\_ofpos.sof.

```

/path_raw/flames_uves_odd_red.fits    FIB_FF_ODD_RED
/path_raw/flames_uves_even_red.fits   FIB_FF_EVEN_RED
/path_raw/flames_uves_all_red.fits    FIB_FF_ALL_RED
/path_pro/masterbias_redl.fits        MASTER_BIAS_REDL (optional)
/path_pro/order_def_redl.fits         FIB_ORDEF_REDL
/path_pro/ordertable_redl.fits        FIB_ORD_TAB_REDL
/path_pro/masterflat_set1_redl.fits   MASTER_SFLAT_REDL
/path_pro/masterflat_set2_redl.fits   MASTER_SFLAT_REDL
/path_pro/masterflat_set2_redl.fits   MASTER_SFLAT_REDL
/path_pro/masterbias_redu.fits        MASTER_BIAS_REDU (optional)
/path_pro/order_def_redu.fits         FIB_ORDEF_REDU
/path_pro/ordertable_redu.fits        FIB_ORD_TAB_REDU
/path_pro/masterflat_set1_redu.fits   MASTER_SFLAT_REDU
/path_pro/masterflat_set2_redu.fits   MASTER_SFLAT_REDU
/path_pro/masterflat_set2_redu.fits   MASTER_SFLAT_REDU

```

Note that the input master bias frames are optional but recommended input.

The command

**esorex flames\_cal\_prep\_sff\_ofpos flames\_cal\_prep\_sff\_ofpos.sof**

will generate the following products:

default recipe filename	format	PRO.CATG	short description
xt_odd_l	table	FIB_FF_ODD_INFO_TAB	info table
xt_even_l	table	FIB_FF_EVEN_INFO_TAB	info table
xt_all_l	table	FIB_FF_ODD_INFO_TAB	info table
slitff_common_redl	imagelist	SLIT_FF_COM_REDL	slitff common frame
slitff_norm_redl	image	SLIT_FF_NOR_REDL	slitff common frame
slitff_dtc_redl	imagelist	SLIT_FF_DTC_REDL	slitff data frame
slitff_sgc_redl	imagelist	SLIT_FF_SGC_REDL	slitff sigma frame
slitff_bpc_redl	imagelist	SLIT_FF_BPC_REDL	slitff badpixel frame

<b>ESO</b>	<b>FLAMES-UVES Pipeline User Manual</b>	Doc:	VLT-MAN-ESO-19500-3016
		Issue:	Issue 5.0
		Date:	Date 2008-11-28
		Page:	35 of 110

slitff_bnc_redl	imagelist	SLIT_FF_BNC_REDL	slitff frames boundaries
orfl	table	FIB_ORDEF_TABLE_REDL	fibre order table
fibreff_common_redl	imagelist	FIB_FF_COM_REDL	fibff common frame
fibreff_norm_redl	imagelist	FIB_FF_NOR_REDL	fibff common frame
fibreff_nsigma_redl	imagelist	FIB_FF_NSIG_REDL	Normalisation sigmas for fibre FF
fibreff_dtc_redl	imagelist	FIB_FF_DTC_REDL	fibff data frame
fibreff_sgc_redl	imagelist	FIB_FF_SGC_REDL	fibff sigma frame
fibreff_bpc_redl	imagelist	FIB_FF_BPC_REDL	fibff badpixel frame
slitff_common_redu	imagelist	SLIT_FF_COM_REDU	slitff common frame
slitff_norm_redu	image	SLIT_FF_NOR_REDU	slitff common frame
slitff_dtc_redu	imagelist	SLIT_FF_DTC_REDU	slitff data frame
slitff_sgc_redu	imagelist	SLIT_FF_SGC_REDU	slitff sigma frame
slitff_bpc_redu	imagelist	SLIT_FF_BPC_REDU	slitff badpixel frame
slitff_bnc_redu	imagelist	SLIT_FF_BNC_REDU	slitff frames boundaries
ordu	table	FIB_ORDEF_TABLE_REDU	fibre order table
fibreff_common_redu	imagelist	FIB_FF_COM_REDU	fibff common frame
fibreff_norm_redu	imagelist	FIB_FF_NOR_REDU	fibff common frame
fibreff_nsigma_redu	imagelist	FIB_FF_NSIG_REDU	Normalisation sigmas for fibre FF
fibreff_dtc_redu	imagelist	FIB_FF_DTC_REDU	fibff data frame
fibreff_sgc_redu	imagelist	FIB_FF_SGC_REDU	fibff sigma frame
fibreff_bpc_redu	imagelist	FIB_FF_BPC_REDU	fibff badpixel frame

and generates the ASCII file masterflat\_redl-0.paf, masterflat\_redu-0.paf which log QC parameters.

**mv \*.fits \*.paf /path\_pro**

7. Then the wavelength calibration is performed. A raw frame illuminated by the long slit through the fibres by an arc lamp are put in the ASCII file flames\_cal\_wavecalf.sof.

```

/path_raw/flames_uves_arc_lamp_red.fits    FIB_ARC_LAMP_RED
/path_pro/ordertable_redl.fits             FIB_ORDEF_TABLE_REDL
/path_pro/lineguesstable_redl.fits         FIB_LIN_GUE_REDL
/path_pro/masterbias_redl.fits             MASTER_BIAS_REDL  (optional)
/path_pro/ordertable_redu.fits             FIB_ORDEF_TABLE_REDU
/path_pro/lineguesstable_redu.fits         FIB_LIN_GUE_REDU
/path_pro/masterbias_redu.fits             MASTER_BIAS_REDU  (optional)
/path_ref/thargood_3.fits                  LINE_REFER_TABLE

```

Note that the input master bias and master flat frames are optional but recommended inputs.

The command

**esorex flames\_cal\_wavecalf flames\_cal\_wavecalf.sof**

will generate the following products:

default recipe filename	format	PRO.CATG	short description
linetable_redl	table	FIB_LINE_TABLE_REDL	line table
linetable_redu	table	FIB_LINE_TABLE_REDU	line table

<b>ESO</b>	<b>FLAMES-UVES Pipeline User Manual</b>	Doc:	VLT-MAN-ESO-19500-3016
		Issue:	Issue 5.0
		Date:	Date 2008-11-28
		Page:	36 of 110

This table contains the solutions for each extracted fibres in several extensions according to the following schema:

```
Line table for fibre 1, saved to extensions 1-3 of 'linetable_redl/redu'
Line table for fibre 2, saved to extensions 4-6 of 'linetable_redl/redu'
Line table for fibre 3, saved to extensions 7-9 of 'linetable_redl/redu'
Line table for fibre 4, saved to extensions 10-12 of 'linetable_redl/redu'
Line table for fibre 5, saved to extensions 13-15 of 'linetable_redl/redu'
Line table for fibre 6, saved to extensions 16-18 of 'linetable_redl/redu'
Line table for fibre 7, saved to extensions 19-21 of 'linetable_redl/redu'
Line table for fibre 8, saved to extensions 22-24 of 'linetable_redl/redu'
Line table for fibre 9, saved to extensions 25-27 of 'linetable_redl/redu'
```

and generates the ASCII file linetable\_redl-0.paf, linetable\_redu-0.paf which logs QC parameters.

**mv \*.fits \*.paf /path\_pro**

- Finally, the raw science frame is reduced. The raw science frame is listed together with master calibration products in the following ASCII file flames\_obs\_scired.sof:

```
/path_raw/flames_uves_science_red.fits FIB_SCI_RED
/path_pro/masterbias_redl.fits MASTER_BIAS_REDL (optional)
/path_pro/slitff_com_redl.fits SLIT_FF_COM_REDL
/path_pro/slitff_nor_redl.fits SLIT_FF_NOR_REDL
/path_pro/slitff_dtc_redl.fits SLIT_FF_DTC_REDL
/path_pro/slitff_bpc_redl.fits SLIT_FF_BPC_REDL
/path_pro/slitff_sgc_redl.fits SLIT_FF_SGC_REDL
/path_pro/fibff_com_redl.fits FIB_FF_COM_REDL
/path_pro/fibff_nor_redl.fits FIB_FF_NOR_REDL
/path_pro/fibff_dtc_redl.fits FIB_FF_DTC_REDL
/path_pro/fibff_bpc_redl.fits FIB_FF_BPC_REDL
/path_pro/fibff_sgc_redl.fits FIB_FF_SGC_REDL
/path_pro/fibff_bnc_redl.fits FIB_FF_BNC_REDL
/path_pro/orfl.fits FIB_ORDEF_TABLE_REDL
/path_pro/linetable_redl.fits FIB_LINE_TABLE_REDL
/path_pro/masterbias_redu.fits MASTER_BIAS_REDU (optional)
/path_pro/slitff_com_redu.fits SLIT_FF_COM_REDU
/path_pro/slitff_nor_redu.fits SLIT_FF_NOR_REDU
/path_pro/slitff_dtc_redu.fits SLIT_FF_DTC_REDU
/path_pro/slitff_bpc_redu.fits SLIT_FF_BPC_REDU
/path_pro/slitff_sgc_redu.fits SLIT_FF_SGC_REDU
/path_pro/fibff_com_redu.fits FIB_FF_COM_REDU
/path_pro/fibff_nor_redu.fits FIB_FF_NOR_REDU
/path_pro/fibff_dtc_redu.fits FIB_FF_DTC_REDU
/path_pro/fibff_bpc_redu.fits FIB_FF_BPC_REDU
/path_pro/fibff_sgc_redu.fits FIB_FF_SGC_REDU
/path_pro/fibff_bnc_redu.fits FIB_FF_BNC_REDU
/path_pro/orfu.fits FIB_ORDEF_TABLE_REDU
/path_pro/linetable_redu.fits FIB_LINE_TABLE_REDU
```

<b>ESO</b>	<b>FLAMES-UVES Pipeline User Manual</b>	Doc:	VLT-MAN-ESO-19500-3016
		Issue:	Issue 5.0
		Date:	Date 2008-11-28
		Page:	37 of 110


The command:

**esorex flames\_obs\_scired flames\_obs\_scired.sof**

will generate the following products:

default recipe file name	format	PRO.CATG	short description
bin_table_info_redl	2d table	FIB_SCI_INFO_TAB_REDL	table with fibre information
fxb_l_raw000i	2d (pix-ord) image	XB_SCI_RAW_REDL	extracted, flatfielded raw frame
fxb_l_rawsig000i	2d (pix-ord) image	ERR_XB_SCI_RAW_REDL	error frame
wfxb_l_raw000i	2d (wav-ord) image	WXB_SCI_RAW_REDL	rebinned, extracted, flatfielded raw frame
wfxb_l_rawsig000i	2d (wav-ord) image	ERR_WXB_SCI_RAW_REDL	error frame
mwfxb_l_rawsig000i	1d (wav) image	MWXB_SCI_RAW_REDL	merged, rebinned, extracted, flat-fielded raw frame
mwfxb_l_raw000i_sigma	1d (wav) image	ERR_MWXB_SCI_RAW_REDL	merged, rebinned, extracted, flat-fielded error frame
fxb_l_000i	2d (pix-ord) image	XB_SCI_REDL	extracted, flatfielded raw frame
fxb_l_sig000i	2d (pix-ord) image	ERR_XB_SCI_REDL	error frame
wfxb_l_000i	2d (wav-ord) image	WXB_SCI_REDL	rebinned, extracted, flatfielded raw frame
wfxb_l_sig000i	2d (wav-ord) image	ERR_WXB_SCI_REDL	error frame
mwfxb_l_000i	1d (wav) image	MWXB_SCI_REDL	merged, rebinned, extracted, flat-fielded raw frame
mwfxb_l_000i_sigma	1d (wav) image	ERR_MWXB_SCI_REDL	merged, rebinned, extracted, flat-fielded error frame
bin_table_info_redu	2d table	FIB_SCI_INFO_TAB_REDU	table with fibre information
fxb_u_raw000i	2d (pix-ord) image	XB_SCI_RAW_REDU	extracted, flatfielded raw frame
fxb_u_rawsig000i	2d (pix-ord) image	ERR_XB_SCI_RAW_REDU	error frame
wfxb_u_raw000i	2d (wav-ord) image	WXB_SCI_RAW_REDU	rebinned, extracted, flatfielded raw frame
wfxb_u_rawsig000i	2d (wav-ord) image	ERR_WXB_SCI_RAW_REDU	error frame
mwfxb_u_rawsig000i	1d (wav) image	MWXB_SCI_RAW_REDU	merged, rebinned, extracted, flat-fielded raw frame
mwfxb_u_raw000i_sigma	1d (wav) image	ERR_MWXB_SCI_RAW_REDU	merged, rebinned, extracted, flat-fielded error frame
fxb_u_000i	2d (pix-ord) image	XB_SCI_REDU	extracted, flatfielded raw frame
fxb_u_sig000i	2d (pix-ord) image	ERR_XB_SCI_REDU	error frame
wfxb_u_000i	2d (wav-ord) image	WXB_SCI_REDU	rebinned, extracted, flatfielded raw frame
wfxb_u_sig000i	2d (wav-ord) image	ERR_WXB_SCI_REDU	error frame
mwfxb_u_sig000i	1d (wav) image	MWXB_SCI_REDU	merged, rebinned, extracted, flat-fielded raw frame
mwfxb_u_000i_sigma	1d (wav) image	ERR_MWXB_SCI_REDU	merged, rebinned, extracted, flat-fielded error frame

and generates the ASCII file resampled\_ff\_blue-0.paf which logs QC parameters.

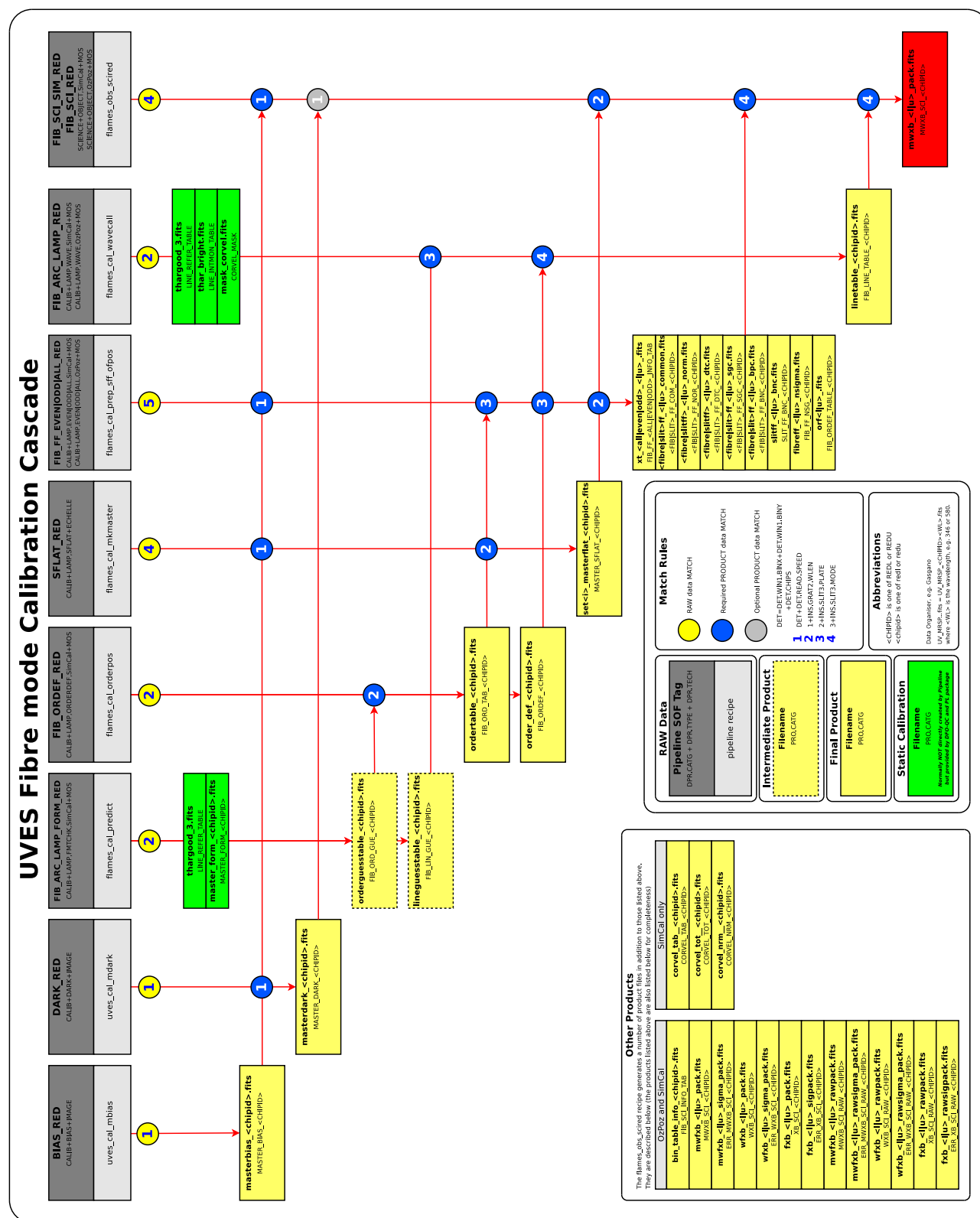

<b>ESO</b>	<b>FLAMES-UVES Pipeline User Manual</b>	Doc:	VLT-MAN-ESO-19500-3016
		Issue:	Issue 5.0
		Date:	Date 2008-11-28
		Page:	39 of <a href="#">110</a>

## 6 Known problems

We suggest the user to execute the data reduction recipes using parameter defaults and all the reference and master calibrations indicated in this manual. The following is a list of currently-known problems with FLAMES-UVES recipes, and workarounds, if available:

- The radial velocity correlation has not yet been verified to give the same results as the ones provided by FLAMES-UVES MIDAS based pipeline. In particular the zero points, are not yet computed during wavelength calibration. The user interested on those results is invited to run the MIDAS based pipeline release.
- Guplots do not work with Gasgano.
- Recipes `flames_cal_prep_sff_ofpos` and `flames_obs_scired` are not portable on 64-bit platforms.

<b>ESO</b>	<b>FLAMES-UVES Pipeline User Manual</b>	Doc:	VLT-MAN-ESO-19500-3016
		Issue:	Issue 5.0
		Date:	Date 2008-11-28
		Page:	40 of <a href="#">110</a>

## 7 Instrument Data Description

To reduce a science frame, the following calibration frames are needed:

- Bias(es)
- Dark(s)
- Fibre formatcheck
- Single Fibre Order definition
- Slit Flat-field(s)
- Fibre flat-field frames having illuminated: odd, even, all fibres
- Fibre wavelength calibration

It is also necessary to have handy a reference (ThAr) line table. The table thargood\_3.tfits is included in the pipeline installation kit. The table thargood\_< *n* >.tfits (where *n* is 2 or 3) was also included in Pre P82 Service mode packages. The file UV\_GLRE\_<YYMMDD>A\_line\_refer\_table.fits is included in Service mode packages for P82 onward.

### 7.1 Bias frames

Bias frames give the read out of the CCD detector of zero integration time with the shutter closed. Usually they are taken as a set of five exposures from which, through stacking, a Master Bias is created thus reducing the read out noise. This needs to be subtracted for example from the science frame to get the signal contribution from the source only.

### 7.2 Dark frames

Dark frames are measured occasionally, with the shutter closed. They are used to measure the dark current. They are measured for 1x1 binning with typical exposure times of 1h. There are also open-shutter DARKs (since December 2001). They include, in addition to the CCD dark current, contributions from the camera enclosure. Typical values are reported on the ESO Website under [www.eso.org/observing/dfo/quality/UVES/qc/dark\\_qc1.html](http://www.eso.org/observing/dfo/quality/UVES/qc/dark_qc1.html). As the contribution of UVES dark exposures may be considered, in first approximation, negligible, they can be excluded from the data reduction chain as we will assume here in the following section.

### 7.3 Slit flat-field frames

Slit flat-field frames are long slit exposures taken with a continuum lamp. They give information on the response of the detector, allowing to measure variations in efficiency at small (pixel-to-pixel), intermediate (fringing, in the far red) and large (the blaze function) scale. Thus they are usually taken as three sets of three frames each at a



<b>ESO</b>	<b>FLAMES-UVES Pipeline User Manual</b>	Doc:	VLT-MAN-ESO-19500-3016
		Issue:	Issue 5.0
		Date:	Date 2008-11-28
		Page:	41 of 110

given cross disperser setting, which after bias subtraction are stacked in a master to reject statistical outliers like cosmic ray events. The final Master flat-field is also background subtracted to eliminate diffused light from the orders in the inter-order regions. Science frames need to be corrected for pixel-to-pixel variations, interference fringes and the blaze function through division by the master flat-field frame.

To properly evaluate the background level and avoid overlapping of orders, each slit flat-field frame has to be taken with a limited slit aperture. To cover all the fibres (allowing also for some variability in their placement) several slit flat-field frames are thus necessary, taken with different offsets with respect to a central position and overlapping between each other.

## 7.4 Formatcheck frames

It is necessary to have a special formatcheck frame taken by illuminating the simultaneous calibration fibre with a ThAr lamp. This is used in combination with a physical model of UVES and the information contained in the FITS header and in a ThAr reference line table, and the known offset of the simultaneous calibration fibre from the UVES slit's center, to find a "guess" solution of the spectral format (order locations and wavelength calibration). This allows the user later on to obtain robust and automatic spectral format solutions.

## 7.5 Single fibre order definition frames

A single fibre order definition frame is a calibration exposure obtained with the calibration fibre illuminated by a continuum lamp. It is a very high signal-to-noise ratio fibre-echelle frame describing accurately the simultaneous calibration fibre order locations.

## 7.6 Fibre Flat-Field (odd-even-all) frames

The fibre flat-fields serve to locate the fibre positions, both relative to one another and absolute, within some tolerance; to know the cross-sectional profile of each fibre at each wavelength; and to determine the relative throughput of different fibres. These fibre flat-fields are only usable if the fibre traces fall in regions illuminated in at least one of the slit flat-fields; fibres whose light falls, even in part, in dark regions of the slit flat-fields are discarded. The data reduction software considers as really flat only a subregion of each slit flat-field frame discarding a few pixels at the flat-field border.

The fibres' images on the detector are so closely packed that the cross-sectional profiles of adjacent fibres overlap to some non-negligible extent. Therefore, a single flat-field frame, would not yield enough information on individual fibres to perform a correct data reduction. Instead, the pipeline requires fibre flat-fields on which the fibres can be clearly separated, such as one containing only odd-numbered fibres and another containing only even-numbered fibres.

The order definition frame, if taken with the simultaneous fibre, could also be used as an additional flat-field frame covering the simultaneous calibration fibre. Such inclusion is necessary to be able to extract the simultaneous calibration fibre. In such a case one would add it into the pool of odd-even fibre flat fields (see below), and this set of frames would thus contain completely separable images of all fibres.

<b>ESO</b>	<b>FLAMES-UVES Pipeline User Manual</b>	Doc:	VLT-MAN-ESO-19500-3016
		Issue:	Issue 5.0
		Date:	Date 2008-11-28
		Page:	42 of 110

In addition to the above, an all-fibre flat-field is needed to correct for relative throughput differences between odd and even fibres. This can be replaced by a Nasmyth screen all-fibre flat-field for higher accuracy.

If it is planned to observe using the simultaneous calibration fibre, it is also necessary to have in the calibration data a ThAr frame with the ThAr lamp feeding all the same fibres used during the night (7OzPoz+1SimCal) as well as an all fibre flat-field with the flat-field lamp feeding all the same fibres used during the night (7OzPoz+1SimCal). This is the minimum fibre flat-field set usable by the pipeline.

## 7.7 Wavelength calibration frames

Wavelength calibration frames are fibre exposures taken (for FLAMES-UVES) illuminating all fibres with a ThAr arc lamp. They are used to find the wavelength calibration solution.

The wavelength calibration frames must be taken through the same set of fibres as the science frame.

We would like to emphasize the importance of having “correct” calibration data to be able to reduce the science observations with best quality. For example, if because of thermal shifts the available slit flat-fields do not cover the position where a (science-frame) fibre is located, this fibre will be discarded and the pipeline will be unable to extract its signal. The same is true if fibre positions in the science-frame are shifted too much with respect to the corresponding positions of the fibre flat-field frames used to do the data reduction. This may occur in case of automatic online pipeline data reduction on Paranal where usually the calibration data base is updated with several months time scale. Moreover, the cross-sectional profiles of the fibres may change e. g. due to slightly different focusing. Since this, in turn, changes the contamination between adjacent fibres, a significantly different profile between calibrations and science data will lead to an incorrect deconvolution of the fibres. The pipeline searches for possible instrument shifts along Y (i. e. perpendicular to dispersion) in the interval [-6,6] pixel. It satisfactorily reduced commissioning data which, due a not completely stable instrument setup at the time, were characterized by shifts between the science observation and the calibration in the range [-3.5,3.5]. It must be said that the source of such Y shift was discovered and fixed during an instrument software upgrade conducted before the beginning of FLAMES operations and so we do not expect the user will receive data affected by shifts between science and calibrations of more than one pixel. The pipeline allows to interactively recover extraction also in cases of shifts wider than 6 pixels. This can be done by executing the science recipe once with parameter defaults; then by plotting the information available in the product `cor_table_chip.fits` (`chip='l'` or `'u'`) the user can asses the proper value of the correlation shift, and set correspondly the correlation offset in order to reduce, within a narrower correlation range centered on the given offset, the given science data in a following data reduction iteration.

The FLAMES positioner uses two plates and two distinct sets of fibres, with two distinct entrance positions in UVES.

A simultaneous calibration fibre for precise radial velocity measurements is also available. Using this fibre the spectral format changes, including (at most) 7 OzPoz fibres and 1 SimCal fibre. Therefore, only calibration data need to be taken with the same FLAMES plate and ( $\leq 8$  OzPoz fibres)/( $\leq 7$  OzPoz fibres+1 SimCal fibre) configuration as the actual science frames to be reduced should be used in the reduction.

Although we introduced means to detect and appropriately treat saturated pixels in the data reduction, it might happen that a saturated frame affects the data reduction. More specifically, while isolated saturated pixels can be easily detected by the pipeline and excluded from the subsequent data reduction, strongly saturated emission lines can produce a very irregular distribution of scattered light on the frame, causing a poor background subtrac-

<b>ESO</b>	<b>FLAMES-UVES Pipeline User Manual</b>	Doc:	VLT-MAN-ESO-19500-3016
		Issue:	Issue 5.0
		Date:	Date 2008-11-28
		Page:	43 of 110

tion and, as a consequence, wrong results in the determination of shifts to be applied to each fibre (correlation step). In particularly pathological cases, this may lead to a failure of the automatic data reduction procedure, and require manual intervention to properly extract the data.

## 7.8 Supported raw frames (keyword identifiers)

In this section we describe all possible types of raw frames for the different observing modes. More information on those data may be found on [www.eso.org/qc/uves](http://www.eso.org/qc/uves). The different frame types can be identified by the values of the DPR keywords of their FITS headers (see [5]). These keywords are generated by the FLAMES-UVES templates (for a description of the FLAMES-UVES templates see [9]). A given frame type can be processed by one or several different dedicated pipeline recipes. The individual pipeline recipes are described in section 11. In most cases, calibration data frames are needed to reduce a given frame. These calibration data have to match the input frame in a number of instrument parameters (e.g. to apply a flat-field correction to a science frame only a flat-field frame taken in the same central wavelength, same slit length, etc. should be used for the correction). These parameters are listed under *relevant instrument parameters*.

The following raw frame types are possible:

- Relevant instrument parameters group 1 common to all raw frames:

Conversion e $\rightarrow$ ADU:                      ESO DET OUT1 CONAD

- Relevant instrument parameters group 2 common to some raw frames:

Grating used:                              ESO INS GRAT2 ID  
Central wavelength used:                ESO INS GRAT2 WLEN  
Observation Plate used:                INS OBSPLATE  
Simultaneous calibration:               OCS SIMCAL

## Bias frames

- Template signature:  
UVES\_x\_cal\_bias  
(x: red, dic1, dic2)
- DPR keywords:  
ESO DPR CATG = CALIB  
ESO DPR TYPE = BIAS ESO DPR TECH = IMAGE
- DO category:  
BIAS\_RED
- Pipeline recipe: `uves_cal_mbias`

<b>ESO</b>	<b>FLAMES-UVES Pipeline User Manual</b>	Doc:	VLT-MAN-ESO-19500-3016
		Issue:	Issue 5.0
		Date:	Date 2008-11-28
		Page:	44 of <a href="#">110</a>

- Relevant instrument parameters: group 1.
- Reference: Figure [7.8.0](#) (a).

### Dark frames

- Template signature:  
UVES\_x\_cal\_dark  
(x: red, dic1, dic2)
- DPR keywords:  
ESO DPR CATG = CALIB  
ESO DPR TYPE = DARK  
ESO DPR TECH = IMAGE
- DO category:  
DARK\_RED
- Pipeline recipe: uves\_cal\_mdark
- Relevant instrument parameters: group 1.

### Order definition flat-field frames

- Template signature:  
UVES\_x\_tec\_orderdef  
(x: red, dic1, dic2)
- DPR keywords:  
ESO DPR CATG = CALIB  
ESO DPR TYPE = LAMP, ORDERDEF, SimCal  
ESO DPR TECH = MOS
- DO category:  
ORDER\_FLAT\_RED
- Pipeline recipe: flames\_cal\_orderpos
- Relevant instrument parameters groups 1, 2.
- Reference: Figure [7.8.0](#) (c).

### Spectroscopic Slit flat-field frames

- Template signature:  
UVES\_x\_cal\_y  
(x: red, dic1, dic2)  
(y: flatatt, flatfree)

<b>ESO</b>	<b>FLAMES-UVES Pipeline User Manual</b>	Doc:	VLT-MAN-ESO-19500-3016
		Issue:	Issue 5.0
		Date:	Date 2008-11-28
		Page:	45 of 110

- DPR keywords:  
ESO DPR CATG = CALIB  
ESO DPR TYPE = LAMP,SFLAT  
ESO DPR TECH = ECHELLE
- DO category:  
FLAT\_RED
- Pipeline recipe: flames\_cal\_mflat
- Relevant instrument parameters groups 1, 2.
- Reference: Figure 7.8.0 (e).

### Format check spectra

- Template signature:  
UVES\_x\_tec\_fmtchk  
(x: red, dic1, dic2)
- DPR keywords:  
ESO DPR CATG = CALIB  
ESO DPR TYPE = LAMP,FMTCHK,SimCal  
ESO DPR TECH = MOS
- DO category:  
FIB\_ARC\_LAMP\_FORM\_RED
- Pipeline recipe: flames\_cal\_predict
- Relevant instrument parameters: groups 1, 2.
- Reference: Figure 7.8.0 (b).

### Odd Fibre flat calibration spectra

- Template signature:  
UVES\_x\_cal\_y  
(x: red, dic1, dic2)  
(y: waveatt, wavefree)
- DPR keywords:  
ESO DPR CATG = CALIB  
ESO DPR TYPE = LAMP,FLAT,ODD  
ESO DPR TECH = MOS
- DO category:  
FIB\_FF\_ODD\_RED

<b>ESO</b>	<b>FLAMES-UVES Pipeline User Manual</b>	Doc:	VLT-MAN-ESO-19500-3016
		Issue:	Issue 5.0
		Date:	Date 2008-11-28
		Page:	46 of 110

- Pipeline recipe: `flames_cal_prep_sff_ofpos`
- Relevant instrument parameters: groups 1, 2.
- Reference: Figure 7.8.0 (d).

### Even Fibre flat calibration spectra

- Template signature:  

```
UVES_x_cal_y
(x: red, dic1, dic2)
(y: waveatt, wavefree)
```
- DPR keywords:  

```
ESO DPR CATG = CALIB
ESO DPR TYPE = LAMP,FLAT,EVEN
ESO DPR TECH = MOS
```
- DO category:  

```
FIB_FF_EVEN_RED
```
- Pipeline recipe: `flames_cal_prep_sff_ofpos`
- Relevant instrument parameters: groups 1, 2.
- Reference: Figure 7.8.0 (d).

### All Fibre flat calibration spectra

- Template signature:  

```
UVES_x_cal_y
(x: red, dic1, dic2)
(y: waveatt, wavefree)
```
- DPR keywords:  

```
ESO DPR CATG = CALIB
ESO DPR TYPE = LAMP,FLAT,ALL
ESO DPR TECH = MOS
```
- DO category:  

```
FIB_FF_ALL_RED
```
- Pipeline recipe: `flames_cal_prep_sff_ofpos`
- Relevant instrument parameters: groups 1, 2.
- Reference: Figure 7.8.0 (d).

<b>ESO</b>	<b>FLAMES-UVES Pipeline User Manual</b>	Doc:	VLT-MAN-ESO-19500-3016
		Issue:	Issue 5.0
		Date:	Date 2008-11-28
		Page:	47 of <a href="#">110</a>

### OzPoz-fibre Wavelength calibration spectra

- Template signature:  

```
UVES_x_cal_y
(x: red, dic1, dic2)
(y: waveatt, wavefree)
```
- DPR keywords:  

```
ESO DPR CATG = CALIB
ESO DPR TYPE = LAMP,WAVE,OzPoz
ESO DPR TECH = MOS
```
- DO category:  

```
FIB_ARC_LAMP_RED
```
- Pipeline recipe: `flames_cal_wavecal`
- Relevant instrument parameters: groups 1, 2.
- Reference: Figure [7.8.0](#) (d).

### Simcal-fibre Wavelength calibration spectra

- Template signature:  

```
UVES_x_cal_y
(x: red, dic1, dic2)
(y: waveatt, wavefree)
```
- DPR keywords:  

```
ESO DPR CATG = CALIB
ESO DPR TYPE = LAMP,WAVE,SimCal
ESO DPR TECH = MOS
```
- DO category:  

```
FIB_ARC_LAMP_RED
```
- Pipeline recipe: `flames_cal_wavecal`
- Relevant instrument parameters: groups 1, 2.
- Reference: Figure [7.8.0](#) (d).



<b>ESO</b>	<b>FLAMES-UVES Pipeline User Manual</b>	Doc:	VLT-MAN-ESO-19500-3016
		Issue:	Issue 5.0
		Date:	Date 2008-11-28
		Page:	48 of <a href="#">110</a>

## Nasmyth screen Fibre flat calibration spectra

- Template signature:  

```
UVES_x_cal_y
(x: red, dic1, dic2)
(y: waveatt, wavefree)
```
- DPR keywords:  

```
ESO DPR CATG = CALIB
ESO DPR TYPE = LAMP,FLAT,NASMYTH
ESO DPR TECH = MOS
```
- DO category:  

```
FIB_FF_NASMYTH_RED
```
- Pipeline recipe: `flames_cal_prep_sff_ofpos`
- Relevant instrument parameters: groups 1, 2.
- Reference: Figure [7.8.0](#) (d).

## Science spectra


- Template signatures:  

```
UVES_x_obs_y
(x: blue, red, dic1, dic2)
(y: exp, expfree)
```
- DPR keywords:  


```
ESO DPR CATG = SCIENCE
ESO DPR TECH = MOS
```
- DO category:  

```
FIB_SCI_RED (ESO DPR TYPE = OBJECT,OzPoz)
FIB_SCI_SIM_RED (ESO DPR TYPE = OBJECT,SimCal )
FIB_SCI_COM_RED (ESO DPR TYPE = OBJECT,OzPoz )
```
- Pipeline recipe: `uves_obs_scired`
- Relevant instrument parameters: groups 1, 2.

<b>ESO</b>	<b>FLAMES-UVES Pipeline User Manual</b>	Doc:	VLT-MAN-ESO-19500-3016
		Issue:	Issue 5.0
		Date:	Date 2008-11-28
		Page:	49 of 110

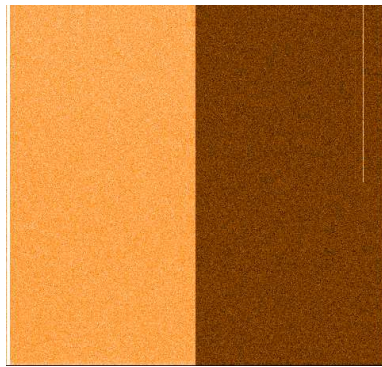

(a raw bias frame)

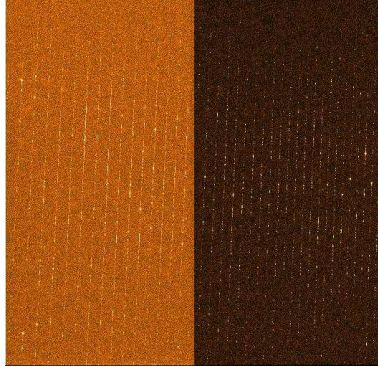

(a raw formatcheck simcal fibre frame)



(a raw order simcal fibre frame)




(a master sflat frame (set1))




(a master sflat frame (set2))



(a master sflat frame (set3))

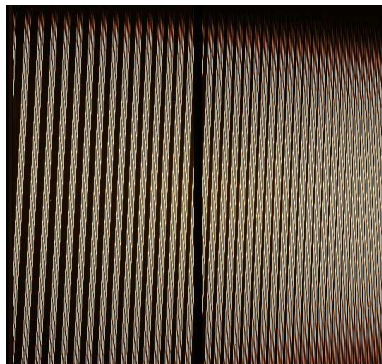


(a raw odd fibres frame)

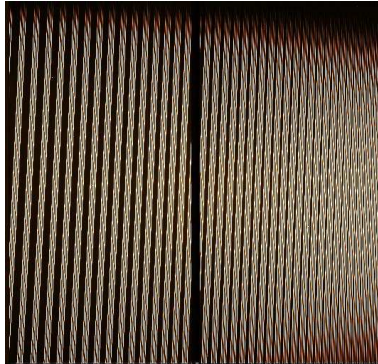


(a raw even fibres frame)



(a raw all fibres frame)

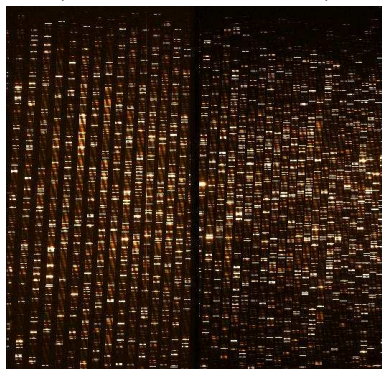

(a raw arc lamp fibre frame)

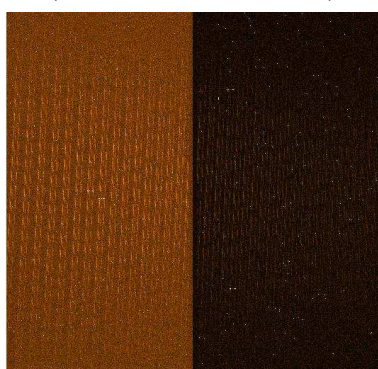

(a raw science frame)

Figure 7.8.0: Several FLAMES-UVES raw frames.

<b>ESO</b>	<b>FLAMES-UVES Pipeline User Manual</b>	Doc:	VLT-MAN-ESO-19500-3016
		Issue:	Issue 5.0
		Date:	Date 2008-11-28
		Page:	50 of 110

## 8 Static Calibration Data

In the following section ancillary data required for FLAMES-UVES data reduction are listed. For each of them we indicate the corresponding value of the HIERARCH ESO PRO CATG, in short PRO.CATG, FITS keyword. This has to be used to identify the frames listed in the *Set of Frames* (see Section 5.2.2, page 24). More information on those data may be found on [www.eso.org/qc/uves](http://www.eso.org/qc/uves).

### 8.1 Line reference table

A reference list of arc lines is necessary to perform the wavelength calibration. Its PRO.CATG is LINE\_REFER\_TABLE. This frame is an input of the recipes `uves_cal_predict` and `uves_cal_wavecal`.

This release of the FLAMES-UVES pipeline provides a new reference line catalog (`thargood_3.fits`)<sup>10</sup> more accurate than the previous one (`thargood_2.fits`). Both catalogs are included but the user is recommended to use the new one.

---

<sup>10</sup>These line reference tables are also included in Service Mode data packages. As of P82 `thargood_3.fits` has been renamed to `UV_GLRE_070222A_line_refer_table.fits`

<b>ESO</b>	<b>FLAMES-UVES Pipeline User Manual</b>	Doc:	VLT-MAN-ESO-19500-3016
		Issue:	Issue 5.0
		Date:	Date 2008-11-28
		Page:	51 of 110

## 9 Data Reduction

In this section, after an overview of the main problems the data reduction needs to solve, we list the required data and the recipes which allow to solve them, giving the data reduction sequence necessary to reduce calibration and science data.

### 9.1 Data reduction overview

In order to fully reduce a set of FLAMES-UVES data, including calibration and science frames, the following data reduction problems need to be solved:

- The detector bias and dark current levels should be measured and subtracted.
- The echelle fiber order traces need to be determined in a robust manner.
- Compute and correct for detector pixel to pixel gain variations and the blaze function.
- For each fibre perform the wavelength calibration in a robust and automatic manner.
- Reduce science fiber data recovering the cross order fibres contamination.

### 9.2 Required input data

To be able to reduce science data one needs to use raw, product and static calibration data together with the pipeline recipes in a given sequence so as to generate all the necessary input to each pipeline recipe at each step along the way. We call this sequence “**data reduction cascade**”. The FLAMES-UVES pipeline involves the following input data:

- Raw frames:
  - Bias frames to determine a master bias.
  - Dark(s) (optional) frames to determine a master dark.
  - Slit flat frames to determine a set of master slit flats covering all the nine fibres.
  - Fibre formatcheck frames to determine guess order and line tables.
  - Calibration fibre flat frames to determine guess order tables and have a flat-field for that fibre.
  - Odd-even-all fibre order definition frames to determine the fibre order traces tables and construct a reference base of non contaminated fibre frames.
  - Fibre arc lamp frames to determine the line table for all fibres.
- calibration data products<sup>11</sup>:
  - A master bias and a master dark to subtract them from master flat, and science frames.

---

<sup>11</sup>These are either created by processing the above data, or they might be provided by another source, e.g. as part of a Service mode data package, in which case the corresponding RAW data listed above is not needed



<b>ESO</b>	<b>FLAMES-UVES Pipeline User Manual</b>	Doc:	VLT-MAN-ESO-19500-3016
		Issue:	Issue 5.0
		Date:	Date 2008-11-28
		Page:	52 of <a href="#">110</a>

- Master slit flat-field frames to correct for different detector pixel
  - Calibration fibre guess order tables to have a stable order tracing.
  - Calibration fibre order flat frames to have a flat frame to flat the calibration fibre and completed the base of fibre flats used to estimate the contamination of each fibre to the others.
  - Calibration line and order guess tables to have an automatic, accurate and stable wavelength calibration. efficiencies, the blaze function, the detector fringing at longer wavelengths.
  - Line tables to calibrate the object and sky spectra in wavelength.
- Reference files:
    - Line table to produce guess and final line tables.

### 9.3 Reduction cascade

Here we outline the logical sequence of steps which are needed to perform a complete spectral extraction of FLAMES-UVES data. All these operations can be carried out using the provided recipes. Examples of data reduction are in chapter 5.

1. Master bias and dark frames are generated.
2. For most frames (fibre and slit flat-field frames, science frames) a variance frame and a bad-pixel mask must be created.
3. A so-called format-check frame is examined. This frame must have been taken with a single fibre illuminated by a (Th-Ar) wavelength-calibration lamp, and all other fibres dark. This frame is compared with a physical model of the CCD illumination, with the help of an appropriate table of lines emitted by the calibration lamp, to derive a first guess of the order positions and of the wavelength-calibration solution.
4. The simultaneous calibration fibre flat frame is used to refine the guess order table and add an otherwise missing flat fibre trace to the pool of odd-even fibres to be used in the science frame extraction.
5. The whole set of fibres is considered (or at least those used in the science exposure), by processing both an odd-numbered- and an even-numbered fibre flat-field frames. In this way a complete (raw) fibre-order position table is obtained.
6. From the latter, the inter-order (background) table is computed, for later computation and subtraction of scattered-light contamination.
7. Slit flat-field frames are processed. At least two half-slit flat-field frames are needed to cover the fibre-illuminated region without overlapping, and these are combined and normalized. They will be used to remove pixel-to-pixel effects later. Associated variance frames and bad-pixel masks are also created.
8. odd/even fibre flat-field frames are processed. These frames contain the required information on the fibre cross-dispersion profile, and need to be input as separate odd/even fibre frames since adjacent fibres have partially overlapping profiles. Fibres for which there is no corresponding slit flat-field information are ignored. The frames are corrected for background light, and associated variance frames and bad-pixel masks are created.

<b>ESO</b>	<b>FLAMES-UVES Pipeline User Manual</b>	Doc:	VLT-MAN-ESO-19500-3016
		Issue:	Issue 5.0
		Date:	Date 2008-11-28
		Page:	53 of 110

9. The algorithm requires also an all-fibre flat-field frame, to compute the relative throughput among the odd and even fibres (that may have different intensity levels). It is treated as if it were a science frame, with the same extraction procedure, to obtain for each fibre and order a wavelength-dependent normalization (instead of a spectrum). Since the FLAMES-UVES data reduction includes two extraction methods (standard and optimal), the same method should be used to reduce the flat field and science frames if possible. Moreover, the all-fibre flat-field frame is also used as a reference for the positions of all fibres, since they are simultaneously lit here. For this purpose, the cross-dispersion shifts between the fibre positions in the odd and even fibre flat-field frames are computed with this common reference, and this additional information is stored in the order-fibre table.
10. Wavelength calibration is done, using a fibre frame with each fibre illuminated by a ThAr lamp, the guess line and fibre-order tables determined from previous data reduction steps. This step extracts each fibre spectrum, looks for ThAr lines and identifies them with help of the guess solution.
11. Eventually, the science spectra are extracted, either using standard or optimal extraction. Because of UVES stability limitations, there may be shifts between the fibre images in the flat-field frames and in the science frames. These shifts are expected to occur in the cross-dispersion direction and to be at most  $\pm 1$  pixel. Since especially optimal extraction is extremely sensitive to such shifts (spectra with large S/N cannot be extracted satisfactorily for shifts larger than 0.1 pix), the extraction routine computes the value of the shift (using a correlation-function method), and applies the opposite shift to the fibre flat-field frames (that is why pixel-to-pixel effects need to be computed separately using slit flat-field frames). After doing that, the optimal extraction is performed using the fibre flat-field frames as a model to fit the science frame at every wavelength (the fitted amplitudes are the extracted spectra at each wavelength), using at the same time the fibre profile information to deconvolve the partial fibre overlap. Standard extraction is somewhat simpler, doing only a sum of the science frame fibre spectra over a pre-defined window across dispersion, without fitting, but including also a deconvolution of fibre cross-contamination.
12. Next, extracted spectra are wavelength-calibrated.
13. Last, wavelength-calibrated, merged spectra are created.

### 9.3.1 Data reduction peculiarities

The data reduction software design to reduce FLAMES-UVES fibre-echelle mode data (described in [2], [3] and [5]) has some peculiarities when compared for example to the standard echelle data reduction, or with respect to the extraction of fibre spectra like those of FEROS.

In contrast to the standard echelle package:

- The fibre flat-field cross dispersion profiles (with pixel resolution) are used as "true" physical model of the fibre light distribution at each pixel on the detector in the optimal extraction of science targets, instead of some analytic approximation (In UVES-echelle mode, the corr-dispersion profile is measured using the science frame itself).
- A very good solution of the adjacent fibre contamination (which is a peculiar problem of fibre-fed multi-object spectrographs with close packed fibres and thus does not occur in echelle mode) is achieved.

<b>ESO</b>	<b>FLAMES-UVES Pipeline User Manual</b>	Doc:	VLT-MAN-ESO-19500-3016
		Issue:	Issue 5.0
		Date:	Date 2008-11-28
		Page:	54 of <a href="#">110</a>

- Spectra corresponding to different fibres, having a different fibre throughput, are corrected giving reliable relative fibre throughputs.
- Thanks to the use of three different kind of flat-field spectra (single fibre, odd-even-all fibres, slit flat-field frames) usually the final merged spectra have a very high quality, in which any residual oscillation artifacts are well below the noise level.

With respect to FEROS we have the same differences as before with an additional one:

- In FEROS the flat-fielding is done on the extracted spectra, first extracting the science and the flat-field and then dividing the first by the latter. Although this procedure has empirically been shown to lead to negligible differences with respect to the "correct" one, we still use the "correct" one (first correct science by flat-field, next extract the spectra).

The adopted data reduction together with the mentioned pros implies a number of cons. Essentially we have two:

- We need a lot of calibrations, compared for example to FEROS.
- The achieved high quality extraction quality is paid for with a quite extensive computational cost.



<b>ESO</b>	<b>FLAMES-UVES Pipeline User Manual</b>	Doc:	VLT-MAN-ESO-19500-3016
		Issue:	Issue 5.0
		Date:	Date 2008-11-28
		Page:	55 of 110

## 10 Pipeline Recipes Interfaces

In this section, for each recipe example, the required input data (and their classification tags) are provided. The column labeled “nr” indicates the number of required input frames with the following convention:

- 1 for a single frame (in case of pipeline calibration products, one for each detector chip of the corresponding arm),
- + for one or more input frames,
- ? for input frames which are optional (in case of pipeline calibration products one for the red1 and one for the redu detector chips).

In the following it is assumed that the input files (in our examples raw data are located in the directory /path\_raw, reference data are in the directory /path\_ref and pipeline products are in the directory /path\_pro) are existing FITS files (e.g. /data1/uves/com2/UVES.2004-08-16T02:54:04.353.fits, and /cal/uves/ech/cal/thargood\_3.fits, /cal/uves/ech/cal/linetable\_red1.fits, /cal/uves/ech/cal/linetable\_redu.fits).

A full list of pipeline products is also provided for each recipe, indicating their default recipe name (optionally replaced by EsoRex to a given standard), the value of the FITS keyword HIERARCH ESO PRO CATG (in short PRO.CATG) and a short description. The keyword PRO.CATG is used to classify each frame. The DPR and other relevant keywords are used to associate to each raw frame the proper calibration frame.

The data may be recognized and organized according to the values of the following FITS keywords:

Association keyword	Information
HIERARCH ESO DPR TYPE	raw data type
HIERARCH ESO DPR CATG	raw data category
HIERARCH ESO DPR TECH	raw data technique
HIERARCH ESO INS GRAT2 WLEN	Instrument setup central wavelength
HIERARCH ESO DET DIT	Integration time
HIERARCH ESO INS OBSPLATE	Observing plate
HIERARCH ESO OCS SIMCAL	Simultaneous calibration observation mode

The pipeline is able to also process pipeline products generated by the MIDAS based pipeline. The viceversa is not true.

For each recipe the input parameters (as they appear in the recipe configuration file), the corresponding parameter aliases (to be set on the command line) and their default values are listed. Each recipe has the following common parameters:

```
--debug           : Whether or not to save intermediate results to local
                    : directory. [FALSE]
--plotter         : Any plots produced by the recipe are redirected to the
                    : command specified by this parameter. The plotting
                    : command must contain the substring 'gnuplot' and must
```

<b>ESO</b>	<b>FLAMES-UVES Pipeline User Manual</b>	Doc:	VLT-MAN-ESO-19500-3016
		Issue:	Issue 5.0
		Date:	Date 2008-11-28
		Page:	56 of 110

be able to parse gnuplot syntax on its standard input. Valid examples of such a command may include 'gnuplot -persist' and 'cat > mygnuplot\$\$gp'. A finer control of the plotting options can be obtained by writing an executable script, e.g. my\_gnuplot.pl, that executes gnuplot after setting the desired gnuplot options (e.g. set terminal pslatex color). To turn off plotting, set this parameter to 'no'. [no]

--process\_chip : For RED arm data proces the redl, redu, or both chip(s). <both | redl | redu | REDL | REDU> [both]

A full description of each parameter is obtained by running the command **esorex -parameters**, or **esorex -help** or **esorex -man-page**, or by looking at the Recipe Input Parameters section of the dedicated Gasgano window. Also, the role of the most important parameters is described in section 11.

Also quality control parameters are computed. They are written to the header of corresponding pipeline products. More and up-to-date information on instrument quality control can be found on [www.eso.org/qc](http://www.eso.org/qc).

## 10.1 uves\_cal\_mbias

The recipe uves\_cal\_mbias computes a master bias frame.

### 10.1.1 Input

frame tag/category	nr	filename example
BIAS_ARM	+	/path_raw/uves_bias_arm.fits

### 10.1.2 Output

default recipe file name	format	PRO.CATG	short description
masterbias_chip.fits	2d (pix-pix) image	MASTER_BIAS_CHIP	master bias frame

### 10.1.3 Quality control

The recipe computes the following QC parameters:

QC.DUTYCYCL	Time to store a frame
QC.OUT1.RON.RAW	Read noise frame in ADU
QC.OUT1.RON.MASTER	Read noise frame in ADU
QC.OUT1.STRUCTY	Structure in Y
QC.OUT1.STRUCTX	Structure in X

<b>ESO</b>	<b>FLAMES-UVES Pipeline User Manual</b>	Doc:	VLT-MAN-ESO-19500-3016
		Issue:	Issue 5.0
		Date:	Date 2008-11-28
		Page:	57 of 110

#### 10.1.4 Parameters

There are no recipe specific input parameter for uves\_cal\_mbias.

```
--clean_traps      : Detector traps clean. If TRUE detector traps are
                    : corrected. The bad pixels are replaced by the average
                    : of nearest good pixels in the same column, or simply
                    : marked as bad. The positions of bad pixels are
                    : hard-coded (as function of UVES chip) their
                    : intensities need to be corrected. [TRUE]
```

### 10.2 uves\_cal\_mdark

The recipe uves\_cal\_mdark computes a master dark frame.

#### 10.2.1 Input

frame tag/category	nr	filename example
DARK_ARM	+	/path_raw/uves_dark_arm.fits
MASTER_BIAS_CHIP	1	/path_raw/masterbias_chip.fits

#### 10.2.2 Output

default recipe file name	format	PRO.CATG	short description
masterdark_chip.fits	2d (pix-pix) image	MASTER_DARK_CHIP	master dark frame

#### 10.2.3 Quality control

DFO trend the dark median level (PRO.DATAMED) for a given PRO.DATANCOM value and checks the values (min, max, median, mean, rms) of the dark on certain rectangular definable regions distributed across the detector.

QC.DATANCOM	Number of coadded frames
PRO.DATAMED	Median frame level
QC.REGij.MIN	Min of region i j of size <b>box_sx</b> × <b>box_sy</b>
QC.REGij.MAX	Max of region i j of size <b>box_sx</b> × <b>box_sy</b>
QC.REGij.AVG	Mean of region i j of size <b>box_sx</b> × <b>box_sy</b>
QC.REGij.MED	Median of region i j of size <b>box_sx</b> × <b>box_sy</b>
QC.REGij.RMS	Rms of region i j of size <b>box_sx</b> × <b>box_sy</b>
QC.REG.MIN.MIN	Min of all region Mins
QC.REG.MIN.MAX	Max of all region Mins
QC.REG.MIN.AVG	Mean of all region Mins
QC.REG.MIN.MED	Median of all region Mins
QC.REG.MIN.RMS	Rms of all region Mins
QC.REG.MAX.MIN	Min of all region Maxs

<b>ESO</b>	<b>FLAMES-UVES Pipeline User Manual</b>	Doc:	VLT-MAN-ESO-19500-3016
		Issue:	Issue 5.0
		Date:	Date 2008-11-28
		Page:	58 of 110

QC.REG.MAX.MAX	Max of all region Maxs
QC.REG.MAX.AVG	Mean of all region Maxs
QC.REG.MAX.MED	Median of all region Maxs
QC.REG.MAX.RMS	Rms of all region Maxs
QC.REG.AVG.MIN	Min of all region Means
QC.REG.AVG.MAX	Max of all region Means
QC.REG.AVG.AVG	Mean of all region Means
QC.REG.AVG.MED	Median of all region Means
QC.REG.AVG.RMS	Rms of all region Means
QC.REG.MED.MIN	Min of all region Medians
QC.REG.MED.MAX	Max of all region Medians
QC.REG.MED.AVG	Mean of all region Medians
QC.REG.MED.MED	Median of all region Medians
QC.REG.MED.RMS	Rms of all region Medians
QC.REG.RMS.MIN	Min of all region Rms
QC.REG.RMS.MAX	Max of all region Rms
QC.REG.RMS.AVG	Mean of all region Rms
QC.REG.RMS.MED	Median of all region Rms
QC.REG.RMS.RMS	Rms of all region Rms

#### 10.2.4 Parameters

```
--qc_dark.reg.num_x   : Number of regions on X direction (where mean/med/rms
                        : are computed). [4]
--qc_dark.reg.num_y   : Number of regions on Y direction (where mean/med/rms
                        : are computed). [4]
--qc_dark.reg.box_sx  : Region X size. [100]
--qc_dark.reg.box_sy  : Region Y size. [100]
--qc_dark.reg.border_x: X distance from the left hand side detector's border
                        : and the left hand side regin's bottom corner. [100]
--qc_dark.reg.border_y: Y distance from the left hand side detector's border
                        : and the left hand side regin's bottom corner. [100]
--qc_dark.reg.when    : When QC analysis is performed. 0: on each raw frame or
                        : 1: on the master frame. [0]
```

### 10.3 flames\_cal\_mkmaster

The recipe flames\_cal\_mkmaster computes the master slit flat frames.

#### 10.3.1 Input

frame tag/category	nr	filename example
SFLAT_RED	+	/path_raw/uves_flat_set1_red.fits
SFLAT_RED	+	/path_raw/uves_flat_set2_red.fits
SFLAT_RED	+	/path_raw/uves_flat_set3_red.fits
MASTER_BIAS_REDL	?	/path_pro/masterbias_redl.fits

<b>ESO</b>	<b>FLAMES-UVES Pipeline User Manual</b>	Doc:	VLT-MAN-ESO-19500-3016
		Issue:	Issue 5.0
		Date:	Date 2008-11-28
		Page:	59 of 110

```
| MASTER_BIAS_REDU | ? | /path_pro/masterbias_redu.fits |
```

### 10.3.2 Output

default recipe file name	format	PRO.CATG	short description
set1_masterflat_redl.fits	2d (pix-pix) image	MASTER_SFLAT_REDL	master flat frame
set1_masterflat_redu.fits	2d (pix-pix) image	MASTER_SFLAT_REDU	master flat frame
set2_masterflat_redl.fits	2d (pix-pix) image	MASTER_SFLAT_REDL	master flat frame
set2_masterflat_redu.fits	2d (pix-pix) image	MASTER_SFLAT_REDU	master flat frame
set3_masterflat_redl.fits	2d (pix-pix) image	MASTER_SFLAT_REDL	master flat frame
set3_masterflat_redu.fits	2d (pix-pix) image	MASTER_SFLAT_REDU	master flat frame

### 10.3.3 Quality control

The pipeline monitors the number of coadded frames (PRO.DATANCOM). Additional quality control information is monitored by DFO and can be found at [www.eso.org/qc](http://www.eso.org/qc).

### 10.3.4 Parameters

```
--backsub.mmethod      : Background measuring method. If equal to 'median' the
                        : background is sampled using the median of a
                        : subwindow. If 'minimum', the subwindow minimum value
                        : is used. If 'no', no background subtraction is done.
                        : <median | minimum | no> [median]
--backsub.npoints      : This is the number of columns in interorder space used
                        : to sample the background. [82]
--backsub.radiusy      : The height (in pixels) of the background sampling
                        : window is (2*radiusy + 1). This parameter is not
                        : corrected for binning. [2]
--backsub.sdegree      : Degree of interpolating splines. Currently only degree
                        : = 1 is supported. [1]
--backsub.smoothx      : If spline interpolation is used to measure the
                        : background, the x-radius of the post-smoothing window
                        : is (smoothx * image_width). Here, 'image_width' is
                        : the image width after binning. If negative, the
                        : default values are used: (25.0/4096) for blue
                        : flat-field frames, (50.0/4096) for red flat-field
                        : frames, (300.0/4096) for blue science frames and
                        : (300.0/4096) for red science frames. [-1.0]
--backsub.smoothy      : If spline interpolation is used to measure the
                        : background, the y-radius of the post-smoothing window
                        : is (smoothy * image_height). Here, 'image_height' is
                        : the image height after binning. If negative, the
                        : default values are used: (100.0/2048) for blue
                        : flat-field frames, (300.0/2048) for red flat-field
                        : frames, (200.0/2048) for blue science frames and
                        : (500.0/2048) for red science frames. [-1.0]
```

## 10.4 flames\_cal\_predict

The recipe flames\_cal\_predict computes the line and order guess tables using a model of UVES and the atmospheric pressure information written in the header of the frame, the known offset between the UVES slit's center

<b>ESO</b>	<b>FLAMES-UVES Pipeline User Manual</b>	Doc:	VLT-MAN-ESO-19500-3016
		Issue:	Issue 5.0
		Date:	Date 2008-11-28
		Page:	60 of 110

position and the one of the simultaneous calibration fibre, temperature and the instrument setting provided by the FITS header of the raw formatcheck frame which is acquired by illuminating the simultaneous fibre with a line emission lamp.

#### 10.4.1 Input

frame tag/category	nr	filename example
FIB_ARC_LAMP_FORM_RED	1	/path_raw/flames_uves_arc_lamp_form_red.fits
LINE_REFER_TABLE	1	/path_ref/thargood_3.fits
MASTER_BIAS_REDL	?	/path_pro/masterbias_redl.fits
MASTER_BIAS_REDU	?	/path_pro/masterbias_redu.fits

#### 10.4.2 Output

default recipe file name	format	PRO.CATG	short description
lineguesstable_redl.fits	table	FIB_LIN_GUE_REDL	Guess line table
lineguesstable_redu.fits	table	FIB_LIN_GUE_REDU	Guess line table
orderguesstable_redl.fits	table	FIB_ORD_REDL_CHIP	Guess order table
orderguesstable_redu.fits	table	FIB_ORD_REDU_CHIP	Guess order table

The guess line table contains the following columns:

X	Position along x
Y	Position along y
PEAK	line peak
Ident	line catalog wavelength
YNEW	Computed predicted line y position
Order	Relative order number
WAVEC	Predicted line wavelength of line peak
Aux	Product of wavelength and order number
XREG	Result of the polynomial regression $XREG=Aux(X)$
Pixel	Local dispersion
RORD	Order location (double precision to use it in a fit)
XPRED	Predicted X line position
YPRED	Predicted Y line position
XDIF	Difference between measured and predicted X line position
YDIF	Difference between measured and predicted Y line position
SELPLOT	selection column

The guess order table contains the following columns:

ABS_ORDER	Absolute order number
ORDER	Relative order number
X	Position along x

<b>ESO</b>	<b>FLAMES-UVES Pipeline User Manual</b>	Doc:	VLT-MAN-ESO-19500-3016
		Issue:	Issue 5.0
		Date:	Date 2008-11-28
		Page:	61 of 110

Y	Position along y
YFIT	Computed predicted order y position
RESIDUAL	Residual (Y-YFIT)

### 10.4.3 Quality control

The pipeline generates the following QC parameters:

QC.MODEL.NLINALL	Total number of detected lines
QC.MODEL.NLINSEL	Number of selected lines
QC.MODEL.DIFFXRMS	RMS difference of predicted and measured line x positions
QC.MODEL.DIFFXAVG	Mean difference of predicted and measured line x positions
QC.MODEL.DIFFXMED	Median difference of predicted and measured line x positions
QC.MODEL.DIFFYRMS	RMS difference of predicted and measured line y positions
QC.MODEL.DIFFYAVG	Mean difference of predicted and measured line y positions
QC.MODEL.DIFFYMED	Median difference of predicted and measured line y positions
QC.MODEL.RESXRMS	Std dev of X difference to physical model
QC.MODEL.RESXAVG	Average of X difference to physical model
QC.MODEL.RESXMED	Median of X difference to physical model
QC.MODEL.RESYRMS	Std dev of Y difference to physical model
QC.MODEL.RESYAVG	Average of Y difference to physical model
QC.MODEL.RESYMED	Median of Y difference to physical model
QC.MODEL.WLENMIN	Minimum predicted lines wavelength
QC.MODEL.WLENMAX	Maximum predicted lines wavelength
QC.MODEL.ORDMIN	Minimum predicted absolute order
QC.MODEL.ORDMAX	Maximum predicted absolute order
QC.WLENMIN	Minimum wavelength of spectral format
QC.WLENMAX	Maximum wavelength of spectral format
QC.ORDMIN	Minimum relative order
QC.ORDMAX	Maximum relative order
QC.FIB1.ABSTRANS	Average transmission countrate measured on the fibre
QC.FIB1.NHOTPIX	Number of found hot pixels and hot columns.
QC.FIB1.PLATENO	This is the Id of the plate to which the fibres were connected

### 10.4.4 Parameters

```
--mbox_x      : Match box X size. [40]
--mbox_y      : Match box Y size. [40]
--trans_x     : Detector translation along X. [0.0]
--trans_y     : Detector translation along Y. [0.0]
--ech_angle_off : Offset on echelle angle. [0.0]
--cd_angle_off : Offset on cross disperser angle. [0.0]
--ccd_rot_angle_off : Offset on CCD rotation angle. [0.0]
--compute_regression_sw : Compute regression?. [TRUE]
--x_axis_scale : Scale X axis. [0.0]
--y_axis_scale : Scale Y axis. [0.0]
--def_poll    : Polynomial X deg. [4]
```



<b>ESO</b>	<b>FLAMES-UVES Pipeline User Manual</b>	Doc:	VLT-MAN-ESO-19500-3016
		Issue:	Issue 5.0
		Date:	Date 2008-11-28
		Page:	62 of 110

```
--def_pol2      : Polynomial Y deg. [5]
--kappa         : Kappa value in kappa sigma clipping on RESIDUAL
                  between YFIT and Y columns. [4.5]
--tol          : Tolerance in kappa sigma clipping on RESIDUAL between
                  YFIT and Y columns. [2.0]
```

## 10.5 flames\_cal\_orderpos

The recipe `flames_cal_orderpos` generates an order table from a raw frame taken with the calibration fibre illuminated by a continuum lamp. This order table is more accurate than the order guess table previously generate. Additional product of this step is an order fibre flat frame, that covers the simultaneous calibration fibre, to be added to the pool of odd-even order flat frame to build a complete base of uncontaminate fibres frame to extract science data.

### 10.5.1 Input

frame tag/category	nr	filename example
FIB_ORDEF_RED	1	/path_raw/flames_uves_order_flat_red.fits
FIB_ORD_GUE_REDL	?	/path_pro/orderguesstable_redl.fits
FIB_ORD_GUE_REDU	?	/path_pro/orderguesstable_redu.fits
MASTER_BIAS_REDL	?	/path_pro/masterbias_redl.fits
MASTER_BIAS_REDU	?	/path_pro/masterbias_redu.fits

### 10.5.2 Output

default recipe file name	format	PRO.CATG	short description
ordertable_redl.fits	table	FIB_ORD_TAB_REDL	order table
ordertable_redu.fits	table	FIB_ORD_TAB_REDU	order table
order_def_redl.fits	table	FIB_ORDEF_REDL	fiber order def rotated frame
order_def_redu.fits	table	FIB_ORDEF_REDU	fiber order def rotated frame

The output table contains the columns

Order	Relative order number
X	Position along x
Y	Order line centroid location
dY	Uncertainty of Y
Residual_Square	Squared residual
OrderRMS	Root mean squared residual of initial one-dimensional linear fit of order
OrderSlope	Slope of order
Yfit	The fitted order location
dYfit_Square	Variance of Yfit
Residual	(Y - Yfit)

<b>ESO</b>	<b>FLAMES-UVES Pipeline User Manual</b>	Doc:	VLT-MAN-ESO-19500-3016
		Issue:	Issue 5.0
		Date:	Date 2008-11-28
		Page:	63 of 110

### 10.5.3 Quality control

QC.ORD.RESIDMIN	Min residuals in order table
QC.ORD.RESIDMAX	Max residuals in order table
QC.ORD.RESIDAVG	Mean residuals in order table
QC.ORD.RESIDRMS	RMS residuals in order table
QC.ORD.NPRED	Predicted number of orders
QC.ORD.NDET	Detected number of orders
QC.ORD.NPOSALL	Number of positions found
QC.ORD.NPOSSEL	Number of positions selected
QC.ORDMIN	Minimum (relative) order value
QC.ORDMAX	Maximum (relative) order value

Where the residuals measure the difference between the order solution obtained by applying the polynomial model and the corresponding order location measurements on the frame.

### 10.5.4 Parameters

```
--use_guess_tab      : If a Guess order table is provided this parameter
                        : defines how it is used:0: No usage, 1: use it to set
                        : lower/upper Y rows where order are searched 2: the
                        : order table tries to fully match the guess. [1]
--radx                : X-radius of median filtering window. [2]
--rady                : Y-radius of median filtering window. [1]
--mmethod              : Background subtraction method. If equal to 'median'
                        : the background is sampled using the median of a
                        : sub-window. If 'minimum', the minimum sub-window
                        : value is used. If 'no', no background subtraction is
                        : done. <median | minimum | no> [median]
--backsubgrid          : Number of grid points (in x- and y-direction) used to
                        : estimate the background (mode=poly). [50]
--backsubradiusy       : The height (in pixels) of the background sampling
                        : window is (2*radiusy + 1). This parameter is not
                        : corrected for binning. [2]
--backsubkappa         : The value of kappa in the one-sided kappa-sigma
                        : clipping used to estimate the background (mode=poly).
                        : [4.0]
--backsubdegx          : Degree (in x) of polynomial used to estimate the
                        : background (mode=poly). [2]
--backsubdegx          : Degree (in y) of polynomial used to estimate the
                        : background (mode=poly). [2]
--samplewidth          : Separation of sample traces (used by Hough transform)
                        : in input image. [50]
--minslope             : Minimum possible line slope. This should be the
                        : 'physical' slope on the chip, i.e. not taking binning
                        : factors into account, which is handled by the recipe.
                        : [0.0]
--maxslope             : Maximum possible line slope. [0.2]
--sloperes             : Resolution (width in pixels) of Hough space. [120]
--norders              : Number of echelle orders to detect. If set to 0 the
                        : predicted number of orders will be read from the
                        : guess order table. If no guess order table is given,
                        : the recipe will try to autodetect the number of
                        : orders. If the raw frame is a red chip and this
                        : parameter is specified, the given value will be used
```

<b>ESO</b>	<b>FLAMES-UVES Pipeline User Manual</b>	Doc:	VLT-MAN-ESO-19500-3016
		Issue:	Issue 5.0
		Date:	Date 2008-11-28
		Page:	64 of 110

```

                                for both red chips. [0]
--pthres                        : In automatic mode, or if the number of orders to
                                detect is read from a guess table, the detection of
                                new lines stops when the intensity of a candidate
                                line drops to less than 'pthres' times the intensity
                                of the previous detection. Otherwise - i.e. if the
                                number of orders to detect was specified by setting
                                the 'norders' parameters - this parameter is ignored.
                                [0.2]
--tracestep                    : The step size used when tracing the orders. [10]
--minthresh                    : The minimum threshold value is (min + minthresh*(max -
                                min)). Here 'min' and 'max' are the lowest and
                                highest pixel values in the central bin of the order.
                                [0.2]
--maxgap                       : If the order line drops below detection threshold, the
                                order tracing algorithm will try to jump a gap of
                                maximum size 'maxgap' multiplied by the image width.
                                [0.2]
--maxrms                       : When fitting the orders with straight lines, this is
                                the maximum allowed RMS relative to the median RMS of
                                all orders. [100.0]
--defpoll                      : The degree of the bivarite fit (cross dispersion
                                direction). If negative, the degree is optimized to
                                give the best fit. [-1]
--defpoll2                     : The degree of the bivarite fit (order number). If
                                negative, the degree is optimized to give the best
                                fit. [-1]
--kappa                       : Used for kappa-sigma clipping of the final polynomial
                                fit. If negative, no clipping is done. [6.0]

```

## 10.6 flames\_cal\_prep\_sff\_ofpos

The recipe `uves_cal_prep_sff_ofpos` is used to determine the fibre order table and construct several frames needed to extract a science fibre frame using tree input fibre frames obtained by illumination with a continuum source first the odd fibres, then the even fibres and, finally all the fibres.

### 10.6.1 Input

frame tag/category	nr	filename example
FIB_FF_ODD_RED	1	/path_raw/flames_uves_odd_red.fits
FIB_FF_EVEN_RED	1	/path_raw/flames_uves_even_red.fits
FIB_FF_ALL_RED	1	/path_raw/flames_uves_all_red.fits
MASTER_BIAS_REDL	?	/path_pro/masterbias_redl.fits
MASTER_BIAS_REDU	?	/path_pro/masterbias_redu.fits
FIB_ORD_TAB_REDL	1	/path_pro/ordertable_redl.fits
FIB_ORD_TAB_REDU	1	/path_pro/ordertable_redu.fits
MASTER_SFLAT_REDL	1	/path_pro/masterflat_set1_redl.fits
MASTER_SFLAT_REDL	1	/path_pro/masterflat_set2_redl.fits
MASTER_SFLAT_REDL	1	/path_pro/masterflat_set3_redl.fits
MASTER_SFLAT_REDU	1	/path_pro/masterflat_set1_redu.fits
MASTER_SFLAT_REDU	1	/path_pro/masterflat_set2_redu.fits
MASTER_SFLAT_REDU	1	/path_pro/masterflat_set3_redu.fits

<b>ESO</b>	<b>FLAMES-UVES Pipeline User Manual</b>	Doc:	VLT-MAN-ESO-19500-3016
		Issue:	Issue 5.0
		Date:	Date 2008-11-28
		Page:	65 of 110

Note: Simultaneous fibre order table frames used to be classified as FIB\_ORD\_GUE\_REDL/REDU. This tag is the same adopted for guess order tables generated by the `uves_cal_predict` and `flames_cal_predict` recipes. In order to remove this product classification name clash, and for backward compatibility (in order to support data reduction of input data generated by the MIDAS based pipeline) both tags are supported. The user is advised to use the new tag (FIB\_ORD\_TAB\_REDL/REDU), in case of data reduced with the CPL based pipeline.

## 10.6.2 Output

default recipe filename	format	PRO.CATG	short description
xt_odd_l	table	FIB_FF_ODD_INFO_TAB	info table
xt_even_l	table	FIB_FF_EVEN_INFO_TAB	info table
xt_all_l	table	FIB_FF_ODD_INFO_TAB	info table
slitff_common_redl	imagelist	SLIT_FF_COM_REDL	slitff common frame
slitff_norm_redl	image	SLIT_FF_NOR_REDL	slitff common frame
slitff_dtc_redl	imagelist	SLIT_FF_DTC_REDL	slitff data frame
slitff_sgc_redl	imagelist	SLIT_FF_SGC_REDL	slitff sigma frame
slitff_bpc_redl	imagelist	SLIT_FF_BPC_REDL	slitff badpixel frame
slitff_bnc_redl	imagelist	SLIT_FF_BNC_REDL	slitff frames boundaries
orfl	table	FIB_ORDEF_TABLE_REDL	fibre order table
fibreff_common_redl	imagelist	FIB_FF_COM_REDL	fibff common frame
fibreff_norm_redl	imagelist	FIB_FF_NOR_REDL	fibff common frame
fibreff_nsigma_redl	imagelist	FIB_FF_NSIG_REDL	Normalisation sigmas for fibre FF
fibreff_dtc_redl	imagelist	FIB_FF_DTC_REDL	fibff data frame
fibreff_sgc_redl	imagelist	FIB_FF_SGC_REDL	fibff sigma frame
fibreff_bpc_redl	imagelist	FIB_FF_BPC_REDL	fibff badpixel frame
slitff_common_redu	imagelist	SLIT_FF_COM_REDU	slitff common frame
slitff_norm_redu	image	SLIT_FF_NOR_REDU	slitff common frame
slitff_dtc_redu	imagelist	SLIT_FF_DTC_REDU	slitff data frame
slitff_sgc_redu	imagelist	SLIT_FF_SGC_REDU	slitff sigma frame
slitff_bpc_redu	imagelist	SLIT_FF_BPC_REDU	slitff badpixel frame
slitff_bnc_redu	imagelist	SLIT_FF_BNC_REDU	slitff frames boundaries
ordu	table	FIB_ORDEF_TABLE_REDU	fibre order table
fibreff_common_redu	imagelist	FIB_FF_COM_REDU	fibff common frame
fibreff_norm_redu	imagelist	FIB_FF_NOR_REDU	fibff common frame
fibreff_nsigma_redu	imagelist	FIB_FF_NSIG_REDU	Normalisation sigmas for fibre FF
fibreff_dtc_redu	imagelist	FIB_FF_DTC_REDU	fibff data frame
fibreff_sgc_redu	imagelist	FIB_FF_SGC_REDU	fibff sigma frame
fibreff_bpc_redu	imagelist	FIB_FF_BPC_REDU	fibff badpixel frame

## 10.6.3 Quality control

## 10.6.4 Parameters

--ext\_method : Extraction method. <std | opt | fst | fop | qst | qop>

<b>ESO</b>	<b>FLAMES-UVES Pipeline User Manual</b>	Doc:	VLT-MAN-ESO-19500-3016
		Issue:	Issue 5.0
		Date:	Date 2008-11-28
		Page:	66 of 110

```

[opt]
--bias_method      : Bias subtraction method, M for master bias frame, N
                    : for no bias subtraction, V to subtract a constant
                    : bias level defined by the parameter bias_value. <M |
                    : V | N> [M]
--bias_value       : Bias value (only if bias_method = V). [200]
--filter_switch    : Filter switch. <none | median> [none]
--sat_thr          : Saturation threshold. [55000]
--fileprep         : Slitff* and Fibreff* file preparation. If fast
                    : extraction method is used it should be set to FALSE.
                    : [TRUE]
--cubify           : Cubify switch. [TRUE]
--save_flat_size   : To be sure to use the flat part of a slit flatsone may
                    : need to subtract this bit. [-2]

```

## 10.7 flames\_cal\_wavecal

The recipe flames\_cal\_wavecal is used to determine the wavelength dispersion coefficients and construct a wavelength calibration table from a frame were all the fibres are illuminated from a arc line calibration lamp.

### 10.7.1 Input

frame tag/category	nr	filename example
FIB_ARC_LAMP_RED	1	/path_raw/flames_uves_arc_lamp_red.fits
MASTER_BIAS_REDL	?	/path_pro/masterbias_redl.fits
MASTER_BIAS_REDU	?	/path_pro/masterbias_redu.fits
FIB_ORDEF_TABLE_REDL	1	/path_pro/ordertable_redl.fits
FIB_ORDEF_TABLE_REDU	1	/path_pro/ordertable_redu.fits
FIB_LIN_GUE_REDL	1	/path_pro/lineguesstable_redl.fits
FIB_LIN_GUE_REDU	1	/path_pro/lineguesstable_redu.fits
LINE_REFER_TABLE	1	/path_ref/thargood_3.fits

### 10.7.2 Output

default recipe file name	format	PRO.CATG	short description
linetable_redl.fits	table	FIB_LINE_TABLE_REDL	line table
linetable_redu.fits	table	FIB_LINE_TABLE_REDU	line table

The output line table(s), LINE\_TABLE\_REDL, LINE\_TABLE\_REDU contains three extensions per fibre (i.e. 6x3=18 or 8x3=24 extensions). The first extensions contain a table with the columns described below

X	Horizontal position (from Gaussian fit) of detected line
dX	Uncertainty (one sigma) of X
XWidth	Width (in pixels) of detected line (from Gaussian fit)
Y	Relative order number of detected line
Peak	Intensity of detected line
Background	Fitted background (ADU) of detected line

<b>ESO</b>	<b>FLAMES-UVES Pipeline User Manual</b>	Doc:	VLT-MAN-ESO-19500-3016
		Issue:	Issue 5.0
		Date:	Date 2008-11-28
		Page:	67 of 110

Slope	Linear background slope (ADU/pixel) of detected line
Order	Relative order number of detected line
AbsOrder	Absolute order number of detected line
Ynew	Vertical position of detected line
WaveC	Wavelength of this line (computed using the resulting dispersion relation)
dLambdaC	Uncertainty (one sigma) of 'WaveC'.
Pixel	The local dispersion in A/pixel.
Residual	Residual (in A) of this line
Residual_pix	Residual (in pixels) of this line
Lambda_candidate	Nearest line in catalogue
dLambda_cat_sq	Squared distance to nearest catalogue line
dLambda_nn_sq	Squared distance to nearest neighbour multiplied by ALPHA
Ident	The wavelength associated with this emission line, or NULL if this line was not identified.
dIdent	Uncertainty of catalogue wavelength
Select	1 if the line was identified, 0 otherwise
NLinSol	1 if the line was identified and accepted for the polynomial fit, 0 otherwise
Intensity	Intensity of detected line scaled to unit exposure time. (This column is present only if a LINE_INTMON_TABLE is provided.)

The 2nd table extension contains the dispersion relation on the form  $p(x, m) = \lambda \cdot m$ , where  $m$  is the order number. The 3rd table extension contains the map from (pixel, pixel)-space to physical order numbers (used internally by the calibration recipe, a 2d polynomial on the form  $p(x, y) = m$ ).

If there is more than one extraction window (default is 3), the results of each calibration is stored in subsequent table extensions of the same FITS file. For example, extensions 4, 5 and 6 would contain the resulting line table (and its two associated polynomials) for the second extraction window. The results for the calibration of the  $n$ 'th extraction window is stored in extensions  $(3*n - 2)$  to  $3*n$ .

The offset of the extraction window is stored in the FITS keyword like "HISTORY OFFSET -7.5". The corresponding window number (counting from 1) is stored in keywords like "HISTORY WINDOW 3"

The polynomials are stored in table extensions as in the example

Order1	Order2	Coeff
-1	-1	a0
-1	-1	a1
-1	-1	a2
-1	-1	b0
-1	-1	b1
-1	-1	b2
0	0	c00
0	1	c01
0	2	c02

<b>ESO</b>	<b>FLAMES-UVES Pipeline User Manual</b>	Doc:	VLT-MAN-ESO-19500-3016
		Issue:	Issue 5.0
		Date:	Date 2008-11-28
		Page:	68 of 110

1	0	c10
1	1	c11
1	2	c12
2	0	c20
2	1	c21
2	2	c22

The third column contains the polynomial coefficients corresponding to the degree defined by the two first columns. The six first table rows defines a linear transformation of the dependent and independent variables. For example the table shown represents the polynomial  $p$  defined by

$$(p(x, y) - a_0)/b_0 = q((x - a_1)/b_1, (y - a_2)/b_2)$$

and

$$q(x, y) = \sum_{i=0}^2 \sum_{j=0}^2 c_{ij} x^i y^j$$

The linear transformation of the three variables was introduced in order to ensure numerical stability in the polynomial fitting routine.

### 10.7.3 Quality control

The recipe computes the following quality control parameters:

QC.FWHMAVG	Average FWHM of lines selected
QC.FWHMRMS	Standard deviation of FWHM of selected lines
QC.FWHMMED	Median FWHM of selected lines
QC.RESOLAVG	Average resolving power of selected lines
QC.RESOLRMS	Standard deviation of the resolving power of selected lines
QC.RESOLMED	Median resolving power of selected lines
QC.LINE.RESIDAVG	Mean of residuals of line positions to fit
QC.LINE.RESIDRMS	Sigma of residuals of line positions to fit
QC.WLENMIN	Minimum wavelength of detected order
QC.WLENMAX	Maximum wavelength of detected order
QC.ORDMIN	Minimum order number detected
QC.ORDMAX	Maximum order number detected
QC.NLINTOT	Total number of lines found on the frame
QC.NLINSEL	Number of selected lines
QC.NLINRES1	Number of lines with residuals < 0.1 nm
INS.SLIT3.WID	Slit width
INS.GRAT2.WLEN	Grating central wavelength [nm].
INS.TEMP2.MEAN	Average temperature [C].



<b>ESO</b>	<b>FLAMES-UVES Pipeline User Manual</b>	Doc:	VLT-MAN-ESO-19500-3016
		Issue:	Issue 5.0
		Date:	Date 2008-11-28
		Page:	69 of 110

QC.NHOTPIX	Number of pixels outside the range [-20,55000]
QC.PLATENO	Plate setting value
QC.FIB8.DRSNO	Sequential number of detected fibre
QC.FIB8.SEQ	Sequential number of detected fibre
QC.FIB8.POS	Detected fibre Y position relative to order slit center
QC.FIB8.MSK	Detected fibre mask value (1 if the fibre was recognised as lit)
QC.FIB8.FWHMAVG	Average FWHM of lines detected for the given fibre
QC.FIB8.FWHMRMS	RMS of lines detected for the given fibre
QC.FIB8.FWHMMED	Median FWHM of lines detected for the given fibre
QC.FIB8.RESOLAVG	Average instrumental resolution for the given fibre
QC.FIB8.RESOLRMS	RMS of instrumental resolution for the given fibre
QC.FIB8.RESOLMED	Median instrumental resolution for the given fibre
QC.FIB8.LINE.RESIDAVG	Average line residuals to reference values for the given fibre
QC.FIB8.LINE.RESIDRMS	RMS of line residuals to reference values for the given fibre
QC.FIB8.WLENMIN	Minimum wavelength [nm] of detected lines for the given fibre
QC.FIB8.WLENMAX	Maximum wavelength [nm] of detected lines for the given fibre
QC.FIB8.ORDMIN	Minimum order value for the given fibre
QC.FIB8.ORDMAX	Maximum order value for the given fibre
QC.FIB8.NLINTOT	Total number of detected lines for the given fibre
QC.FIB8.NLINSEL	Number of selected lines for the given fibre
QC.FIB8.NLINSOL	Number of lines used for the final wavelength solution for the given fibre
QC.FIB8.NLINRES1	Number of lines with residuals less than 1 Å for the given fibre

The line FWHMs and the corresponding resolving power are measured along the dispersion direction (see also Figure 11.2.2). The residuals measures the differences between the solution from the polynomial fit model and the corresponding line positions found on the arc lamp frame (see also Figure 11.2.1).

#### 10.7.4 Parameters

--nwindows	: Number of extraction windows per trace. The windows will be aligned (i.e. no overlap and no spacing between adjacent windows). Unless an offset is specified, the middle window(s) is centered on the trace. [1]
--length	: Length (in pixels) of each extraction window. This parameter is also equal to the separation of adjacent window centers, causing the extraction windows to always be aligned. The parameter is automatically adjusted according to the binning of the input raw frame. If negative, the extraction window length is determined automatically to cover the full slit. [7.0]
--offset	: A global offset (in pixels) of all extraction windows. [0.0]
--extract.method	: Extraction method. <average   linear   2d   weighted   optimal> [average]
--extract.kappa	: In optimal extraction mode, this is the threshold for bad (i.e. hot/cold) pixel rejection. If a pixel deviates more than kappa*sigma (where sigma is the uncertainty of the pixel flux) from the inferred spatial profile, its weight is set to zero. If this

<b>ESO</b>	<b>FLAMES-UVES Pipeline User Manual</b>	Doc:	VLT-MAN-ESO-19500-3016
		Issue:	Issue 5.0
		Date:	Date 2008-11-28
		Page:	70 of 110

```

parameter is negative, no rejection is performed.
[10.0]
--extract.chunk      : In optimal extraction mode, the chunk size (in pixels)
                        used for fitting the analytical profile (a fit of the
                        analytical profile to single bins would suffer from
                        low statistics). [32]
--extract.profile    : In optimal extraction mode, the kind of profile to
                        use. 'gauss' gives a Gaussian profile, 'moffat' gives
                        a Moffat profile with beta=4 and a possible linear
                        sky contribution. 'virtual' uses a virtual resampling
                        algorithm (i.e. measures and uses the actual object
                        profile). 'constant' assumes a constant spatial
                        profile and allows optimal extraction of wavelength
                        calibration frames. 'auto' will automatically select
                        the best method based on the estimated S/N of the
                        object. For low S/N, 'moffat' or 'gauss' are
                        recommended (for robustness). For high S/N, 'virtual'
                        is recommended (for accuracy). In the case of virtual
                        resampling, a precise determination of the order
                        positions is required; therefore the order-definition
                        is repeated using the (assumed non-low S/N) science
                        frame. <constant | gauss | moffat | virtual | auto>
                        [auto]
--extract.skymethod  : In optimal extraction mode, the sky subtraction method
                        to use. 'median' estimates the sky as the median of
                        pixels along the slit (ignoring pixels close to the
                        object), whereas 'optimal' does a chi square
                        minimization along the slit to obtain the best
                        combined object and sky levels. The optimal method
                        gives the most accurate sky determination but is also
                        a bit slower than the median method. <median |
                        optimal> [optimal]
--extract.oversample : The oversampling factor used for the virtual
                        resampling algorithm. If negative, the value 5 is
                        used for S/N <=200, and the value 10 is used if the
                        estimated S/N is > 200. [-1]
--extract.best       : (optimal extraction only) If false (fastest), the
                        spectrum is extracted only once. If true (best), the
                        spectrum is extracted twice, the second time using
                        improved variance estimates based on the first
                        iteration. Better variance estimates slightly improve
                        the obtained signal to noise but at the cost of
                        increased execution time. [TRUE]
--range              : Width (in pixels) of search window is 2*range + 1.
                        This parameter is automatically adjusted according to
                        binning. [8]
--minlines           : Minimum number of lines to detect. If zero, the
                        default value (2000 for BLUE/REDL chips; 1000 for
                        REDU chip) is used. [0]
--maxlines           : Maximum number of lines to detect. If zero, the
                        default value (2400 for BLUE/REDL chip; 1400 for REDU
                        chip) is used. [0]
--centeringmethod    : Line centering method. <gaussian | gravity> [gaussian]
--shiftmax           : The maximum shift (in pixels) compared to guess
                        solution. This parameter is automatically corrected
                        for binning. [10.0]
--shiftstep          : The step size (in pixels) used when searching for the
                        optimum shift. This parameter is automatically
                        corrected for binning. [0.1]
--shiftttoler        : Tolerance (in pixels) when matching shifted lines.
                        This parameter is not adjusted according to binning.
                        [0.05]
--alpha              : The parameter that controls the distance to the
                        nearest neighbours. [0.1]
--maxerror           : This parameter controls the graceful exit of the

```

<b>ESO</b>	<b>FLAMES-UVES Pipeline User Manual</b>	Doc:	VLT-MAN-ESO-19500-3016
		Issue:	Issue 5.0
		Date:	Date 2008-11-28
		Page:	71 of 110

```

identification loop. If the RMS of the global fit
exceeds this value (in pixels) the iteration stops.
[20.0]
--degree      : Degrees of the global 2d dispersion polynomial. If a
                negative number is specified, the polynomial degrees
                are automatically selected by starting from (1, 1)
                and increasing the degrees as long as the RMS residual
                decreases significantly. [5]
--tolerance    : Tolerance of fit. If positive, the tolerance is in
                pixel units. If negative, abs(tolerance) is in
                wavelength units. Lines with residuals worse than the
                tolerance are excluded from the final fit. Unlike in
                previous versions, this parameter is not corrected
                for CCD binning. This rejection based on the absolute
                residual in pixel can be effectively disabled by
                setting the tolerance to a very large number (e.g.
                9999). In that case outliers will be rejected using
                only kappa sigma clipping. [0.6]
--kappa        : Lines with residuals more then kappa stdev are
                rejected from the final fit. [4.0]
--rebin.wavestep : The bin size (in w.l.u.) in wavelength space. If
                negative, a step size of 2/3 * ( average pixel size )
                is used. [-1.0]
--rebin.scale   : Whether or not to multiply by the factor dx/dlambda
                (pixels per wavelength) during the rebinning. This
                option is disabled as default in concordance with the
                method used in the MIDAS pipeline. This option should
                be set to true to convert the observed flux (in
                pixel-space) to a flux per wavelength (in
                wavelength-space). [FALSE]

```

## 10.8 flames\_obs\_scired

This recipe reduces a science fiber frame.

### 10.8.1 Input

frame tag/category	nr	filename example
FIB_SCI_RED	1	/path_raw/flames_uves_science_red.fits
SLIT_FF_COM_REDL	1	/path_pro/slitff_com_redl.fits
SLIT_FF_COM_REDU	1	/path_pro/slitff_com_redu.fits
SLIT_FF_NOR_REDL	1	/path_pro/slitff_nor_redl.fits
SLIT_FF_NOR_REDU	1	/path_pro/slitff_nor_redu.fits
SLIT_FF_DTC_REDL	1	/path_pro/slitff_dtc_redl.fits
SLIT_FF_DTC_REDU	1	/path_pro/slitff_dtc_redu.fits
SLIT_FF_BPC_REDL	1	/path_pro/slitff_bpc_redl.fits
SLIT_FF_BPC_REDU	1	/path_pro/slitff_bpc_redu.fits
SLIT_FF_SGC_REDL	1	/path_pro/slitff_sgc_redl.fits
SLIT_FF_SGC_REDU	1	/path_pro/slitff_sgc_redu.fits
FIB_FF_COM_REDL	1	/path_pro/fibreff_com_redl.fits
FIB_FF_COM_REDU	1	/path_pro/fibreff_com_redu.fits
FIB_FF_NOR_REDL	1	/path_pro/fibreff_nor_redl.fits
FIB_FF_NOR_REDU	1	/path_pro/fibreff_nor_redu.fits

<b>ESO</b>	<b>FLAMES-UVES Pipeline User Manual</b>	Doc:	VLT-MAN-ESO-19500-3016
		Issue:	Issue 5.0
		Date:	Date 2008-11-28
		Page:	72 of 110

FIB_FF_DTC_REDL	1	/path_pro/fibreff_dtc_redl.fits
FIB_FF_DTC_REDU	1	/path_pro/fibreff_dtc_redu.fits
FIB_FF_BPC_REDL	1	/path_pro/fibreff_bpc_redl.fits
FIB_FF_BPC_REDU	1	/path_pro/fibreff_bpc_redu.fits
FIB_FF_SGC_REDL	1	/path_pro/fibreff_sgc_redl.fits
FIB_FF_SGC_REDU	1	/path_pro/fibreff_sgc_redu.fits
FIB_FF_BNC_REDL	1	/path_pro/fibreff_bnc_redl.fits
FIB_FF_BNC_REDU	1	/path_pro/fibreff_bnc_redu.fits
FIB_ORDEF_TABLE_REDL	1	/path_pro/ordertable_redl.fits
FIB_ORDEF_TABLE_REDU	1	/path_pro/ordertable_redu.fits
FIB_LINE_TABLE_REDL	1	/path_pro/linetable_redl.fits
FIB_LINE_TABLE_REDU	1	/path_pro/linetable_redu.fits
MASTER_BIAS_REDL	?	/path_pro/masterbias_redl.fits
MASTER_BIAS_REDU	?	/path_pro/masterbias_redu.fits

## 10.8.2 Output

default recipe file name	format	PRO.CATG	short description
bin_table_info_redl	table	FIB_SCI_INFO_TAB_REDL	table with information on fibres setup
fxb_l_raw000i	2d (pix-ord) image	XB_SCI_RAW_REDL	extracted, flatfielded raw frame
wfxb_l_raw000i	2d (wav-ord) image	WXB_SCI_RAW_REDL	rebinned, extracted, flatfielded raw frame
mwfxb_l_raw000i	1d (wav) image	MWXB_SCI_RAW_REDL	merged,rebinned, extracted, flatfielded raw frame
fxb_l_rawsig000i	2d (pix-ord) image	ERR_XB_SCI_RAW_REDL	error frame
wfxb_l_rawsig000i	2d (wav-ord) image	ERR_WXB_SCI_RAW_REDL	error frame
mwfxb_l_raw000i_sigma	1d (wav) image	ERR_MWXB_SCI_RAW_REDL	error merged, rebinned, extracted, flat fielded raw frame
fxb_l_000i	2d (pix-ord) image	XB_SCI_REDL	extracted, flatfielded raw frame
wfxb_l_000i	2d (wav-ord) image	WXB_SCI_REDL	rebinned, extracted, flatfielded raw frame
mwfxb_l_000i	1d (wav) image	MWXB_SCI_REDL	merged, rebinned, extracted, flatfielded raw frame
fxb_l_sig000i	2d (pix-ord) image	ERR_XB_SCI_REDL	error frame
wfxb_l_sig000i	2d (wav-ord) image	ERR_WXB_SCI_REDL	error frame
mwfxb_l_000i_sigma	1d (wav) image	ERR_MWXB_SCI_REDL	merged, rebinned, extracted, flat fielded raw frame
fxb_u_raw000i	2d (pix-ord) image	XB_SCI_RAW_REDU	extracted, flatfielded raw frame
wfxb_u_raw000i	2d (wav-ord) image	WXB_SCI_RAW_REDU	rebinned, extracted, flatfielded raw frame
mwfxb_u_raw000i	1d (wav) image	MWXB_SCI_RAW_REDU	merged,rebinned, extracted, flatfielded raw frame
fxb_u_rawsig000i	2d (pix-ord) image	ERR_XB_SCI_RAW_REDU	error frame
wfxb_u_rawsig000i	2d (wav-ord) image	ERR_WXB_SCI_RAW_REDU	error frame
mwfxb_u_raw000i_sigma	1d (wav) image	ERR_MWXB_SCI_RAW_REDU	error merged, rebinned, extracted, flat fielded raw frame
fxb_u_000i	2d (pix-ord) image	XB_SCI_REDU	extracted, flatfielded raw frame
wfxb_u_000i	2d (wav-ord) image	WXB_SCI_REDU	rebinned, extracted, flatfielded raw frame
mwfxb_u_000i	1d (wav) image	MWXB_SCI_REDU	merged,rebinned, extracted, flatfielded raw frame
fxb_u_sig000i	2d (pix-ord) image	ERR_XB_SCI_REDU	error frame
wfxb_u_sig000i	2d (wav-ord) image	ERR_WXB_SCI_REDU	error frame
mwfxb_u_000i_sigma	1d (wav) image	ERR_MWXB_SCI_REDU	error merged, rebinned, extracted, flat fielded raw frame

## 10.8.3 Quality control

The recipe computes the folowing quality control parameters in the FITS header of the reduced frame.

## 10.8.4 Quality control

The recipe computes the following quality control parameters:

ESO	FLAMES-UVES Pipeline User Manual	Doc:	VLT-MAN-ESO-19500-3016
		Issue:	Issue 5.0
		Date:	Date 2008-11-28
		Page:	73 of 110

QC.FWHMAVG	Average FWHM of lines selected
------------	--------------------------------

### 10.8.5 Parameters

```

--ext_method      : Extraction method. <std | opt | fst | fop | qst | qop>
                   [opt]
--cor_max_fnd     : Find correlation maximum?. <N | Y> [Y]
--cor_def_rng     : Correlation range size?. [6.0]
--cor_def_pnt     : Correlation sampling points?. [25]
--cor_def_off     : Correlation center offset?. [0.0]
--bias_method     : Bias subtraction method. <M | V | N> [M]
--bias_value      : Bias value (only if bias_method = V). [200]
--cubify_sw       : Cubify switch. <Y | N> [N]
--filt_sw         : Filter switch. <none | median> [none]
--bkg_max_io_win  : Background window number in each full inter order.
                   [500]
--bkg_xy_win_sz_x : x maximum size of each background window:. [6]
--bkg_xy_win_sz_y : y maximum size of each background window:. [2]
--pixel_thresh_max : Pixel saturation threshold max. [55000]
--pixel_thresh_min : Pixel saturation threshold min. [-20]
--input_fmt_cube  : Input data format. [TRUE]
--output_fmt_cube : Output data format. [FALSE]
--drs_k_s_thre    : Kappa sigma threshold. [10.0]
--drs_base_name    : Base name for science products. [fxb]
--drs_maxyshift    : Half width of the interval to scan for correlation,
                   when determining y shift. [3.0]
--drs_ext_w_siz   : Integration window size good: 10 (if fibre
                   deconvolution works fine). [10.0]
--rebin.wavestep  : The bin size (in w.l.u.) in wavelength space. If
                   negative, a step size of 2/3 * ( average pixel size )
                   is used. [-1.0]
--rebin.scale     : Whether or not to multiply by the factor dx/dlambda
                   (pixels per wavelength) during the rebinning. This
                   option is disabled as default in concordance with the
                   method used in the MIDAS pipeline. This option should
                   be set to true to convert the observed flux (in
                   pixel-space) to a flux per wavelength (in
                   wavelength-space). [FALSE]
--merge          : Order merging method. If 'optimal', the flux in the
                   overlapping region is set to the (optimally computed,
                   using the uncertainties) average of single order
                   spectra. If 'sum', the flux in the overlapping region
                   is computed as the sum of the single order spectra.
                   If flat-fielding is done, method 'optimal' is
                   recommended, otherwise 'sum'. <optimal | sum>
                   [optimal]

```

## 10.9 flames\_obs\_redchain

### 10.9.1 Input

frame tag/category	nr	filename example
BIAS_RED	+	/path_raw/uves_bias_red.fits
DARK_RED	+	/path_raw/uves_dark_red.fits
MASTER_SFLAT_CHIP	1	/path_pro/set1_master_flat_chip.fits
MASTER_SFLAT_CHIP	1	/path_pro/set2_master_flat_chip.fits

ESO	FLAMES-UVES Pipeline User Manual	Doc:	VLT-MAN-ESO-19500-3016
		Issue:	Issue 5.0
		Date:	Date 2008-11-28
		Page:	74 of 110

MASTER_SFLAT_CHIP	1	/path_pro/set3_master_flat_chip.fits
FIB_ARC_FORM_RED	1	/path_raw/flames_uves_arc_form_red.fits
FIB_ORDER_FLAT_RED	1	/path_raw/flames_uves_oflat_red.fits
FIB_FF_ODD_RED	1	/path_raw/flames_uves_odd_red.fits
FIB_FF_EVEN_RED	1	/path_raw/flames_uves_even_red.fits
FIB_FF_ALL_RED	1	/path_raw/flames_uves_all_red.fits
FIB_ARC_LAMP_RED	1	/path_raw/flames_uves_arc_lamp_red.fits
FIB_SCI_RED	1	/path_raw/flames_uves_science_red.fits
LINE_REFER_TABLE	1	/path_ref/thargood_3.fits

Note: in present release flames\_obs\_redchain cannot directly handle raw SFLAT frames.

## 10.9.2 Output

This recipe generates all the products described for the previous recipes.

## 10.9.3 Quality control

This recipe generates all the quality control parameters described for the previous recipes.

## 10.9.4 Parameters

This recipe accepts all the same parameters as the individual recipes except with the recipe name prefixing the parameter command line option, e.g. - uves\_cal\_mbias.clean\_traps to set the clean\_traps parameter for the mbias part of the recipe. The additional parameter **scired** (default value is TRUE) switches on/off the execution of the last step (science data reduction).

```
--scired          : Whether or not to do science reduction. If false, only
                    master calibration frames are created. If false,
                    either zero or all necessary calibration frames must
                    be provided for each arm. [TRUE]
--uves_cal_mbias.clean_traps : Detector traps clean. If TRUE detector traps
                    are corrected. The bad pixels are replaced by
                    the average of nearest good pixels in the
                    same column, or simply marked as bad. The
                    positions of bad pixels are hard-coded (as
                    function of UVES chip). [TRUE]
--uves_cal_mdark.qc_dark.reg.num_x :
                    Number of regions on X direction (where mean/med/rms
                    are computed). [4]
--uves_cal_mdark.qc_dark.reg.num_y :
                    Number of regions on Y direction (where mean/med/rms
                    are computed). [4]
--uves_cal_mdark.qc_dark.reg.box_sx :
                    Region X size. [100]
--uves_cal_mdark.qc_dark.reg.box_sy :
                    Region Y size. [100]
--uves_cal_mdark.qc_dark.reg.border_x :
                    X distance from the left hand side detector's border
```

<b>ESO</b>	<b>FLAMES-UVES Pipeline User Manual</b>	<b>Doc:</b>	<b>VLT-MAN-ESO-19500-3016</b>
		<b>Issue:</b>	<b>Issue 5.0</b>
		<b>Date:</b>	<b>Date 2008-11-28</b>
		<b>Page:</b>	<b>75 of 110</b>

```

        and the left hand side regin's bottom corner. [100]
--uves_cal_mdark.qc_dark.reg.border_y :
        Y distance from the left hand side detector's border
        and the left hand side regin's bottom corner. [100]
--uves_cal_mdark.qc_dark.reg.when : When QC analysis is performed. 0: on each
        raw frame or 1: on the master frame. [0]
--flames_cal_predict.mbox_x : Match box X size. [40]
--flames_cal_predict.mbox_y : Match box Y size. [40]
--flames_cal_predict.trans_x : Detector translation along X. [0.0]
--flames_cal_predict.trans_y : Detector translation along Y. [0.0]
--flames_cal_predict.ech_angle_off :
        Offset on echelle angle. [0.0]
--flames_cal_predict.cd_angle_off : Offset on cross disperser angle. [0.0]
--flames_cal_predict.ccd_rot_angle_off :
        Offset on CCD rotation angle. [0.0]
--flames_cal_predict.compute_regression_sw :
        Compute regression?. [TRUE]
--flames_cal_predict.x_axis_scale : Scale X axis. [0.0]
--flames_cal_predict.y_axis_scale : Scale Y axis. [0.0]
--flames_cal_predict.def_pol1 : Polynomial X deg. [4]
--flames_cal_predict.def_pol2 : Polynomial Y deg. [5]
--flames_cal_predict.kappa : Kappa value in kappa sigma clipping on RESIDUAL
        between YFIT and Y columns. [4.5]
--flames_cal_predict.tol : Tolerance in kappa sigma clipping on RESIDUAL
        between YFIT and Y columns. [2.0]
--flames_cal_orderpos.use_guess_tab :
        If a Guess order table is provided this parameter
        defines how it is used:0: No usage, 1: use it to set
        lower/upper Y rows where order are searched 2: the
        order table tries to fully match the guess. [1]
--flames_cal_orderpos.radx : X-radius of median filtering window. [2]
--flames_cal_orderpos.rady : Y-radius of median filtering window. [1]
--flames_cal_orderpos.mmmethod : Background subtraction method. If equal to
        'median' the background is sampled using the
        median of a sub-window. If 'minimum', the
        minimum sub-window value is used. If 'no',
        no background subtraction is done. <median |
        minimum | no> [median]
--flames_cal_orderpos.backsubgrid : Number of grid points (in x- and
        y-direction) used to estimate the
        background (mode=poly). [50]
--flames_cal_orderpos.backsubradiusy :
        The height (in pixels) of the background sampling
        window is (2*radiusy + 1). This parameter is not
        corrected for binning. [2]
--flames_cal_orderpos.backsubkappa :
        The value of kappa in the one-sided kappa-sigma
        clipping used to estimate the background (mode=poly).
        [4.0]
--flames_cal_orderpos.backsubdegx : Degree (in x) of polynomial used to
        estimate the background (mode=poly). [2]
--flames_cal_orderpos.backsubdegy : Degree (in y) of polynomial used to
        estimate the background (mode=poly). [2]
--flames_cal_orderpos.samplewidth : Separation of sample traces (used by
        Hough transform) in input image. [50]
--flames_cal_orderpos.minslope : Minimum possible line slope. This should be
        the 'physical' slope on the chip, i.e. not
        taking binning factors into account, which
        is handled by the recipe. [0.0]
--flames_cal_orderpos.maxslope : Maximum possible line slope. [0.2]
--flames_cal_orderpos.sloperes : Resolution (width in pixels) of Hough space.
        [120]
--flames_cal_orderpos.norders : Number of echelle orders to detect. If set to
        0 the predicted number of orders will be
        read from the guess order table. If no guess

```



<b>ESO</b>	<b>FLAMES-UVES Pipeline User Manual</b>	Doc:	VLT-MAN-ESO-19500-3016
		Issue:	Issue 5.0
		Date:	Date 2008-11-28
		Page:	76 of 110

```

order table is given, the recipe will try to
autodetect the number of orders. If the raw
frame is a red chip and this parameter is
specified, the given value will be used for
both red chips. [0]
--flames_cal_orderpos.pthres : In automatic mode, or if the number of orders
to detect is read from a guess table, the
detection of new lines stops when the
intensity of a candidate line drops to less
than 'pthres' times the intensity of the
previous detection. Otherwise - i.e. if the
number of orders to detect was specified by
setting the 'norders' parameters - this
parameter is ignored. [0.2]
--flames_cal_orderpos.tracestep : The step size used when tracing the orders.
[10]
--flames_cal_orderpos.minthresh : The minimum threshold value is (min +
minthres*(max - min)). Here 'min' and
'max' are the lowest and highest pixel
values in the central bin of the order.
[0.2]
--flames_cal_orderpos.maxgap : If the order line drops below detection
threshold, the order tracing algorithm will
try to jump a gap of maximum size 'maxgap'
multiplied by the image width. [0.2]
--flames_cal_orderpos.maxrms : When fitting the orders with straight lines,
this is the maximum allowed RMS relative to
the median RMS of all orders. [100.0]
--flames_cal_orderpos.defpol1 : The degree of the bivarite fit (cross
dispersion direction). If negative, the
degree is optimized to give the best fit.
[-1]
--flames_cal_orderpos.defpol2 : The degree of the bivarite fit (order
number). If negative, the degree is
optimized to give the best fit. [-1]
--flames_cal_orderpos.kappa : Used for kappa-sigma clipping of the final
polynomial fit. If negative, no clipping is
done. [6.0]
--uves_cal_mflat.backsub.mmethod : Background measuring method. If equal to
'median' the background is sampled using
the median of a subwindow. If 'minimum',
the subwindow minimum value is used. If
'no', no background subtraction is done.
<median | minimum | no> [median]
--uves_cal_mflat.backsub.npoints : This is the number of columns in
interorder space used to sample the
background. [82]
--uves_cal_mflat.backsub.radiusy : The height (in pixels) of the background
sampling window is (2*radiusy + 1). This
parameter is not corrected for binning.
[2]
--uves_cal_mflat.backsub.sdegree : Degree of interpolating splines. Currently
only degree = 1 is supported. [1]
--uves_cal_mflat.backsub.smoothx : If spline interpolation is used to measure
the background, the x-radius of the
post-smoothing window is (smoothx *
image_width). Here, 'image_width' is the
image width after binning. If negative,
the default values are used: (25.0/4096)
for blue flat-field frames, (50.0/4096)
for red flat-field frames, (300.0/4096)
for blue science frames and (300.0/4096)
for red science frames. [-1.0]
--uves_cal_mflat.backsub.smoothy : If spline interpolation is used to measure
the background, the y-radius of the

```

ESO	FLAMES-UVES Pipeline User Manual	Doc:	VLT-MAN-ESO-19500-3016
		Issue:	Issue 5.0
		Date:	Date 2008-11-28
		Page:	77 of 110

```

        post-smoothing window is (smoother *
        image_height). Here, 'image_height' is
        the image height after binning. If
        negative, the default values are used:
        (100.0/2048) for blue flat-field frames,
        (300.0/2048) for red flat-field frames,
        (200.0/2048) for blue science frames and
        (500.0/2048) for red science frames.
        [-1.0]

--flames_cal_wavecal.nwindows : Number of extraction windows per trace. The
        windows will be aligned (i.e. no overlap and
        no spacing between adjacent windows). Unless
        an offset is specified, the middle window(s)
        is centered on the trace. [1]

--flames_cal_wavecal.length : Length (in pixels) of each extraction window.
        This parameter is also equal to the separation
        of adjacent window centers, causing the
        extraction windows to always be aligned. The
        parameter is automatically adjusted according
        to the binning of the input raw frame. If
        negative, the extraction window length is
        determined automatically to cover the full
        slit. [7.0]

--flames_cal_wavecal.offset : A global offset (in pixels) of all extraction
        windows. [0.0]

--flames_cal_wavecal.extract.method :
        Extraction method. <average | linear | 2d | weighted |
        optimal> [average]

--flames_cal_wavecal.extract.kappa :
        In optimal extraction mode, this is the threshold for
        bad (i.e. hot/cold) pixel rejection. If a pixel
        deviates more than kappa*sigma (where sigma is the
        uncertainty of the pixel flux) from the inferred
        spatial profile, its weight is set to zero. If this
        parameter is negative, no rejection is performed.
        [10.0]

--flames_cal_wavecal.extract.chunk :
        In optimal extraction mode, the chunk size (in pixels)
        used for fitting the analytical profile (a fit of the
        analytical profile to single bins would suffer from
        low statistics). [32]

--flames_cal_wavecal.extract.profile :
        In optimal extraction mode, the kind of profile to
        use. 'gauss' gives a Gaussian profile, 'moffat' gives
        a Moffat profile with beta=4 and a possible linear
        sky contribution. 'virtual' uses a virtual resampling
        algorithm (i.e. measures and uses the actual object
        profile). 'constant' assumes a constant spatial
        profile and allows optimal extraction of wavelength
        calibration frames. 'auto' will automatically select
        the best method based on the estimated S/N of the
        object. For low S/N, 'moffat' or 'gauss' are
        recommended (for robustness). For high S/N, 'virtual'
        is recommended (for accuracy). In the case of virtual
        resampling, a precise determination of the order
        positions is required; therefore the order-definition
        is repeated using the (assumed non-low S/N) science
        frame. <constant | gauss | moffat | virtual | auto>
        [auto]

--flames_cal_wavecal.extract.skymethod :
        In optimal extraction mode, the sky subtraction method
        to use. 'median' estimates the sky as the median of
        pixels along the slit (ignoring pixels close to the
        object), whereas 'optimal' does a chi square
        minimization along the slit to obtain the best

```

<b>ESO</b>	<b>FLAMES-UVES Pipeline User Manual</b>	<b>Doc:</b>	<b>VLT-MAN-ESO-19500-3016</b>
		<b>Issue:</b>	<b>Issue 5.0</b>
		<b>Date:</b>	<b>Date 2008-11-28</b>
		<b>Page:</b>	<b>78 of 110</b>

```

combined object and sky levels. The optimal method
gives the most accurate sky determination but is also
a bit slower than the median method. <median |
optimal> [optimal]
--flames_cal_wavecals.extract.oversample :
The oversampling factor used for the virtual
resampling algorithm. If negative, the value 5 is
used for S/N <=200, and the value 10 is used if the
estimated S/N is > 200. [-1]
--flames_cal_wavecals.extract.best : (optimal extraction only) If false
(fastest), the spectrum is extracted
only once. If true (best), the spectrum
is extracted twice, the second time
using improved variance estimates based
on the first iteration. Better variance
estimates slightly improve the obtained
signal to noise but at the cost of
increased execution time. [TRUE]
--flames_cal_wavecals.range : Width (in pixels) of search window is 2*range +
1. This parameter is automatically adjusted
according to binning. [8]
--flames_cal_wavecals.minlines : Minimum number of lines to detect. If zero,
the default value (2000 for BLUE/REDL chips;
1000 for REDU chip) is used. [0]
--flames_cal_wavecals.maxlines : Maximum number of lines to detect. If zero,
the default value (2400 for BLUE/REDL chip;
1400 for REDU chip) is used. [0]
--flames_cal_wavecals.centeringmethod :
Line centering method. <gaussian | gravity> [gaussian]
--flames_cal_wavecals.shiftmax : The maximum shift (in pixels) compared to
guess solution. This parameter is
automatically corrected for binning. [10.0]
--flames_cal_wavecals.shiftstep : The step size (in pixels) used when
searching for the optimum shift. This
parameter is automatically corrected for
binning. [0.1]
--flames_cal_wavecals.shifftoler : Tolerance (in pixels) when matching shifted
lines. This parameter is not adjusted
according to binning. [0.05]
--flames_cal_wavecals.alpha : The parameter that controls the distance to the
nearest neighbours. [0.1]
--flames_cal_wavecals.maxerror : This parameter controls the graceful exit of
the identification loop. If the RMS of the
global fit exceeds this value (in pixels)
the iteration stops. [20.0]
--flames_cal_wavecals.degree : Degrees of the global 2d dispersion polynomial.
If a negative number is specified, the
polynomial degrees are automatically selected
by starting from (1, 1) and increasing the
degrees as long as the RMS residual decreases
significantly. [4]
--flames_cal_wavecals.tolerance : Tolerance of fit. If positive, the tolerance
is in pixel units. If negative,
abs(tolerance) is in wavelength units.
Lines with residuals worse than the
tolerance are excluded from the final fit.
Unlike in previous versions, this parameter
is not corrected for CCD binning. This
rejection based on the absolute residual in
pixel can be effectively disabled by
setting the tolerance to a very large
number (e.g. 9999). In that case outliers
will be rejected using only kappa sigma
clipping. [0.6]
--flames_cal_wavecals.kappa : Lines with residuals more than kappa stdev are

```

<b>ESO</b>	<b>FLAMES-UVES Pipeline User Manual</b>	Doc:	VLT-MAN-ESO-19500-3016
		Issue:	Issue 5.0
		Date:	Date 2008-11-28
		Page:	79 of 110

```

rejected from the final fit. [4.0]
--flames_cal_wavecal.rebin.wavestep :
    The bin size (in w.l.u.) in wavelength space. If
    negative, a step size of 2/3 * ( average pixel size )
    is used. [-1.0]
--flames_cal_wavecal.rebin.scale : Whether or not to multiply by the factor
    dx/dlambda (pixels per wavelength) during
    the rebinning. This option is disabled as
    default in concordance with the method
    used in the MIDAS pipeline. This option
    should be set to true to convert the
    observed flux (in pixel-space) to a flux
    per wavelength (in wavelength-space).
    [FALSE]
--flames_cal_prep_sff_ofpos.ext_method :
    Extraction method. <std | opt | fst | fop | qst | qop>
    [opt]
--flames_cal_prep_sff_ofpos.bias_method :
    Bias subtraction method, M for master bias frame, N
    for no bias subtraction, V to subtract a constant
    bias level defined by the parameter bias_value. <M |
    V | N> [M]
--flames_cal_prep_sff_ofpos.bias_value :
    Bias value (only if bias_method = V). [200]
--flames_cal_prep_sff_ofpos.filter_switch :
    Filter switch. <none | median> [none]
--flames_cal_prep_sff_ofpos.sat_thr :
    Saturation threshold. [55000]
--flames_cal_prep_sff_ofpos.fileprep :
    Slitff* and Fibreff* file preparation. If fast
    extraction method is used it should be set to FALSE.
    [TRUE]
--flames_cal_prep_sff_ofpos.cubify :
    Cubify switch. [TRUE]
--flames_cal_prep_sff_ofpos.save_flat_size :
    To be sure to use the flat part of a slit flatsone may
    need to subtract this bit. [-2]
--flames_obs_scired.ext_method : Extraction method. <std | opt | fst | fop |
    qst | qop> [opt]
--flames_obs_scired.cor_max_fnd : Find correlation maximum?. <N | Y> [Y]
--flames_obs_scired.cor_def_rng : Correlation range size?. [6.0]
--flames_obs_scired.cor_def_pnt : Correlation sampling points?. [25]
--flames_obs_scired.cor_def_off : Correlation center offset?. [0.0]
--flames_obs_scired.bias_method : Bias subtraction method. <M | V | N> [M]
--flames_obs_scired.bias_value : Bias value (only if bias_method = V). [200]
--flames_obs_scired.cubify_sw : Cubify switch. <Y | N> [N]
--flames_obs_scired.filt_sw : Filter switch. <none | median> [none]
--flames_obs_scired.bkg_max_io_win :
    Background window number in each full inter order.
    [500]
--flames_obs_scired.bkg_xy_win_sz_x :
    x maximum size of each background window:. [6]
--flames_obs_scired.bkg_xy_win_sz_y :
    y maximum size of each background window:. [2]
--flames_obs_scired.pixel_thresh_max :
    Pixel saturation threshold max. [55000]
--flames_obs_scired.pixel_thresh_min :
    Pixel saturation threshold min. [-20]
--flames_obs_scired.input_fmt_cube :
    Input data format. [TRUE]
--flames_obs_scired.output_fmt_cube :
    Output data format. [FALSE]
--flames_obs_scired.drs_k_s_thre : Kappa sigma threshold. [10.0]
--flames_obs_scired.drs_base_name : Base name for science products. [fxb]
--flames_obs_scired.drs_maxyshift : Half width of the interval to scan for

```

<b>ESO</b>	<b>FLAMES-UVES Pipeline User Manual</b>	Doc:	VLT-MAN-ESO-19500-3016
		Issue:	Issue 5.0
		Date:	Date 2008-11-28
		Page:	80 of <a href="#">110</a>

```

correlation, when determining y shift.
[3.0]
--flames_obs_scired.drs_ext_w_siz : Integration window size good: 10 (if
fibre deconvolution works fine). [10.0]
--flames_obs_scired.rebin.wavestep :
The bin size (in w.l.u.) in wavelength space. If
negative, a step size of 2/3 * ( average pixel size )
is used. [-1.0]
--flames_obs_scired.rebin.scale : Whether or not to multiply by the factor
dx/dlambda (pixels per wavelength) during
the rebinning. This option is disabled as
default in concordance with the method
used in the MIDAS pipeline. This option
should be set to true to convert the
observed flux (in pixel-space) to a flux
per wavelength (in wavelength-space).
[FALSE]
--flames_obs_scired.merge : Order merging method. If 'optimal', the flux in
the overlapping region is set to the (optimally
computed, using the uncertainties) average of
single order spectra. If 'sum', the flux in the
overlapping region is computed as the sum of the
single order spectra. If flat-fielding is done,
method 'optimal' is recommended, otherwise
'sum'. <optimal | sum> [optimal]

```

<b>ESO</b>	<b>FLAMES-UVES Pipeline User Manual</b>	Doc:	VLT-MAN-ESO-19500-3016
		Issue:	Issue 5.0
		Date:	Date 2008-11-28
		Page:	81 of 110

## 11 Algorithms and recipe details

In this section we describe the main algorithms implemented in the FLAMES-UVES pipeline recipes. Relevant data reduction parameters are typed in **bold** face. Parameters that are meant to be critical for data reduction quality and recommended not to be changed unless the user really knows what is the correspondent effect, are in upper case. They are also written in the FITS header of relevant products as FITS keywords. In this release they are hardcoded. This was also true in the MIDAS based release. The chosen values demonstrated in years of operations to provide a good accuracy and robustness. They are defined in `flames_def_drs_par.c` module and the user may change them at own risk. Those are listed in Appendix B.

### 11.1 Algorithms

#### 11.1.1 Error model

The pipeline does full error propagation using the error propagation formula and making the usual assumption about Gaussian error bars. The error of the  $i$ 'th pixel of any input raw frame is defined as

$$\sigma_i = \sqrt{g(C_{i,ADU} - b_i) + ron^2 + (g^2 - 1)/12}, \quad (1)$$

where  $C_{i,ADU}$  is the flux in analog-to-digital units,  $b_i$  is the master bias level and  $g$  and  $ron$  are the detector gain and read-out noise as defined in the FITS header.

#### 11.1.2 Physical model offset setting

On the basis of the plate number, wavelength setting, CCD chip, appropriate values are set to describe the physical model offsets between the simultaneous calibration fibre and the UVES slit center. This allows to use the same physical model implementation as the one describing the echelle mode spectral format.

#### 11.1.3 Fibre frames preparation

To each input fibre frame are associated an error and a bad pixel frame. The descriptor indicating the maximum number of fibres and which fibre is on/off are set.

This function takes as input the same odd/even fibre flat-field frames which are used for the fibre-order tracing (see 11.2.6), normalizes them (each fibre is normalized to 1 integrated flux at each x position), and selects only fibres falling within the region illuminated in the slit flat-field frames, to create a new set of odd/even fibre flat-field frames plus other auxilliary frames (normalization, fibre boundaries) used by subsequent DRS commands. Any bad pixels in the input fibre flat-fields are filled with interpolated values, if possible. Each fibre is then correlated with a Gaussian centered at the position traced by the order-fibre table, with a variable offset.

The variable y offset is allowed to span the intervall between **-MAXYSHIFT** and **+MAXYSHIFT**. The half width at half maximum of the Gaussian pseudo-fibres is read from the **GUASSFIBRESIGMA** keyword, and their total half width is read from the **GAUSSHALFWIDTH** keyword. The offset corresponding to the maximum correlation is saved in a keyword (**YOFFSET**), to be used as a “zero-point” correction for fibre asymmetry

<b>ESO</b>	<b>FLAMES-UVES Pipeline User Manual</b>	Doc:	VLT-MAN-ESO-19500-3016
		Issue:	Issue 5.0
		Date:	Date 2008-11-28
		Page:	82 of 110

when measuring actual fibre offsets on the frames. For this correction to be effective, subsequent correlations are forced to use the same Gaussian parameters (i.e. **GAUSSFIBRESIGMA** and **GAUSSHALFWIDTH**) which were used here.

#### 11.1.4 Slit frames preparation

To each input slit frame are associated an error and a bad pixel map frames.

This function also orders the set of input slit flats by Y position and rescale them in intensity so that each frame intensity is the same as the one of the previous adjacent one where the two overlap, and creates a new, minimal set of slit flat-field frames plus other auxiliary frames (normalization, boundary data) used by subsequent data reduction commands. On each input slit flat frame, the **HALFWIDTH** and **YSHIFT** FITS keywords are set, defining respectively the half width size of the flat part of the orders and the offset between the order centres in this frame and the order centres as traced by the order-fibre table. See also 11.1.13. Of course, for relative equalisation to be possible, each frame *must* have some overlap with the neighbouring ones in the set.

#### 11.1.5 Lit fibres flag setting

This algorithm fills the **MAXFIBRES**, **FIBREMASK** descriptors in the input frame, describing how many fibres are illuminated, and their fibre numbers. This step is executed as soon as the input fibre frames are re-oriented. The descriptors values are assigned according to the frame classification obtained from DPR keywords (ESO DPR TECH, ESO DPR CATG, ESO DPR TYPE). For a proper automatic setting of the descriptors, which is vital for data reduction, it is mandatory that the DPR keyword set is as specified in Table 5.3. Otherwise it is very probable that data reduction (in particular the correlation step) fails.

#### 11.1.6 Inter-order background determination

This command takes as input the order table to create a table of positions where the scattered-light background may be computed, in regions not illuminated by any fibre. The number of windows in each inter-order space is set by keyword **BKG\_MAX\_IO\_WIND** (set by default to 500). The x,y widths of each window are set by keyword **BKG\_XY\_WIN\_SZ** (6 pixels in x, and 2 pixels in y).

- First, the order centres in the input frame are computed as the middle point between the first and last lit fibres (as read from the **FIBREMASK** FITS keyword). If this is offset with respect to the order positions mapped by the order table, the positions of the windows in the background table are corrected accordingly. Then, all the windows in the background table are trimmed to exclude any position in the input frame where light from the fibres is expected to be present. Windows which happen to be trimmed to zero size are obviously discarded at this stage.
- Step two: if the **BKGBADSCAN** keyword is set to “fraction” or “absolute”, a neighbourhood of size determined by the **BKGBADWIN** keyword of each pixel in each window left from the previous step is checked for bad pixels. Bad pixels in each neighbourhood are counted, and if they exceed the fraction set in **BKGBADMAXFRAC** (if **BKGBADSCAN**=“fraction”) or the total number set in **BKGBADMAXTOT** (if **BKGBADSCAN**=“absolute”), that pixel is removed from the window. Windows which turn



<b>ESO</b>	<b>FLAMES-UVES Pipeline User Manual</b>	Doc:	VLT-MAN-ESO-19500-3016
		Issue:	Issue 5.0
		Date:	Date 2008-11-28
		Page:	83 of 110

out to be empty after this additional trimming are discarded. This second trimming step is skipped if **BKGBADSCAN** is set to “none”.

- Step three: a list of points to be fit is built from the input frame and the remaining windows from the steps above, in one out of four possible ways, selected by the value of the **BKGFITMETHOD** keyword. If **BKGFITMETHOD**=“all”, the list consists simply of all the pixels contained in the windows. If **BKGFITMETHOD**=“median”, for each window the pixel with the median value is found and added to the list. If **BKGFITMETHOD**=“minimum”, for each window the pixel with the minimum value is found and added to the list. Finally, if **BKGFITMETHOD**=“average”, for each window the median value is found first, all pixels in the window differing from it more than  $\kappa_{\text{sigma\_threshold}}$  times their standard deviation are discarded, and the weighted average of the remaining pixels in the window is added to the list.

This data reduction step generates the fitted background image.

- As a fourth step, a bivariate polynomial is actually fit to the list of pixels.
- Then, the values in the list are compared to the fitted values: if any pixel in the list deviates from the fitted values more than  $\kappa_{\text{sigma\_threshold}}$  (**DRS\_K\_S\_THRE**) times the standard deviation, at most **MAXDISCARDFRACT** times the number of pixels in the list are discarded in one pass, in order of worst fit, and the fit repeated on the remaining ones. This loop is repeated until no more pixels are rejected from the list, or until it has been repeated **MAXBACKITERS** times, whichever comes first.

The scattered light background is estimated in a FLAMES-UVES frame, fitting a bivariate polynomial on the regions of the input frame listed in the background table. The degrees of the polynomial are read from the **BKGPOL** keyword.

- The resulting bivariate polynomial is then evaluated over the full input frame, and saved.

### 11.1.7 Fibre shift determination

This algorithm measures the shift of each fibre of an input frame with respect to the corresponding fibre of the odd-even base, using a set of normalised fibre flat-field frames, as produced by the fibre frame preparation (11.1.15) and updated by fibre frames normalization, (11.1.16), and using a set of equalised slit flat-field frames, as produced by the slit frames preparation (11.1.4) and an order-fibre table, as produced by `flames_cal_prep_ofpos` (11.2.6) and updated by the fibre frame normalization 11.1.16). The **MAXYSHIFT** keyword sets the value of the maximum  $y$  shift (in world coordinates) to be searched between the input frame and the flat-field frames. The default value of the maximum allowed  $y$  shift is set by the keyword **MAXYSHIFT**.

The measured **YSHIFT** is written into the input frame’s FITS header, and measures the  $y$  shift values of each fibre of the input frame with respect to the fibres in each normalised flat-field frame.

The descriptor **NFLATS** written into the input frame containing the size of the **YSHIFT** descriptor.

The  $y$ -axis shift of the orders/fibres in the input frame is found by correlating the frame with a synthetic one composed of Gaussian-shaped pseudo-fibres centred on the positions traced by the order-fibre table, with a variable offset in the direction perpendicular to dispersion.

The size of the Gaussian shapes is forced to be the same as the one which was used to measure the “zero-point” shift measured at the fibre frames preparation (11.1.15) stage.

<b>ESO</b>	<b>FLAMES-UVES Pipeline User Manual</b>	Doc:	VLT-MAN-ESO-19500-3016
		Issue:	Issue 5.0
		Date:	Date 2008-11-28
		Page:	84 of 110

The correlation is computed using one pixel every **CORRELXSTEP**; thus **CORRELXSTEP**=1 means that all pixels are used, **CORRELXSTEP**=2 means that every other pixels is used and so on. The safest bet (and the default) is **CORRELXSTEP**=1. Increasing **CORRELXSTEP** cuts the correlation times significantly, but **CORRELXSTEP** values larger than the size of the resolution element are risky: on a line emission spectrum it may happen that only the dark parts of the frame are used for the correlation, resulting in numerical instability and useless results.

The variable  $y$  offset is allowed to span the interval between **-MAXYSHIFT** and **+MAXYSHIFT**, and is referred to the fibre positions as measured on an all fibres flat-field frame by the Fiber frames normalization (see [11.1.16](#) command), and stored in the order table.

Note that during the science data reduction, due to the fibre illumination pattern peculiar to the given observation, the cross correlation step may not lead to a single maximum. At this purpose we allow the user to control manually that step, through parameters **cor\_max\_fnd**, **cor\_def\_rng**, **cor\_def\_pnt**, **cor\_def\_off**. The user should check the cross correlation profile in the file **cor\_shape\_suffix.fits** (where suffix='l' or 'u'), and then set proper values of **cor\_def\_off** (centered to the expected correlation maximum) and **cor\_def\_rng** (allowing to scan correlation values near the expected maximum) and **cor\_def\_pnt** (this set the accuracy in the determination of the correlation maximum position).

The maximum of the correlation is determined using a Brent search with accuracy better than **CORRELTOL** and the corresponding  $y$  offset corrected for the “zero-point” shift measured at the fibre frames preparation (see [11.1.15](#)) stage.

This final  $y$  offset is saved as a descriptor called **YSHIFT** in the input frame. Since it is referred to a unique reference position, it is duplicated as many times as the number of normalised fibre flat-field frames.

### 11.1.8 Generation of UVES related fibre information

This algorithm extracts from the table extension of an input raw frame the fibre information corresponding to the UVES fibre link.

This auxiliary command is used to extract important information relative to FLAMES-UVES observations contained in the binary table extensions in the output table **out\_tab** (what filename??).

### 11.1.9 Frames preparation

This data reduction step prepares frames read from disk for usage by the following steps. Its input may be a single frame or a set of frames. For each indicated frame, two new frames are created, a variance frame and a bad-pixel mask frame, respectively with suffixes **\_sigma.fits** and **\_mask.fits**. Their names are written into the FITS keywords **SIGMAFRAME** and **BADPXFRAME** of the input frame.

The variance and bad pixel frames are updated during the following data reduction steps.

Moreover, two additional methods to detect and flag bad pixels are provided: recursive median filtering and clipping of values out of a given validity interval.

<b>ESO</b>	<b>FLAMES-UVES Pipeline User Manual</b>	Doc:	VLT-MAN-ESO-19500-3016
		Issue:	Issue 5.0
		Date:	Date 2008-11-28
		Page:	85 of 110

**Recursive median filtering** If recursive median filtering is enabled setting the parameter **filter\_switch** (of the recipe `flames_cal_prep_sff_ofpos`) to `median`, for each pixel the median of the good pixels in a neighbourhood defined by the **DRS\_FILT\_HW\_X** and **DRS\_FILT\_HW\_Y** keywords is computed and if the pixel value differs by more than **DRS\_K\_S\_THRE** times the computed sigma, the pixel is flagged as bad. This procedure is iterated until no new bad pixels are detected. While this procedure is very effective at pinpointing bad pixels, it is also computationally intensive, and thus very time consuming. For this reason it is disabled by default (**filter\_switch**=none).

**Threshold filtering** If the **DRS\_FILT\_SAT\_SW** is set to YES, each pixel is compared with the validity interval provided by the **DRS\_PTHRE\_MIN** and **DRS\_PTHRE\_MAX** parameters and any pixels outside this interval are flagged as bad. This is meant to quickly catch saturated pixels and/or bad values due to electronics. The cuts must be chosen keeping in mind whether the input frame was already bias-subtracted or not, and allows for normal noise, in order not to cause a spurious mass rejection of pixels.

The final bad pixel mask will thus be the merged set of all bad pixels detected by any means, or previously marked as such in a user-defined bad pixel frame whose name(s) was (were) set in the **BADPXFRAME** keyword (in a FITS file called like `'*_mask.fits'`, for example `b_sci_prefix_mask.fits`, or `b_odd_prefix_mask.fits`, or `b_even_prefix_mask.fits`, or `b_all_prefix_mask.fits` etc..., where prefix='l' or 'u' for the lower or upper detector chip).

#### 11.1.10 Bad pixel map preparation

The location of the traps and dead columns on the CCD image are mapped according to the instrument setting (wavelength, bin, CCD chip), as determined by examining some reference image. Those positions are hardwired in the code. In addition, all pixels above the saturation threshold set by the keyword **SAT\_THRESHOLD** are marked as bad. The user may change this value modifying the parameter **sat\_thr** of `flames_cal_prep_sff_ofpos`, and **pixel\_thresh\_max** in `flames_obs_scired`. With this information a bad-pixel image is created, that will be used in the following data reduction to mask appropriately the bad pixels.

#### 11.1.11 Fibre order tracing

This data reduction step finds all “traces” in a set of fibre flat-field frames, in which lit fibres do not overlap, typically an even and an odd fibres frame. The **NBTRACES** keyword, if set to a value different from 0, defines the number of fibres to be detected.

The traces are found using the Hough transform technique. To improve robustness of the wavelength calibration and of the extraction the user may set with the FITS keywords **DRS\_SCAN\_MIN** and **DRS\_SCAN\_MAX**. Those keywords are defined by three values, **DRS\_SCAN\_MIN**=55,73,73 and **DRS\_SCAN\_MAX**=1993,1975,1975, respectively for the 520, 580 and 860 wavelength setting, and are used to specify the y extension of detector region where the fibre traces are searched.

These traces are then correlated with the order positions from the first guess order table, and labelled according to order and fibre number as deduced from the **FIBREMASK** keyword in the fibre flat field frames.

The orders are then fit with a polynomial of degrees as defined in the **DEFPOL** keyword, whose component are defined, during the single fibre order tracing, assuming constant relative positions of the fibres. The polynomial

<b>ESO</b>	<b>FLAMES-UVES Pipeline User Manual</b>	Doc:	VLT-MAN-ESO-19500-3016
		Issue:	Issue 5.0
		Date:	Date 2008-11-28
		Page:	86 of 110

coefficients are stored in the **COEFFD** FITS keyword of the output order-fibre table, while the relative position of the fibres are stored in each component of the **FIBREPOS** FITS keyword of the same table.

#### 11.1.12 Iter order background determination

To determine the interorder background the information about the fibre location contained in the input order table is used to create a table of positions where the scattered-light background may be computed, in regions not illuminated by any fibre.

#### 11.1.13 Measure of flat part size of Slit frames

The pipeline finds the **FWHM** across  $y$  of the illuminated part of a slit flat-field frame near frame center, finds the offset ( $y$  shift) of the order centres in this frame with respect to order positions as traced by the order table, and writes these data as keywords **YSHIFT** and **HALFWIDTH** of the input slit flat-field frame.

A cross-order trace cut is determined at the image X center. This profile is smoothed (radius=5pix) and normalized to its maximum. Then the average master flat size is determined. To this value is subtracted the value **save\_flat\_size**. This value is critical to the final accuracy. We set this parameter to -2 (thus we actually slightly increase the size of the flat as computed on the frame). The user may use also values in the range [-3,+1] in case of not satisfying quality. The value of **HALFWIDTH** is equal to half of the computed flat width. Then is determined the flat offset with respect to ther order table trace. To obtain this value we build a synthetic flat frame of flat size equal to the one previously computed and centered on the order trace. Then cross correlating the synthetic frame and each slit flat frame we determine the corresponding shift to the order trace and store it in the **YSHIFT** keyword.

#### 11.1.14 Slit frames preparation

This algorithm uses a set of calibration slit flat-field frames (at least 2 frames), orders them by Y shift and rescales them in intensity so that each frame intensity is the same as the one of the previous adjacent frame where the two overlap, and creates a new, minimal set of slit flat-field frames plus other auxiliary frames (containing information related to frame normalization, and slit flats boundaries) that is later on used during data reduction. On each input slit flat field frame, the **HALFWIDTH** and **YSHIFT** FITS keyword should be set as described in 11.1.13. Those define respectively the half width of the flat part of the orders and the offset between the order centres of the corresponding frame and the order centres as traced by the order-fibre table (see 11.1.13).

Of course, for relative equalisation to be possible, each frame *must* have some overlap with the neighbouring ones in the set.

#### 11.1.15 Fiber frames preparation

The input odd/even fibre flat-field frames are normalized (each fibre is normalised to 1 integrated flux at each  $x$  position), the fibres falling within the region illuminated in the slit flat-field frames are selected, and a new set of odd/even fibre flat-field frames and other auxiliary frames (normalization, fibre boundaries) used by subsequent

<b>ESO</b>	<b>FLAMES-UVES Pipeline User Manual</b>	Doc:	VLT-MAN-ESO-19500-3016
		Issue:	Issue 5.0
		Date:	Date 2008-11-28
		Page:	87 of 110

data reduction commands are created. Any bad pixel in the input fibre flat-fields is filled with interpolated values, if possible. Each fibre is then correlated with a syntetic Gaussian profile which peaks at the position traced by the order-fibre table plus a variable offset.

The variable  $y$  offset is allowed to span the interval between **-MAXYSHIFT** and **+MAXYSHIFT**.

The half width at half maximum of the syntetic Gaussian pseudo-fibres is read from the **GAUSSFIBRESIGMA** keyword, and their total half width is read from the **GAUSSHALFWIDTH** keyword.

The offset corresponding to the maximum correlation is saved in FITS keywords, to be used as a “zero-point” correction for fibre asymmetry when measuring actual fibre offsets on other frames. For this correction to be effective, subsequent correlations are forced to use the same Gaussian parameters (i. e. **GAUSSFIBRESIGMA** and **GAUSSHALFWIDTH**) which were used here.

### 11.1.16 Fiber frames normalization

The all-fibre flat-field is (according to the value chosen for the parameter **method**) optimally or standard extracted.

**Background determination** First, if the **BKGFITINLINE** keyword is set to “yes” (as it is by default), a polynomial background is fit to the positions in the background table `bkg_prefix.fits` (prefix='l' or 'u'), and subtracted from the all-fibre flat frame (see ?? for details).

Setting **BKGFITINLINE** to “no” disables inline background fit and subtraction.

**Y shift determination** Then the y-axis shift of the orders/fibres in the all fibre flat-field frame is found by correlating the frame with a synthetic one composed of Gaussian-shaped pseudo-fibres centred on the positions traced by the order-fibre table, with a variable offset in the direction perpendicular to dispersion. The maximum of this correlation is found, and the corresponding offset is corrected for the “zero-point” shift computed by the fibre preparation (see 11.1.15) algorithm.

For robustness, this step perform a preliminar correlation function shape determination and a rough search for its maximum. This feature can be controlled by the keyword(s) **DRS\_COR\_MAX\_FND**, **DRS\_COR\_DEF\_RNG**, **DRS\_COR\_DEF\_PNT** that are set by the corresponding flames\_obs\_scired parameters `cor_max_fnd`, `cor_def_rng`, and `cor_def_pnt`.

The  $y$  offsets of each fibre are saved as a FITS keyword called **YSHIFT** in the science frame. The correlation step can be disabled by setting the **max\_shift** parameter to zero or, equivalently, setting to zero the **MAXYSHIFT** keyword and leaving the **max\_shift** parameter blank. In this case, the values already contained in the **YSHIFT** descriptor are used.

**Y shift correction** The normalised fibre flat-field frames are then shifted to bring them to coincidence with the above offset and are multiplied by the slit flat-field frames to approximate fibre flat-field frames exactly matching the all-fibre flat-field frame. These frames are used to perform an optimal or standard extraction on the science frame, including a deconvolution of adjacent fibres (the deconvolution coefficients are computed directly from the odd/even fibre flat-fields).

<b>ESO</b>	<b>FLAMES-UVES Pipeline User Manual</b>	Doc:	VLT-MAN-ESO-19500-3016
		Issue:	Issue 5.0
		Date:	Date 2008-11-28
		Page:	88 of 110

**Fibre spectra saving** The extracted “spectra” are saved to be later used as relative normalisation factors between fibres in the subsequent extraction of science frames.

#### 11.1.17 Fiber frames extraction

This data reduction step performs an optimal or a standard spectrum extraction. Three are the supported modes: “normal” (**extract**=opt/sta), the one suggested to the user; “fast” (**extract**=fop/fst), which uses input slit flat fields to clean the fibre PSFs but skips the steps of Y-shift determination and compensation (see below), a bit faster than “normal”; “quick” (**extract**=qop/qst) which does not uses the slit flat fields and does not correct for fibre shifts, the fast and less accurate mode.

Then are determined the background and the Y shift of each fibre of the science frame to the odd-even base, and each fibre position is corrected as described in the fibre frames normalization algorithm (see 11.1.16) for the all flat frame.

**Extraction.** The frames determined as explained above are then used to perform an optimal extraction on the science frame, including a deconvolution of adjacent fibres (the deconvolution coefficients are computed directly from the odd/even fibre flat-fields).

Notes:

- Should the correlation step fail to detect a reasonable  $y$  shift, or if speed is of utmost importance, it is advisable to use standard extraction, disabling the correlation and reducing the integration window to the core of the fibres only, rather than optimal extraction, to minimise fibre to fibre contamination, at the price of losing a considerable fraction of the signal.
- In cases of well-behaved set of frames, however, optimal extraction gives the best signal/noise ratio, especially for faint objects. The calculated spectra are finally written on the disk on one set of files for each illuminated fibre.
- The standard extraction includes the deconvolution of adjacent fibres, being them neighbouring fibres of the same order or the first and last fibres of adjacent orders. The deconvolution coefficients are computed directly from the shifted fibre flat-field frames.
- Since in the case of fast optimal extraction no attempt at all is made to compensate for any  $y$  offset of the science frame with respect to the normalised fibre flat-field frames, this method should be selected with some caution.
- In cases in which the  $y$  offset determination is a problem, the fast extraction should be used instead, limiting the integration window to the cores of the fibres to minimize ill effects at the cost of losing a considerable fraction of the signal.

#### 11.1.18 Fiber frames merging

First the extracted spectrum, its variance and its bad-pixel mask are rebinned.



<b>ESO</b>	<b>FLAMES-UVES Pipeline User Manual</b>	Doc:	VLT-MAN-ESO-19500-3016
		Issue:	Issue 5.0
		Date:	Date 2008-11-28
		Page:	89 of 110

Then the orders are optimally merged, computing every pixel in the merged spectrum as a weighted average of all good pixels available covering its wavelength bin, with weights equal to the inverse of the respective variances.

Important information is contained in the binary table extensions.

### 11.1.19 Interorder background subtraction

Two different methods are used to subtract the scattered light background:

- Spline method. The input image is sampled at half-integer order locations at **backsub.npoints** equally spaced sample points. According to the user defined value of the parameter **backsub.mmethod**, the median or the minimum values computed in a subwindow of height ( $2 * \text{backsub.radiusy} / \text{biny\_size} + 1$ ). The window width is given by the distance between the sample positions.

After the spline interpolation, the interorder background image is filtered using an average filter with radius (**backsub.smoothx**, **backsub.smoothy**).

The spline degree is set by using **backsub.sdegree**. Currently, only splines of degree 1 are supported (i.e. linear interpolation). If the **backsub.mmethod** parameter is set to 'no', no background subtraction is done.

- Polynomial method. The polynomial method is used in the order definition recipe because the order locations, required for the spline method, are not known at this initial stage.

A low degree 2d polynomial is fit to a subset of the image pixels and outlier points (such as the orders themselves) that have large positive residuals are continuously rejected (one-sided kappa-sigma clipping).

The input image is sampled on a regular grid with mask size (**image\_width / backsub.npoints** , **image\_height / backsub.npoints**).

After the initial order line detection, the interorder background is sampled (at locations separated by **image\_width / backsub.npoints**) between the order lines.

### 11.1.20 Hough transform

The Hough transform is computed by sampling the input image at every column separated by **samplewidth**.

Each echelle order maps to a peak in the Hough image. After detecting a peak, the peak itself and the area around it are deleted to avoid redetecting the same feature.

The accurate peak locations are calculated as the centroid of the area around the local maximum in the Hough space.

The number of order lines to detect is specified by the parameter **norders**. If this parameter is set to zero, the function will detect lines until the intensity of the next candidate drops to below a fraction **pthres** of the dimmest line.

An important parameter for the peak removal to work is the (approximate) interorder spacing. This parameter is estimated as the raw image height divided by the predicted number of orders (**norders**). In automatic mode it



<b>ESO</b>	<b>FLAMES-UVES Pipeline User Manual</b>	Doc:	VLT-MAN-ESO-19500-3016
		Issue:	Issue 5.0
		Date:	Date 2008-11-28
		Page:	90 of 110

is estimated as the first minimum of the auto-correlation function along the column in the Hough image which contains the global maximum. This fully automatic way of detecting the orders assumes only that the interorder spacing does not vary too much as function of order number.

Possible order line slopes range from **minslope** to **maxslope**. The resolution of the slope-axis in Hough space is defined by **sloperes**.

### 11.1.21 Order tracing

This algorithm, shared with the echelle-slit mode order tracing, is applied on fibre mode to trace the simultaneous calibration fibre order flat.

The parameter **use\_guess\_tab** value, defaulted to 1, allow the user to benefit of the information contained in the guess order table in a different way:

- 0: No usage.
- 1: Use the guess order table to set the lower/upper Y rows where order are searched.
- 2: The order table try to fully match the guess one.

The order tracing is performed as follows:

- First, all orders are traced in both directions starting from the center which is inferred from the solution of the Hough detection (if **use\_guess\_tab** is set to 0), or from the “guess” order table (if **use\_guess\_tab** is set to 1 or 2).

A Gaussian is fit to the order line at x-positions separated by the parameter **tracstep**. The trace stops if the intensity of the order line drops below the threshold defined by the **minthresh** value in an x-range determined by the parameter **maxgap**.

- Then each order is fit with a straight line, and the entire order is rejected if the RMS is large compared with the average RMS.
- A global polynomial of automatic degree is fit to all orders, and individual points are rejected using kappa-sigma clipping. Alternatively, the user can define polynomial degrees using the parameters (**defpol1**, **defpol2**).

### 11.1.22 Line Search

In fiber mode is used the same algorithm as in echelle-slit mode.

A number of emission lines defined by the range **minlines-maxlines** (both inclusive) is searched for in the extracted arc lamp spectrum. This is achieved by adjusting the detection threshold level until an appropriate number of lines is detected.

A 5-parameters Gaussian fit, including the continuum slope, is made if the line peak is above the threshold with respect to the local interorder background level, which is defined as the median of a window of width ( $2 * \text{range}$

<b>ESO</b>	<b>FLAMES-UVES Pipeline User Manual</b>	Doc:	VLT-MAN-ESO-19500-3016
		Issue:	Issue 5.0
		Date:	Date 2008-11-28
		Page:	91 of 110

+ 1) centered on the current position. Finally, doublets (defined as lines with positions within 2.0 pixels) are removed from the set of detected lines.

### 11.1.23 Wavelength calibration first solution determination

In fiber mode is used the same algorithm as in echelle-slit mode.

An initial dispersion relation is obtained by fitting the relation

$$\lambda * m = f(x, m),$$

to a guess line table containing associations from  $(x, m) = (\text{pixel}, \text{order})$  to wavelengths. Here  $f$  is a 2d polynomial.

A systematic x-shift up to **shiftmax** pixels is recovered by finding the maximum position of the cross-correlation function and applying this shift to the initial dispersion solution. The resolution of the cross-correlation function is defined by **shiftstep**. The parameter **shiftole** defines the tolerance in pixels for the line match. The default polynomial **degree** is 5.

### 11.1.24 Line identification

In fiber mode is used the same algorithm as in echelle-slit mode.

The wavelength calibration starts from a first guess dispersion solution.

Lines are iteratively identified and a dispersion solution is fit, until no new identifications can be made. After the first convergence, all identifications are cleared (to remove possible false identifications), and the loop repeats, but this time ignoring lines with residuals worse than **tolerance** and lines with residuals worse than **kappa** sigma. If set to a negative value, the polynomial degree is automatically adjusted based on the line residuals.

Identifications are made based on a match between the detected line predicted wavelength,  $\lambda_{com}$  and a catalogue wavelength,  $\lambda_{cat}$ . An identification is made if

- The nearest catalogue wavelength is within two linewidths of the predicted wavelength:  $|\lambda_{cat} - \lambda_{com}| < 2 * \sigma$ , where  $\sigma$  is the detected line width,
- The distance to the 2nd nearest neighbours (in the spectrum as well as in the catalogue) is much larger than the residual of the match  $|\lambda_{cat} - \lambda_{nn}| * ALPHA > |\lambda_{cat} - \lambda_{com}|$ , (where **ALPHA** is a "safety parameter" less than one.)
- The nearest neighbour (in the spectrum and in the catalogue) is farther away than the average tolerance distance, which measures the precision of the identifications: **tolerance** < **ALPHA** \*  $|\lambda_{cat} - \lambda_{nn}|$ . Refer to the source code for the exact definition of the **tolerance**.

The purpose of the first criterion is to make the correct identifications. The purpose of the latter two criterions is to avoid making incorrect identifications.

<b>ESO</b>	<b>FLAMES-UVES Pipeline User Manual</b>	Doc:	VLT-MAN-ESO-19500-3016
		Issue:	Issue 5.0
		Date:	Date 2008-11-28
		Page:	92 of 110

## 11.2 Recipes

In the following sections we are going to describe the recipes provided by the FLAMES-UVES pipeline. Recipe common parameters are: **debug**, to activate the debug mode which saves intermediate results to disk and **plotter** to activate the plotting facility (gnuplot). Additionally each recipe can be customized by modifying specific data reduction parameters which are described in the following section or in the referred algorithm description.

### 11.2.1 Master bias generation

We refer the reader to the UVES context documentation for the recipe `uves_cal_mbias`. In fact FLAMES data reduction uses that command (see [1]).

### 11.2.2 Master dark generation

We refer the reader to the UVES context documentation for the recipe `uves_cal_mdark`. In fact FLAMES data reduction uses that command (see [1]).

### 11.2.3 Spectral format guess

The recipe `flames_cal_predict` expects as input a fibre formatcheck frame which is obtained by illuminating the calibration fibre with a ThAr lamp, and a reference ThAr lines table. Optional input is a reference formatcheck frame.

This recipe runs the UVES physical model as for ECHELLE data. For details we remind the user to [1]. The only difference is that in case of FIBER mode data, to take into account of the offset between the calibration fiber and the central position of the order, X and Y offsets are set depending on the given instrument setting (wavelength, plate number).

### 11.2.4 Single fibre order table generation

To get an order table from a simultaneous calibration (single) fibre flat field one can use the recipe `flames_cal_orderpos` providing in input a guess order table the one generated by `flames_cal_predict`.

### 11.2.5 Master slit flat field generation

The recipe `flames_cal_mkmaster` generates a master flat frame from a set of slit flats.

To have full coverage of all 9 fibres and wide enough windows to determine the inter order background are taken three sets of three slit flats each, each set taken at different Y positions. From each set it is computed a master flat field as for echelle mode data. The recipe divides the raw slit flat fields in groups each having the same Y position (having values of the FITS keyword ESO.INS.SLIT3.X1ENC and ESO.INS.SLIT3.X2ENC

<b>ESO</b>	<b>FLAMES-UVES Pipeline User Manual</b>	Doc:	VLT-MAN-ESO-19500-3016
		Issue:	Issue 5.0
		Date:	Date 2008-11-28
		Page:	93 of 110

within 5 encoder steps) and next on each group applies a standard master flat field generation (without inter-order background subtraction as this operation is performed later in the data reduction). The output are a set of (usually three) master slit flat fields each corresponding to a given Y position.

### 11.2.6 Slit flats preparation, odd-even-all fibre flats preparation, normalization and fibre-order table determination

The recipe `flames_cal_prep_sff_ofpos` prepares the input slit flat-field and fibre flat-field frames for inclusion in the calibration database. First are processed the lower-chip data and next the upper-chip data. The slit FF frames are normalised to be equal where they overlap, then a minimal set of frames covering the largest  $y$  interval is produced.

This recipe executes the following algorithm in sequence:

- A bad pixel map is generated (see [11.1.10](#)),
- A variance and a bad pixel map are associated to the input frames (see [11.1.9](#)),
- All fibres are located and an order position table is determined (see [11.2.6](#)),
- Suitable positions for background-light extraction are determined (see [11.1.12](#)),
- Input slit flats are prepared (see [11.1.4](#)),
- Input fibre flats are prepared (see [11.1.15](#)),
- Fibre to fibre normalizations are computed (see [11.1.16](#)).
- All output products are classified.

To improve robustness (of this command and `flames_cal_wavec`, `flames_obs_scired` commands) this step checks that for each detected order all the fibres are traced and if not a fibre tracing is forced. Moreover the recipe `flames_cal_prep_sff_ofpos`, to prevent problems in the wavelength calibrations, excludes from the order tracing the bottom right and upper left hand side corners. This is done setting appropriate values of keywords **DRS\_SCAN\_MIN** and **DRS\_SCAN\_MAX**.

### 11.2.7 Wavelength Calibration

The recipe `flames_cal_wavec` performs the wavelength calibration using previously determined solutions for the first guess line table and the order table. This recipe is a modification of the **uves\_cal\_wavec** recipe to reduce UVES echelle data, in the sense that instead of determining the solution corresponding to tree order traces, one at the center and the other two at 15 pixels up and down, now are determined the solutions corresponding to each of the fiber traces as indicated by the **FIBREMASK** FITS keyword. For details we remind the user to [1].

ESO	FLAMES-UVES Pipeline User Manual	Doc:	VLT-MAN-ESO-19500-3016
		Issue:	Issue 5.0
		Date:	Date 2008-11-28
		Page:	94 of 110

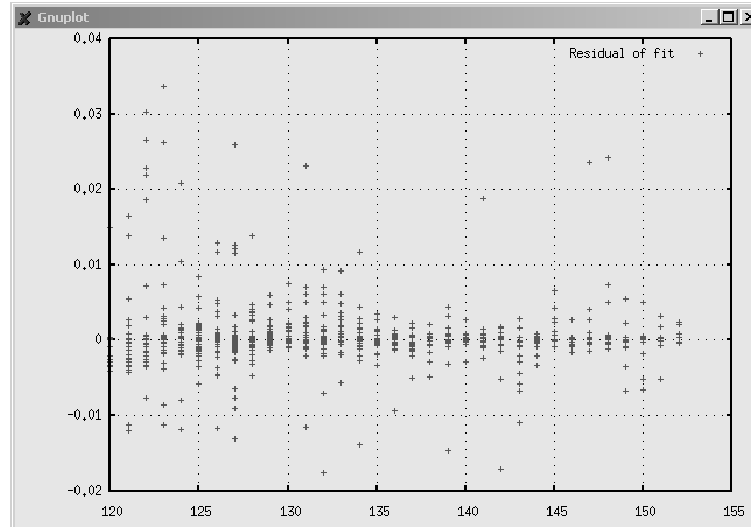


Figure 11.2.1: Wavelength calibration residuals (A) as function of order number.

### 11.2.8 Science data reduction

This recipe performs the actual spectrum extraction, either optimal or standard. The normalized odd/even fibre flat-field frames output of previous commands are compared to the science frame, any y-axis shift between them is computed and compensated by shifting the odd/even fibre flat-fields; an optimal/standard extraction is performed on the science frame, and the resulting spectra are created, one set of files for each illuminated fibre. Next, using the line table created by **flames\_cal\_wavecalf**, the extracted spectra are rebinned and merged. Finally, the products are classified. This commands supports different extraction modes. The normal extraction mode (the one just described, using optimal or standard extraction, setting the parameter **extract** to “opt” or “sta”) is the one we suggest to the user. Fast extraction mode (**extract** set to “fop” or “fst”) is faster as skip the correlation step (but as such it may have also lower quality). Quick extraction mode, which does not use at all the slit flat fields in the data reduction chain, and does not do any correlation step, was implemented in answer to a Paranal Operation request, and if used, requires that the complete data reduction (recipes **flames\_cal\_prep\_sff\_ofpos** is performed in a coherent way (using quick data reduction mode).

In case the correlation step fails (usually due to a large Y shift between the observation and the calibrations) we suggest the user to do the following. The shape of the correlation function is contained in tables **cor\_shape\_x.fits** (x=l or x=u respectively for EEV or MIT chips). From this table one can get the Y offset at which the correlation function has a maximum.

Next repeat the extraction appropriately setting the correlation function parameters (through parameter **extract**). In this case usually **drs\_cor\_def\_rng** can be decreased to 2 or even 1.

For example to do optimal extraction with correlation (OPT) function search (Y) in the range [-2,2], using 2\*25+1 pixels and having measured a -4.5 pixel shift, one may give the command:

```
esorex flames_obs_reduce --cor_pnt=25 --cor_rng=2 --cor_off=4.5 flames_obs_reduce.sof
```

ESO	FLAMES-UVES Pipeline User Manual	Doc:	VLT-MAN-ESO-19500-3016
		Issue:	Issue 5.0
		Date:	Date 2008-11-28
		Page:	95 of 110

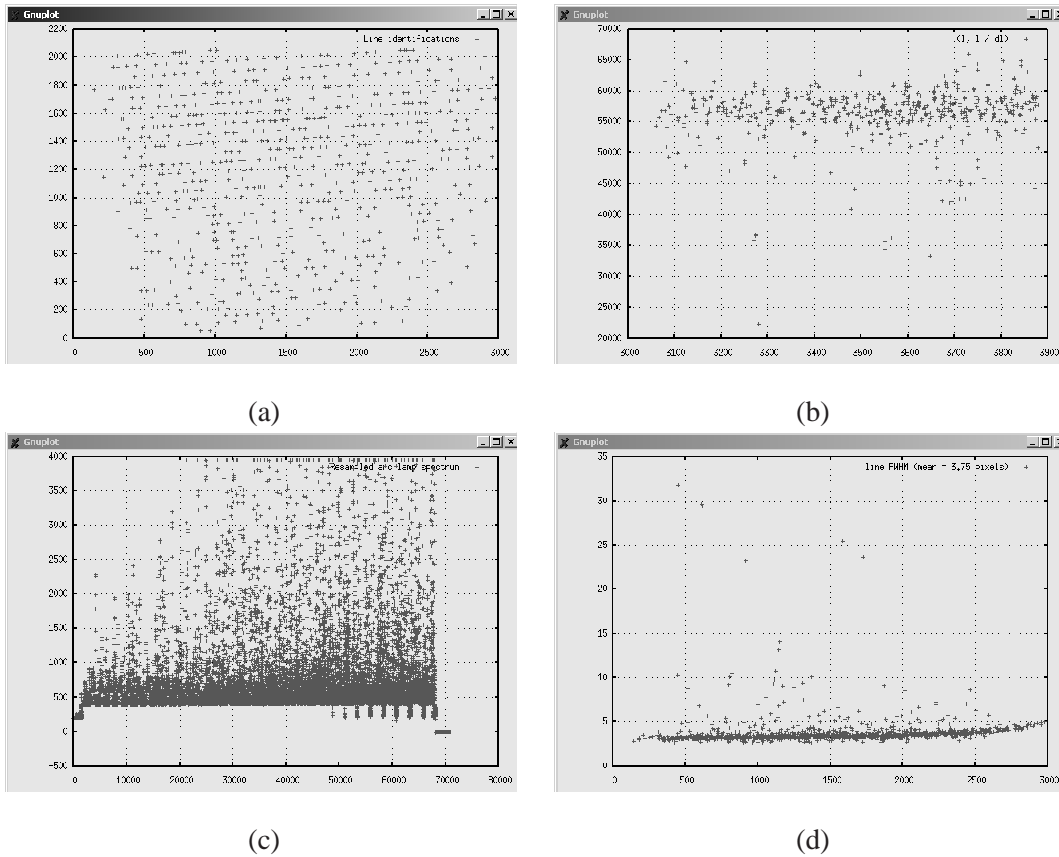


Figure 11.2.2: These plots display: (a) the distribution of identified lines across the detector, (b) the resolving power as a function of wavelength, (c) the resampled ThAr spectrum, (d) the line FWHM as a function of wavelength.

### 11.2.9 Full UVES-FIBRE data reduction

The recipe `flames_obs_redchain` runs the full FLAMES-UVES reduction chain. It runs in a cascade the following recipes:

- `uves_cal_mbias` (if no master bias is provided)
- `uves_cal_mdark` (if darks are provided)
- `flames_cal_predict` (if no guess order and line tables are provided)
- `flames_cal_orderpos` (if no order table is provided)
- `flames_cal_mflat` (if no master slit flat is provided)
- `flames_cal_prep_sff_ofpos` (if no `slitff*` and `fibreff*` are provided)
- `flames_cal_wavecal` (if no dispersion solution line table is provided)
- `flames_obs_scired` (unless the option **scired** is set to FALSE)

<b>ESO</b>	<b>FLAMES-UVES Pipeline User Manual</b>	Doc:	VLT-MAN-ESO-19500-3016
		Issue:	Issue 5.0
		Date:	Date 2008-11-28
		Page:	96 of 110

## A Installation

This chapter gives instructions on how to obtain, build and install the UVES pipeline (ECHELLE and FIBRE modes). Even if this chapter is kept as up-to-date as possible, it may not be fully applicable to a particular release. This might especially happen for patch releases. One is therefore advised to read the installation instructions delivered with the UVES pipeline distribution kit. These release-specific instructions can be found in the file `README` located in the top-level directory of the unpacked UVES pipeline source tree. The supported platforms are listed in Section A.1. It is recommended reading through Section A.3 before starting the installation.

A bundled version of the UVES pipeline with all the required tools and an installer script is available from [www.eso.org/pipelines](http://www.eso.org/pipelines).

### A.1 Supported platforms

The UVES pipeline has been verified to install and execute correctly with EsoRex on the VLT target platforms:

- Linux (glibc 2.1) (tested Fedora core 6 and Scientific Linux 4.0)

and on

- Mac Darwin 8.9.0,

using the GNU C compiler (version 3.2 or newer). Correct execution using Gasgano has been verified on

- Scientific Linux 4.0,
- Linux Fedora core 6,
- Mac Darwin 8.9.0.

### A.2 Requirements

To compile and install the UVES pipeline one needs:

- the GNU C compiler (version 3.2 or later),
- the GNU `gzip` data compression program,
- a version of the `tar` file-archiving program and
- the GNU `make` utility.

For Gasgano support one needs in addition

- the Java Development Kit (version 1.5 or newer)



<b>ESO</b>	<b>FLAMES-UVES Pipeline User Manual</b>	Doc:	VLT-MAN-ESO-19500-3016
		Issue:	Issue 5.0
		Date:	Date 2008-11-28
		Page:	97 of 110

### A.3 Building the UVES pipeline

The UVES pipeline distribution kit contains:

uves-manual-5.0.pdf	The UVES pipeline manual
install_pipeline	Install script
qfits-6.2.0.tar.gz	QFITS 6.2.0
cfitsio-2510.tar.gz	CFITSIO 2510
cpl-4.1.0.tar.gz	CPL 4.1.0
esorex-3.6.8.tar.gz	esorex 3.6.8
gasgano-2.2.7.tar.gz	GASGANO 2.2.7
uves-4.3.0.tar.gz	UVES 4.3.0
uves-calib-4.3.0.tar.gz	UVES static calibration files 4.3.0

Here is a description of the installation procedure:

1. Change directory to where you want to retrieve the UVES pipeline 4.3.0 package. It can be any directory of your choice but not:

```
$HOME/gasgano
$HOME/.esorex
```

2. Download from the ESO ftp server, [www.eso.org/pipelines](http://www.eso.org/pipelines), the latest release of the UVES pipeline distribution.
3. Verify the checksum value of the tar file with the cksum command. **cksum uves-kit-4.3.0.tar.gz**
4. Unpack using the following commands:  
**gunzip uves-kit-4.3.0.tar.gz tar -xvf uves-kit-4.3.0.tar**

Note that the size of the installed software (including *Gasgano*) together with the static calibration data is about 27Mb.

5. Set the environment variable JAVA\_HOME to the directory where you have the JDK 1.5 installed. If this value is not set, the installation script will try to guess it, but if no JDK is found, the gasgano distribution will not be installed; QFits, CPL, EsoRex and the pipeline will be installed anyway.
6. Install: after moving to the top installation directory,  
**cd uves-kit-4.3.0**  
it is possible to perform a simple installation using the available installer script (*recommended*):  
**./install\_pipeline**
7. Check the installation log: probably this will suggest you to set the environment variable CPLDIR and to extend your PATH.

<b>ESO</b>	<b>FLAMES-UVES Pipeline User Manual</b>	Doc:	VLT-MAN-ESO-19500-3016
		Issue:	Issue 5.0
		Date:	Date 2008-11-28
		Page:	98 of <a href="#">110</a>

By default the script will install the UVES recipes, *Gasgano*, *EsoRex*, all the necessary libraries, and the static calibration tables, into a directory tree rooted at `$HOME`. A different path may be specified as soon as the script is run.

The only exception to all this is the *Gasgano* tool, that, if you have the required proper installation of the JDK (version 1.5), will always be installed under the directory `$HOME/gasgano`. Note that the installer will move an existing `$HOME/gasgano` directory to `$HOME/gasgano.old` before the new *Gasgano* version is installed.

Important: the installation script would ensure that any existing *Gasgano* and *EsoRex* setup would be inherited into the newly installed configuration files (avoiding in this way any conflict with other installed instrument pipelines).

<b>ESO</b>	<b>FLAMES-UVES Pipeline User Manual</b>	Doc:	VLT-MAN-ESO-19500-3016
		Issue:	Issue 5.0
		Date:	Date 2008-11-28
		Page:	99 of 110

## B FLAMES-UVES pipeline data reduction parameters

The FLAMES-UVES pipeline gives the user the possibility to customize the data reduction. This can be done modifying the settings in the file

`${PIPE_HOME}/flames/pipe/proc/flames_def_drs_par.prg`

We suggest to do it with care, being sure of what one does and at his own risk. Here we describe the tuning parameters, and suggest also whether each MIDAS keyword (data reduction parameter) should or should not be changed.

- **NICE\_CREA**: this keyword affect the display (Y) or not (N) of some frames. Default is N. Can be modified. If set to Y the data reduction is slightly slowed.
- **DECENTSNR**: How large should the SNR on a fibre be in a calibration frame, at a given order and x, for that slice to be considered “good”? Default: 10. Suggested to not change it. Since it is used while preprocessing the flat field calibration data, for changes to take effect one ought to rerun some or all of the relevant steps (e. g. PREPSLIT/FLAMES, PREPFIBRE/FLAMES, PRENOR/FLAMES. . . ).
- **MATCHTHRESH**: How close should a detected fibre be to the position predicted by the zero order approximation, so that it is recognized and labeled by the matchorder step. Default 4. Suggested to not change it. Setting it too low will cause orders/fibres to be skipped, setting it too large will cause traces to be assigned to the wrong order/fibre. Only used in OFPOS/FLAMES.
- **HALFIBREWIDTH**: Half size of each fibre on the detector measured so that the entire size contains 100% of the flux. Default: 7.5. This parameter is only read in OFPOS/FLAMES, and its value hence stored in a descriptor and used throughout data reduction, for consistency. If one wants to change it, he should probably repeat the whole data reduction procedure for the new value to take effect (don’t do it).
- **MAXFIBRES**: Maximum number of fibres. Default: 9. Keep it. If decreased the DRS will probably crash.
- **SAV\_BORD\_SZ**: Amount of pixels from the 50% intensity threshold of the SlitFF size which is discarded on both sides of the slitFF to be sure to be in its flat part. A large value ( $\sim 5$ ) makes you sure that only the flat part of the slitFF is used, and probably the final quality will be better; on the other hand, this will shrink the effective usable size of the used slitFF so that you need to have enough overlapping between them and enough slitFF to cover all the fibres you like to extract. A small value (2) may introduce some quality defect if the flat part of the slitFFs is not really flat. Default: 3. It may be changed if really necessary with the given advices. Only used in MSFFSZ/FLAMES. If you *really* want to tune the YSHIFT and HALFWIDTH descriptors of the slitFFs, measure them with GET/GCURLSOR on a cross-dispersion cut of the frames and set the desired values manually.
- **X\_WIND\_SIZE**. Half X window size for median filter used in determining the half size of each slitFF and its offset. Default: 3. Leave it as it is. Only used in MSFFSZ/FLAMES. See the comments for SAV\_BORD\_SZ.
- **Y\_WIND\_SIZE**. Half Y window size for median filter used in determining the half size of each slitFF and its offset. Default: 5. Leave it as it is. Only used in MSFFSZ/FLAMES. See the comments for SAV\_BORD\_SZ.

<b>ESO</b>	<b>FLAMES-UVES Pipeline User Manual</b>	Doc:	VLT-MAN-ESO-19500-3016
		Issue:	Issue 5.0
		Date:	Date 2008-11-28
		Page:	100 of 110

- **Y\_SEARCH\_WIND.** Search window size to filter image used in determining the half size of each slitFF and its offset. Default: 100. Leave it as it is. Only used in MSFFSZ/FLAMES. See the comments for SAV\_BORD\_SZ.
- **ORD\_TRESH.** Number of orders cut either on max and min detected orders. Used in determining the half size of each slitFF and its offset. We cut some orders to be sure to get Y intercepts on the cross-order cut at the detector center within the detector. Default: 2. It may be slightly increased, but it is preferable to leave it as it is. Only used in MSFFSZ/FLAMES. See the comments for SAV\_BORD\_SZ.
- **N\_CLIP\_MED.** K-S clipping iterations over median. Used in determining the half size of each slitFF and its offset. Default: 4. Leave it as it is. Only used in MSFFSZ/FLAMES. See the comments for SAV\_BORD\_SZ.
- **N\_CLIP\_AVG.** K-S clipping iterations over average. Used in determining the half size of each slitFF and its offset. Default: 2. Leave it as it is. Only used in MSFFSZ/FLAMES. See the comments for SAV\_BORD\_SZ.
- **INT\_TRESH.** Signal fraction of the slitFF top value used to define the borders of the flat part of the slitFF. A value as small as 0.5 is good to be sure to detect slitFF borders. Default: 0.5. Leave it as it is. Only used in MSFFSZ/FLAMES. See the comments for SAV\_BORD\_SZ.
- **MAXYSHIFT.** This is the halfwidth of the interval of y shifts which are sampled when searching the local maximum of the correlation. Default: 3. Making it smaller may make the search for the maximum slightly faster, but it *will* cause it to fail if the actual y shift is outside this interval. Making it larger may bracket more than one maximum and cause a spurious y shift determination. Don't touch this if you don't need to.
- **NBTRACES.** 0. Number of traces (order times fibres). Default 0.
- **MINFIBREFRAC.** Minimum fibre fraction coverage for extraction. Default: 0.3. Leave it as it is.
- **BKGFITNLINE.** Inline background fitting. Default yes. Keep it.
- **BKGFITMETHOD.** Background fitting method. Possible values are ALL, MEDIAN, MINIMUM, AVERAGE. Default: AVERAGE.
- **BKGBADSCAN.** Background table bad pixel frame scanning switch. Possible values are FRACTION, ABSOLUTE and NONE. Default: NONE.
- **BKGBADWIN.** Background table bad pixel frame scanning window size. Default: 50,50
- **BKGBADMAXFRAC.** Background table bad pixel frame scanning threshold fraction. Default: 0.02
- **BKGBADMAXTOT.** Background table bad pixel scanning threshold number. Default: 200.
- **BKG\_MAX\_IO\_WIN.** Number of background windows in each full inter order. Default: 500.
- **DRS\_BKG\_FIT\_POL.** Polynomial degree used for BKG fit order. Default 2,2
- **BKG\_XY\_WIN\_SZ.** x,y maximum size of each background window. Default 6,2.
- **DRS\_FILT\_HW\_XY.** Size of X half window and Y half window. Default 2,1.

<b>ESO</b>	<b>FLAMES-UVES Pipeline User Manual</b>	Doc:	VLT-MAN-ESO-19500-3016
		Issue:	Issue 5.0
		Date:	Date 2008-11-28
		Page:	101 of 110

- **DRS\_FILT\_IMAX.** Maximum filtering iterations in frame preparation. Default 300.
- **DRS\_FILT\_KS.** Threshold above which pixel are flagged as bad in filtering in frame preparation. Default 10.
- **DRS\_FILT\_SAT\_SW.** Do you want mask saturated pixels in frame preparation ([yes]/No)?. Default is YES.
- **DRS\_FILT\_MASK.** Do you want a filter/generated bad pixel mask: ([NONE]/MEDIAN)?. Default is NONE.
- **MAXYSHIFT.** Half width of the interval to scan for correlation, when determining y shift. Default is 3.
- **GAUSSFIBRESIGMA.** Gaussian pseudofibre total halfwidth for correlation. Default: 1.5. Only read in PREPFIBRE/FLAMES, then the value is stored in a table descriptor and used throughout, for consistency. Don't make it much smaller or it will make the correlation numerically fragile due to increased pixelisation effects. Don't make it too large or signals outside the fibres will spuriously affect the correlation. As a rule of the thumb,  $3 \cdot \text{GAUSSFIBRESIGMA} \lesssim \text{GAUSSHALFWIDTH} \lesssim \text{HALFIBREWIDTH}$ .
- **GAUSSHALFWIDTH** Gaussian pseudofibre total halfwidth for correlation. Default: 6. Only read in PREPFIBRE/FLAMES, then the value is stored in a table descriptor and used throughout, for consistency. See comments above for GAUSSFIBRESIGMA.
- **MAXBACKITERS** This is the maximum number of kappa-sigma clipping iterations which we are willing to perform in background fitting. Default 20.
- **MAXDISCARDFRACT** This is the maximum fraction of windows/pixels which we are willing to discard by kappa-sigma clipping in each iteration of the background fitting loop. Default: 0.1. Making it larger will cause the background fitting to converge in less iterations, but can make it more fragile. Making it smaller will render the kappa-sigma clipping more cautious, hence it will require more iterations for all bad pixels to be discarded and convergence reached. The default chosen works well in most situations, leave it as it is.
- **MAXCORRITERS.** This is the maximum number of iterations which we are willing to perform in correlation. Default: 30.
- **CORRELTOL.** This is the absolute accuracy with which the maximum in the correlation function must be determined. Default 0.005. Leave it as it is.
- **CORRELXSTEP.** This is the X step to use while computing the correlation: it must be a positive integer, 1 means "use all pixels" 2 means "use every other pixel", 3 means "use one every three" etc.. Default: 1. A larger value can make the correlation much faster, but it will also make it more fragile in case of spectra dominated by strong and narrow emission lines.
- **GAUSSCORRELSCL.** Obsolete (to be removed in future releases).
- **GAUSSCORRELWND.** Obsolete (to be removed in future releases).
- **MAXCLEANITERS** This is the maximum number of cleaning iterations to be tried on a given slice in flames\_prep\_fibre, before giving up on its normalisability and falling back to the second cleaning strategy. Default 10. Only used in PREPFIBRE/FLAMES.

<b>ESO</b>	<b>FLAMES-UVES Pipeline User Manual</b>	Doc:	VLT-MAN-ESO-19500-3016
		Issue:	Issue 5.0
		Date:	Date 2008-11-28
		Page:	102 of 110

- **MAXSINGLEPXFRC** This is the maximum acceptable fraction of the flux in a given slice in one single pixel; a higher fraction than this means that there was some numerical instability and/or unmasked bad pixel, and that it must be discarded. Default: .3. Only used in PREPFIBRE/FLAMES. Leave it as it is.
- **MAXOPTITERSINT** This is the maximum number of iterations which we are willing to perform in optimal extraction. Default: 25.
- **MINOPTITERSINT** This is the minimum number of iterations we are willing to perform in optimal extraction. Default: 2.
- **XKILLSIZE**
- **YKILLSIZE** When performing sigma-clipping in the optimal extraction, how many other adjacent pixels in the x and/or y direction(s) should be discarded along with the one exceeding threshold? A cosmic or cosmetic problem is likely to affect a spot larger than 1 pixel. Default: XKILLSIZE=1, YKILLSIZE=1. Larger values will cause a slightly faster convergence of optimal extractions, as clusters of bad pixels are removed in one go. However, some good pixels may get thrown away along with the bad ones. The default is the safe (and slow) choice.
- **MAXORDER.** In flames\_tracing, a few arrays are statically defined, and need to be large enough to contain all the orders which could be possibly found. We defined it to a hefty one thousand. Default 1000. You will hardly ever want to touch this.
- **DYRANGE** Width of the y half window to be used when performing order/fibre labeling. Default: 300. It works well, leave it as it is.
- **DYSTEP** Step of the scan in y while performing order/fibre labelling. Default 0.1. It works well, leave it as it is.
- **MASKTHRES** Defines the minimum value a rebinned mask must have in order to consider the corresponding pixel of the rebinned spectrum “good”. Since a pixel in the rebinned frame may be computed using also bad pixels, we want to throw out pixels which contain even a very small fraction of “badness”. Default: 0.99. It must be a number close enough to 1 to ensure that only good pixels were used to derive a rebinned value, but not too close, to avoid throwing away good pixels just due to numerical rounding errors. The default works well, leave it as it is or, at most, tune it slightly. Used only in MERGE/FLAMES
- **FRACSLICESTHRES** Defines the minimum fraction of slices that must be good, for a fibre to be considered covered enough at the end of fillholes, in flames\_prep\_fibreff, to avoid numerical instabilities in gaussselfcorrel if a fibre is only very barely covered by the slitFF frames. Default 0.3. Leave it as it is.
- **MAXROWS** The maximum permitted sizes for the order tables created by the DRS: Default 300000. You will hardly ever want to touch this.
- **MAXCOLS** The maximum permitted sizes for the order tables created by the DRS: Default 10. You will hardly ever want to touch this.
- **DRS\_PTHRE\_MAX** Maximum threshold for pixel saturation detection: Default  $6 \cdot 10^4$ . Good  $5 \cdot 6 \cdot 10^4$
- **DRS\_PTHRE\_MIN** Minimum threshold for pixel saturation detection: Default -20. Good -20

<b>ESO</b>	<b>FLAMES-UVES Pipeline User Manual</b>	Doc:	VLT-MAN-ESO-19500-3016
		Issue:	Issue 5.0
		Date:	Date 2008-11-28
		Page:	103 of 110

- DRS\_EXT\_MTD Extraction method. Default opt. Possible values opt: optimal; sta: standard; fop: fast-optimal; fst: fast-standard; qop: quick-optimal; qst: quick-standard.
- DRS\_K\_S\_THRE Kappa-sigma Threshold used in extraction. Default 10.
- DRS\_KSIGMA\_THRE Kappa sigma Threshold. Default 10.
- DRS\_EXT\_W\_SIZ Integration window size. Default 10. 10 is a good value if fibre deconvolution works fine.
- DRS\_MIDAS2FITS MIDAS to FITS conversion switch. At the end of a procedure it allow MIDAS to FITS format conversion. Default N. Possible value Y/N.
- DRS\_BIAS\_MTD Bias subtraction Method: M: Master Bias subtraction; <num> constant value subtraction; N: no bias subtraction.
- DRS\_FILT\_SW Filter switch. Default none. Possible values: none/average
- DRS\_SFF\_FIBFF Slitff\*-Fibreff frames preparation. Y: yes, slower, for DFO N: no, faster, for PSO once CDB is ready. Default Y
- DRS\_CREA\_BP\_TAB Switch to create (Y) or not (N) a bad pixels table. Default N
- DRS\_PHYSMOD\_REC Physical model auto recover switch: Y: yes the physical model autorecovers N: No the physical model does not autorecover. Default Y.
- DRS\_USE\_ORDEF Use (Y) or not (N) of the ORDERDEF as part of the odd/even fibre FF Keep it at Y to support SimCal mode. Default is N.
- DRS\_MER\_MTD What merging method are we using: FLAMES or ECHELLE. If ECHELLE metod REDUCE/FLAMES uses REBIN/ECHELLE and MERGE/ECHELLE comands. If FLAMES method REDUCE/FLAMES uses MERGE/FLAMES. ECHELLE method is more robust.  
Default is ECHELLE.
- DRS\_SCI\_RAW Produce or not raw science data. Default is Y.
- DRS\_VERBOSITY DRS MIDAS C Library verbosity level. Default LOW.
- DRS\_MES\_LEV MIDAS procedure main messaging level. Default is 4
- DRS\_CUBIFY Switch to activate creation and usage of slitff\* and fibreff\* cubes. Default is N.
- DRS\_COR\_MAX\_FND DRS strengthened correlation shape definition and pre search of maximum switch: Y do pre search of correlation's maximum. N don't do it (do only search of max with modified Brent method starting from points -3,0,+3). Default is Y.
- DRS\_COR\_DEF\_RNG DRS Correlation function's range: [-DRS\_COR\_DEF\_RNG,+DRS\_COR\_DEF\_RNG] pix 6.
- DRS\_COR\_DEF\_PNT DRS Correlation function's number of definition points. Effective number of points is 2\*DRS\_COR\_DEF\_PNT+1. Default is 25.



<b>ESO</b>	<b>FLAMES-UVES Pipeline User Manual</b>	Doc:	VLT-MAN-ESO-19500-3016
		Issue:	Issue 5.0
		Date:	Date 2008-11-28
		Page:	104 of 110

- DRS\_COR\_DEF\_OFF DRS Correlation function's center of definition. Default is 0.
- DRS\_WCAL\_MODE. Wavecal mode. Default is Auto.
- DRS\_WCAL\_RPLT. Default for wavecal resolution plots creations. Y.
- DRS\_WCAL\_FITS. Defaults for FITS table production in wavecal. Default N.
- DRS\_BASE\_NAME. Default base name for extracted science file prefix. Default is science.
- DRS\_SCI\_CUBE. For packing science products on planes and cubes. Default is N.
- DRS\_WCAL\_DC. Default for wavecal solution polynomial degree. Originally was 5, more precise, next to have better robustness we switched to 4.
- DRS\_WCAL\_TOL. Default for wavelength calibration polynomial solution parameter TOL. This sets a threshold. Points which discard by the found solution more than TOL pix are outliers. Current defaults is 0.6. A user willing to achieve a very accurate wavelength calibration solution, using the improved ThAr table (thargood\_3.tfits), may decide (at the risk of a decreased wavecal recipe robustness), to decrease this value down to 0.07.
- DRS\_WCAL\_PAR. DRS\_WCAL\_PAR=DRS\_WCAL\_MODE, DRS\_WCAL\_RPLT, DRS\_WCAL\_FITS, DRS\_WCAL\_DC, DRS\_WCAL\_TOL
- DRS\_SCI\_PAR1=DRS\_EXT\_MTD, DRS\_COR\_MAX\_FND, DRS\_COR\_DEF\_RNG, DRS\_COR\_DEF\_PNT, DRS\_COR\_DEF\_OFF
- DRS\_SCI\_PAR2=DRS\_BIAS\_MTD, DRS\_FILT\_SW, DRS\_PTHRE\_MAX
- DRS\_P8\_OFPOS Parameter P8 of hough/echelle. Default “.10,-.1,.1,.005,1.”. Reasonable values are:
  - 520 “.05,-.1,.1,.005,1.”
  - 580 “.10,-.1,.1,.005,1.”
  - 860 “.10,-.1,.1,.005,1.”
  - DRS\_P8\_OFPOS\_S1 First component of parameter P8 used in hough/echelle. Default 0.05,.10,.10.
  - DRS\_SCAN\_MIN
  - DRS\_SCAN\_MAX This MIDAS keyword set the parameter SCAN(1) and SCAN(2) respectively used by the command hough/echelle. The chosen values, allow to skip a CCD portion approximately equal to the Y span covered by eight fibres, in the preliminary determination of fibre traces done by the command OFPOS/FLAMES. This to prevent detection of very small portion of orders in the lower right hand side and upper left hand side part of a splitted frame. This help to make wavelength calibration and science extraction more robust for 860 setting. For other settings possible values are also 1,2048 which allows to cover the entire detector. We suggest the user to keep the chosen defaults: DRS\_SCAN\_MIN= 55,73,73 and DRS\_SCAN\_MAX= 1993,1975,1975 respectively for settings 520, 580, 860.
  - DRS\_MER\_DELTA Cutting values using in merging frames. Default 5,5.
  - DRS\_DEL\_SW Choice switch to set DRS\_MER\_DELTA values. Used in UVES echelle extraction. Possible values are D (Default: no setting), A (Automatic setting), and U (User defined setting). Default A.

<b>ESO</b>	<b>FLAMES-UVES Pipeline User Manual</b>	Doc:	VLT-MAN-ESO-19500-3016
		Issue:	Issue 5.0
		Date:	Date 2008-11-28
		Page:	105 of <a href="#">110</a>

- DRS\_SCI\_SIM Switch to activate use of SimCal setting.
- DRS\_SFF\_HW\_MIN Threshold to minimum value accepted in Slit Flat Width Half Width. Default is 10.

<b>ESO</b>	<b>FLAMES-UVES Pipeline User Manual</b>	Doc:	VLT-MAN-ESO-19500-3016
		Issue:	Issue 5.0
		Date:	Date 2008-11-28
		Page:	106 of <a href="#">110</a>

## C Full list of changes

List of differences between UVES-CPL pipeline releases 4.3.0 and 3.9.0

Changes in 4.3.0

Fixed compiler warnings

Changes in 4.2.9

Fixed bug in setting delt1/2 parameters. Updated their man page.

Set a corresponding proper range.

Changes in 4.2.8

Removes xxx\_ functions (originally used in place of cpl\_table)

fixed DFS05803 (seg fault due to missing check on array boundary)

Changes in 4.2.7

Fixed some problems from NRI and recipe unit tests

Changes in 4.2.6

Fixed mem leaks

Changes in 4.2.5

Add back uves\_cal\_lingain, uves\_cal\_ronbias

Changed uves\_cal\_ronbias to reflect detmon function API change

Added recipes/tests

Changes in 4.2.4

Fixed robustness problem corvel computation. (5110)

Fixed robustness problem (DFS5458)

Changes in 4.2.3

Added merge\_delt1, merge\_delt2 parameters to allow to chop off some of the

Fixed problem misalignment order-line guess tables on single fiber frames

Now user can set trans\_x/y also in fiber pmode to correct physical model

Changes in 4.2.2

DFS05346 Problem wavecal/merge: absorption lines may appear as doublets in merged spectra. Solved by decreasing wavecal poly degree default to 4.

Changes in 4.2.1

Added uves\_cal\_ronbias

DFS05302 Added PRO.REC1 keys in uves\_cal\_prep\_sff\_ofpos products

DFS05106 Now header of slit\_ff\_products is the one of the 1st frame of the

DFS05102 Implemented flames\_utl\_decubify

Changes in 4.2.0

DFS05140: changed LF.DO.CATG in DO.CATG for CD\_ALIGN\_RED to uniform to CD

DFS05139: Updated gasgano rule files to support classification of input raw data for uves\_cal\_lingain recipe

Added products of uves\_cal\_lingain to oca rules.

Changes in 4.1.0

Removed uves\_cal\_ccdtest from oca file

Added flames\_cal\_mkmaster in UVES.prefs

Disables uves\_cal\_ronbias

<b>ESO</b>	<b>FLAMES-UVES Pipeline User Manual</b>	Doc:	VLT-MAN-ESO-19500-3016
		Issue:	Issue 5.0
		Date:	Date 2008-11-28
		Page:	107 of <a href="#">110</a>

Updated cpl.version and esorex.version dependency IDs  
Changes in 4.0.0  
As 3.5.7 but updated release ID for PSO release (as includes  
fiber mode)

<b>ESO</b>	<b>FLAMES-UVES Pipeline User Manual</b>	Doc:	VLT-MAN-ESO-19500-3016
		Issue:	Issue 5.0
		Date:	Date 2008-11-28
		Page:	108 of <a href="#">110</a>

## D Abbreviations and acronyms

ANSI	American National Standards Institute
ASCII	American Standard Code for Information Interchange
CalibDB	Calibration Database
CPL	Common Pipeline Library
DFO	Data Flow Operations department
DFS	Data Flow System department
DMD	Data Management and Operations Division
DRS	Data Reduction System
ESO	European Southern Observatory
EsoRex	ESO-Recipe Execution tool
FITS	Flexible Image Transport System
FOV	Field Of View
FPN	Fixed Patter Noise
GUI	Graphical User Interface
OB	Observation Block
PSO	Paranal Science Operations
QC	Quality Control
RON	Read Out Noise
UVES	Ultraviolet Visual Echelle Spectrograph
SOF	Set Of Frames
UT	Unit Telescope
VLT	Very Large Telescope

## Index

ANSI, [69](#)

ASCII, [69](#)

Background extraction

    Description, [60](#)

CalibDB, [69](#)

CPL, [69](#)

    web page, [9](#)

Data reduction cascade, [41](#), [43](#)

Definition

    ANSI, [69](#)

    ASCII, [69](#)

    CalibDB, [69](#)

    CPL, [69](#)

    DFS, [69](#)

    DMD, [69](#)

    DRS, [69](#)

    ESOREX, [69](#)

    FITS, [69](#)

    FOV, [69](#)

    FPN, [69](#)

    GUI, [69](#)

    OB, [69](#)

    PSO, [69](#)

    QC, [69](#)

    RON, [69](#)

    SOF, [69](#)

    UT, [69](#)

    UVES, [69](#)

    VLT, [69](#)

Description

    Background extraction, [60](#)

    Efficiency determination, [64](#)

    First guess solution, [56](#)

    Instrument response determination, [62](#)

DFS, [69](#)

DMD, [69](#)

DRS, [69](#)

Efficiency determination, [64](#)

ESO, [69](#)

ESOREX, [69](#)

EsoRex

    Configuration file, [20](#)

    Example

        uves\_cal\_mkmaster, [23](#)

        uves\_cal\_orderpos, [22](#)

        uves\_cal\_predict, [22](#)

        uves\_cal\_response, [25](#)

        uves\_cal\_wavecal, [23](#)

        uves\_obs\_scired, [25](#)

    Help, [19](#)

    Introduction, [19](#)

    Parameters listing, [20](#)

    Pipeline recipe configuration, [20](#)

    Pipeline recipe execution, [20](#)

    Pipeline recipe help, [20](#)

    Pipeline recipe man page, [20](#)

    Pipeline recipes listing, [20](#)

    Set of frames definition, [19](#)

    Syntax, [19](#)

    Usage, [19](#)

    web page, [9](#)

FITS, [69](#)

FOV, [69](#)

FPN, [69](#)

Gasgano

    File classification, [15](#)

    File selection, [15](#)

    Introduction, [14](#)

    Pipeline Recipe execution, [15](#)

    Usage, [15](#)

    web page, [9](#)

GUI, [69](#)

IFS, [41](#)

Input data, [42](#)

Instrument response determination, [62](#)

Integral Field Spectroscopy, [41](#)

Master frame generation, [60](#)

OB, [69](#)

<b>ESO</b>	<b>FLAMES-UVES Pipeline User Manual</b>	Doc:	VLT-MAN-ESO-19500-3016
		Issue:	Issue 5.0
		Date:	Date 2008-11-28
		Page:	110 of 110

Order Definition, [59](#)

Physical Model

Diagnostic plots, [56](#)

Physical Model, [56](#)

Parameter corrections, [57](#)

Pipeline

Cookbook, [14](#)

Example

uves\_cal\_mkmaster, [23](#)

uves\_cal\_orderpos, [22](#)

uves\_cal\_predict, [22](#)

uves\_cal\_response, [25](#)

uves\_cal\_wavecal, [23](#)

uves\_obs\_scired, [25](#)

Input data, [42](#)

Overview, [9](#), [11](#)

Purposes, [11](#)

Recent improvements, [12](#)

Recipe execution, [20](#)

Recipes

sinfo\_rec\_detlin, [43](#)

sinfo\_rec\_distortion, [44](#)

sinfo\_rec\_jitter, [44](#)

sinfo\_rec\_mdark, [44](#)

sinfo\_rec\_mflat, [44](#)

sinfo\_rec\_wavecal, [44](#)

uves\_cal\_mbias, [14](#)

uves\_cal\_mdark, [14](#)

uves\_cal\_mflat, [14](#)

uves\_cal\_mkmaster, [14](#), [48](#)

uves\_cal\_mkmaster\_sflat, [14](#)

uves\_cal\_orderpos, [14](#), [47](#)

uves\_cal\_predict, [14](#), [46](#)

uves\_cal\_response, [14](#), [51](#)

uves\_cal\_tflat, [14](#)

uves\_cal\_wavecal, [14](#), [50](#)

uves\_obs\_cd\_align, [14](#)

uves\_obs\_redchain, [14](#)

uves\_obs\_scired, [14](#), [53](#)

uves\_obs\_spatred, [14](#)

PSO, [69](#)

QC, [69](#)

QFITS

web page, [9](#)

Raw data, [28](#)

Arc Lamp Frames, [37](#)

Bias frames, [38](#)

Flat frames, [38](#)

Formatcheck, [37](#)

Order Frame, [37](#)

Science frames, [39](#)

STD star frames, [39](#)

Raw Frames

Characteristic FITS keywords, [30](#)

Supported, [30](#)

Reference data

Atmospheric extinction table, [40](#)

Line reference table, [40](#)

Standard stars flux table, [40](#)

RON, [69](#)

Set of frames

Example

uves\_cal\_mkmaster, [21](#), [22](#)

uves\_cal\_orderpos, [21](#)

uves\_cal\_predict, [21](#)

uves\_cal\_response, [22](#)

uves\_cal\_wavecal, [21](#)

uves\_obs\_scired, [22](#)

SOF, [69](#)

UT, [69](#)

UVES, [13](#), [69](#)

Instrument, [13](#)

Instrument description, [13](#)

Pipeline publication, [9](#)

Pipeline User Manual, [9](#)

Quality Control web page, [9](#)

Raw data, [28](#)

Science Verifications web page, [9](#)

User Manual, [9](#)

web page, [9](#)

uves\_cal\_predict, [56](#)

VLT, [69](#)

Wavelength calibration, [61](#)

Diagnostic plots

Residuals, [62](#)

Resolution, [62](#)